

Star^{HYGIENE} DENTALEZ

BETTER HYGIENE MADE EASY WITH iSTAR™ PROPHY HANDPIECE

Our iStar prophy features a simple, compact, cordless design that dramatically simplifies the prophylaxis procedure. No more stress from the pull of handpiece tubing. The iStar is ideal for the modern dental operator and features thoughtful design elements that make work more efficient and comfortable.

Sealed system with custom disposable barriers facilitate disinfection between patients

Five speed options from 500-2,500 rpm allow speed adjustment based on user preference

On/Off finger tip control for intuitive operation

Compatible with most Disposable Prophy Angles

Part number: 265993



PROPHY STAR® 3

Constructed of anodized aluminum and weighing just 62 grams, the lightweight Prophy Star 3 helps minimize operator fatigue. The 360-degree swivel provides better control for greater precision and comfort.



Air-driven control up to 5,000 rpm

LubeFree design is maintenance free which saves time and money

Use with your preferred disposable angle (available separately)


Part number: 264422

Specifications

Speed	0-5,000 rpm adjustable
Air Pressure	35-40 psi
Dimensions	0.64" in diameter
Weight	62 grams
Length	4.3"
Connection	4-Line Configuration

TITAN® BLIS-SONIC™ SCALER

Sonic scaler to effectively and gently remove hard calculus deposits and stains. Conveniently connects to existing air-driven tubings and foot control. Ergonomic silicone grips with comfortable positioning reduces hand fatigue.

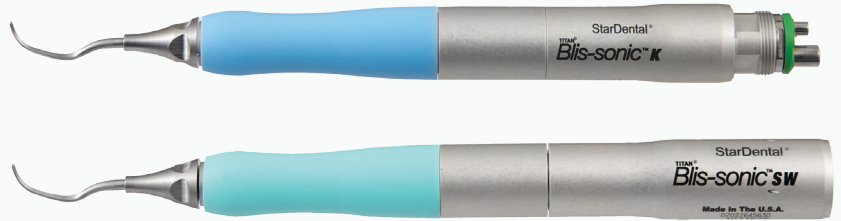

Swivel (SW) and Fixed (K) backend configurations

Lightweight design

Includes five scaling tip designs

SW part number: 264563

K part number: 264565

Star® Nomenclature Key

SW = Swivel (coupler) back end

K = Fixed Back End



Universal

- Use for all Surfaces, medium and heavy calculus



Sickle

- Recommended for sensitive patients
- Less powerful than the Universal Tip
- Less suitable for deep root scaling



Perio

- Recommended for lingual supragingival scaling of mandibular incisors and deep local pockets
- Provides same energy as Universal Tip
- Also includes left and right perio tips. (not shown)

Blis-sonic SW Specifications

Air Pressure	40 psi
Weight	54.6 grams
Length	6.8 inches (17.3cm)

Blis-sonic K Specifications

Air Pressure	40 psi
Weight	50.7 grams
Length	6.7 inches (17.0cm)