



DENTALEZ

EQUIPMENT

2024


TABLE OF CONTENTS

EQUIPMENT

Pivot Chair Mount Packages	8	Rear Mount Side Delivery Systems	51
Standard IC Control Head	8	Cabinetry Mounting Configurations	53
Euro Control Head	9	7020PRO Duo Swing Delivery Systems	54
Fixed Chair Mount with Sidebox Packages	10	7021PRO Duo Swing Delivery Systems	56
Standard IC Control Head	10	7010 Duo Swing Delivery Systems	58
Euro Control Head	11	3486 Swing Delivery Systems	60
3900 Package Ordering Guide	12	Storage Organizer, Trays and Monitor Mount	62
Fusion Packages	24	Ancillary Ordering and Options	63
Pivot Chair Mount	25	Standard IC Delivery Unit	65
Fixed Mount with Sidebox	26	Euro Delivery Unit	66
Rear Mount Side Delivery	27	Fusion Delivery Unit	67
Fusion Package Ordering Guide	28	Handpiece Bundles	68-69
Surface Treatments	32	Control Head Accessories	70
Designer Friendly Colors and Surfaces	33	Ancillary Products	70
Chairs and Stools	35	Simplified Cuspidors and Assistant's Instrumentation	73
3900 Hydraulic Patient Chair	37	Fixed Mount with Sidebox, External Umbilical	76
3900 Chair Upholstery Options	38	Fixed Chair Post Mount, External Umbilical	80
3800 Hydraulic Patient Chair	39	Rear Beneath Shelf or Work Surface Mount	81
Chair Accessories	40	Simplified Rear and Side Delivery Systems	
Economy and Deluxe Stools	41	Side Cabinet and Wall Mount	82
Comfort Stools	42	Rear or Side Beneath Counter Mount	83
Lights and Monitor Mounts	43	Rear Beneath Counter Mount	84
Post Mount Lights	45	Panel Mounts	86
Track Mount Lights	46	Rear or Side Mount Compact Controls	87
Ceiling Mount Lights	47	Mobile Carts	88
Cabinet and Wall Mount Lights	48	Factory Options and Parts	89
Monitor Mount Assemblies	49	General Information	95
		Dimensions	95
		Equipment Maintenance	110
		Equipment Terms and Warranty	111

Ramvac

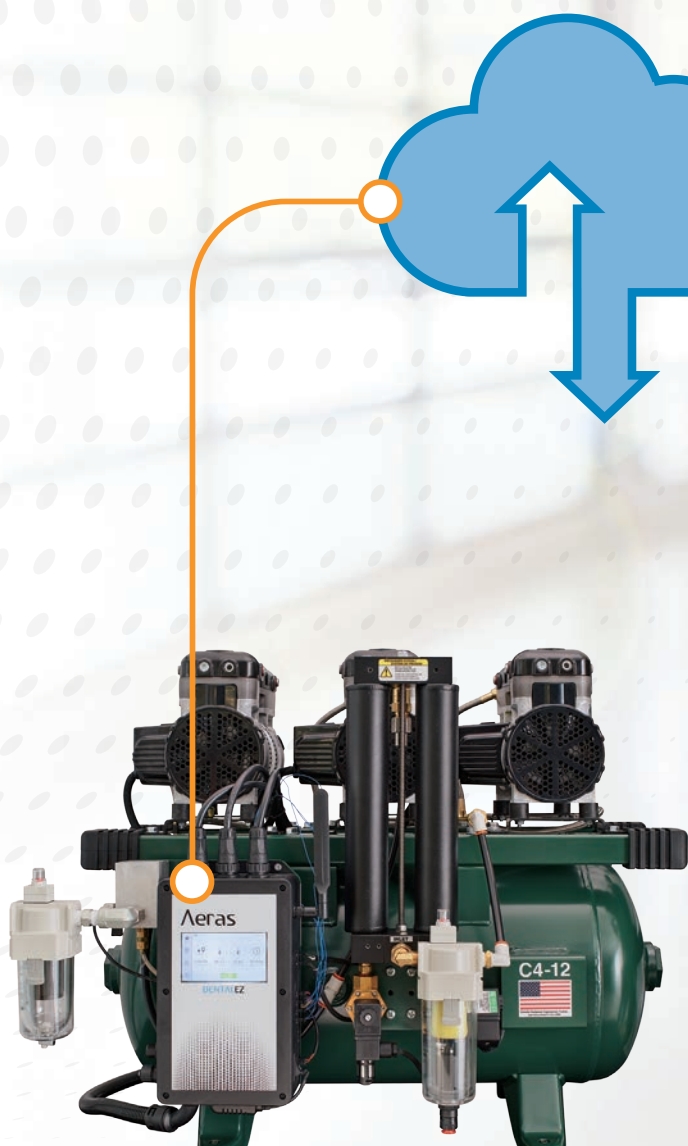DENTALEZ



RAMVAC® WITH AERAS INTEL

In an increasingly complex world, smart devices simplify how we integrate new technology into everyday life. From smartphones to smart homes, intelligent devices are simplifying the way we live. And now, they're doing some of the work for us at work too. With Aeras Intel, the smart revolution has reached the dental office, offering interconnectivity, insight, and control.

AERAS, Inc. is a separate corporation from DENTALEZ, Inc.

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Forest

PACKAGES



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Standard Chair Color Options

Powder Coat Paint Forest Grey

Shown with Diplomat Blue Upholstery

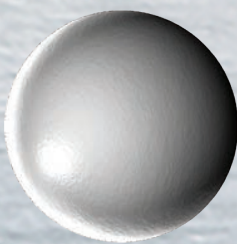


Powder Coat Paint Shadow


Shown with Raven Wing Upholstery




Designer Friendly Chair Color Options




Arctic Silver

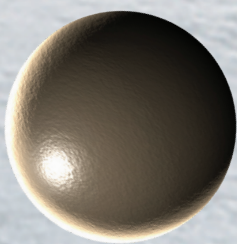


Bronze

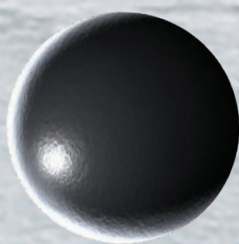


Charcoal

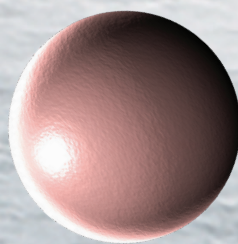


Copper

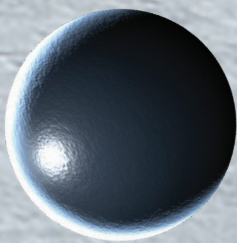

Midnight

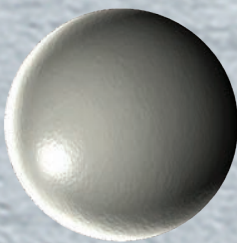


Pewter

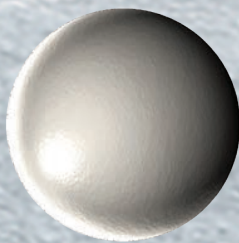

White


Pivot Chair Mount

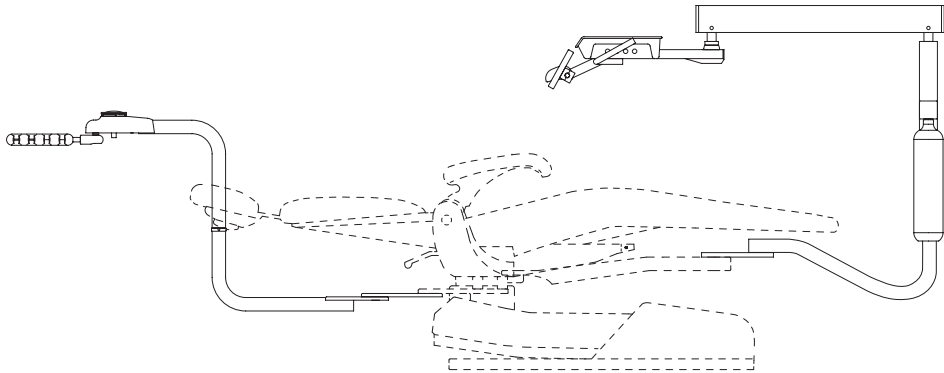Pivot Chair Mounts easily move left to right on large roller bearings to accommodate left and right-handed practitioners. Unlike a Fixed Chair Mount delivery system, the mounting arm comes from underneath the patient rather than over the top, making it less obtrusive and more patient-friendly. Pivot Assistant's Instrumentation, Operatory Lights and Monitor Mounts are popular options to support your style of practice.



Fixed Chair Mount With Sidebox

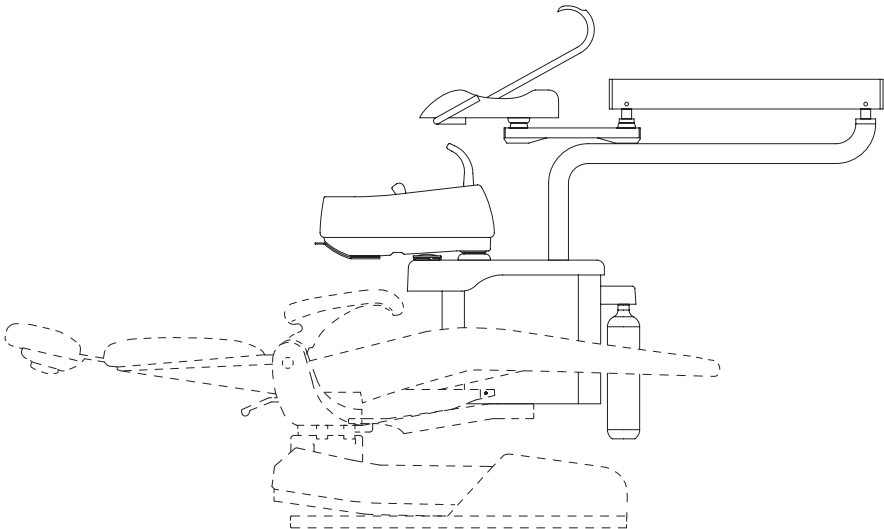Fixed Chair Mounts With Sidebox allow for a Cuspidor, Assistant's Instrumentation, a Monitor and Light to mount on a sidebox, which is large enough to accommodate a host of ancillary support products. Finish the aluminum side covers in one of our nine powder coat paint colors and choose a complementary carbon fiber trim to finish the top cap. If you use a cuspidor, select our beautiful high-gloss ceramic black for a customized and beautiful patient-friendly delivery system.



Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Base Package



Note: Package shown in shadow paint with raven wing ultraleather.

Popular Options



Dr's Deluxe Stool,
Ultraleather Upholstery

Part number: 6159C



E- Torque Electric
Handpiece System

Part number:
STAR-ETORQUE



LED Operatory Light, 9073

Part number: 9073



Assistants Instrumentation
rear Pivot with touchpad

Part number: 1096

Base Package As Shown

Description	Part Number
3900 Patient Chair	3900
Delivery System w/Standard IC, Pivot	4486PI
Assistants Instrumentation, Rear Pivot	1095
LED Operatory Light, 9078	9078
Light Post, Pivot	1113-654
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Popular Options

Description	Part Number
Dr. Deluxe Stool, Ultraleather Upholstery	6159C
Assistants Deluxe Stool No Back, Ultraleather Upholstery	6160C
E-Torque Electric HP System	STAR-ETORQUE
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Quick Disconnect for Water	0016-454-CS
Quick Disconnect for Air	0016-456-CS
Routing Valve Kit	1115-315-CS
Additional HVE w/Dry Oral Cup	1400-113-CS
Soft Touch Back Massage	3900-100
LED Operatory Light, 9073	9073
Assistants Instrumentation, Rear Pivot with Touchpad	1096

Equipment dimensions, see page 95
Shipping weight estimate:
Delivery unit & chair 450 lb / 204 kg
Base package 580 lb / 264 kg
As shown package 630 lb / 286 kg



Base Package



Popular Options



Dr's Deluxe Stool,
Ultraleather Upholstery

Part number: 6159C



E- Torque Electric
Handpiece System

Part number:
STAR-ETORQUE



LED Operatory Light, 9073

Part number: 9073



Assistants Instrumentation
rear Pivot with touchpad

Part number: 1096

Base Package As Shown

Description	Part Number
3900 Patient Chair	3900
Delivery System w/Euro, Pivot	4988PI
Assistants Instrumentation, Rear Pivot	1095
LED Operatory Light, 9078	9078
Light Post, Pivot	1113-654
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Popular Options

Description	Part Number
Dr. Deluxe Stool, Ultraleather Upholstery	6159C
Assistants Deluxe Stool No Back, Ultraleather Upholstery	6160C
E-Torque Electric HP System	STAR-ETORQUE
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Quick Disconnect for Water	0016-454-CS
Quick Disconnect for Air	0016-456-CS
Routing Valve Kit	1115-315-CS
Additional HVE w/Dry Oral Cup	1400-113-CS
Soft Touch Back Massage	3900-100
LED Operatory Light, 9073	9073
Assistants Instrumentation, Rear Pivot with Touchpad	1096

Equipment dimensions, see page 95
Shipping weight estimate:
Delivery unit & chair 445 lb / 203 kg
Base package 533 lb / 264 kg
As shown package 630 lb / 286 kg



Base Package



Popular Options



Dr's Deluxe Stool,
Ultraleather Upholstery
Part number: 6159C



E- Torque Electric
Handpiece System
Part number:
STAR-ETORQUE



Swivel Tray
Part number: 1110-468



LED Operatory Light, 9073
Part number: 9073

Base Package As Shown

Description	Part Number
3900 Patient Chair	3900
Delivery System, Sidebox w/Standard IC	4485SI
Routing Valve	Included
Quick Disconnect for Water	Included
Quick Disconnect for Air	Included
Deluxe Cuspidor	1577
LED Operatory Light, 9078	9078
Light Post, Sidebox	1110-900
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Popular Options

Description	Part Number
Dr. Deluxe Stool, Ultraleather Upholstery	6159C
Assistants Deluxe Stool No Back, Ultraleather Upholstery	6160C
E-Torque Electric HP System	STAR-ETORQUE
Satelec SP	SATELEC-S
Satelec SP, LED lit	SATELEC-S-LED
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Swivel Tray	1110-468
Additional Handpiece Position, Wet	1100-189-CS
Deluxe Cuspidor Remote Activation	1100-480-CS
LED Operatory Light, 9073	9073
Soft Touch Back Massage	3900-100

Equipment dimensions, see page 96
Shipping weight estimate:
Delivery unit & chair 440 lb / 202 kg
Base package 556 lb / 253 kg
As shown package 561 lb / 255 kg



Base Package



Popular Options



Dr's Deluxe Stool,
Ultraleather Upholstery
Part number: 6159C



Assistants Deluxe Stool,
No Back,
Ultraleather Upholstery
Part number: 6160C



E- Torque Electric
Handpiece System
Part number:
STAR-ETORQUE



LED Operator Light, 9073
Part number: 9073

Base Package As Shown

Description	Part Number
3900 Patient Chair	3900
Delivery System, Sidebox w/Euro	4987SI
Routing Valve	Included
Quick Disconnect for Water	Included
Quick Disconnect for Air	Included
Deluxe Cuspidor	1577
LED Operator Light, 9078	9078
Light Post, Sidebox	1110-900
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Popular Options

Description	Part Number
Dr. Deluxe Stool, Ultraleather Upholstery	6159C
Assistants Deluxe Stool No Back, Ultraleather Upholstery	6160C
E-Torque Electric HP System	STAR-ETORQUE
Satelec SP	SATELEC-S
Satelec SP, LED lit	SATELEC-S-LED
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Additional Handpiece Position, Wet	1100-189-CS
Deluxe Cuspidor Remote Activation	1100-480-CS
LED Operator Light, 9073	9073
Soft Touch Back Massage	3900-100

Equipment dimensions, see page 96
Shipping weight estimate:
Delivery unit & chair 445 lb / 203 kg
Base package 551 lb / 251 kg
As shown package 551 lb / 251 kg





SURFACE FINISH OPTIONS

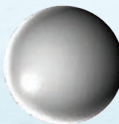
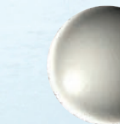
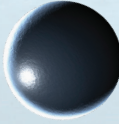

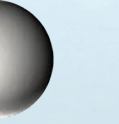
ULTRALEATHER® 12 CLASSIC COLORS

					
Autumn Leaf	Baltic	Boysenberry	Charcoal	Diplomat Blue	Doe
					
Dove Grey	Papyrus	Parrot	Raven Wing	St. John	Stone




ULTRALEATHER 5 PEARLIZED COLORS

		
Midnight	Oz	Pewter
		
Pharoah	Tin Man	

POWDER COAT PAINT 9 COLORS

				
Arctic Silver	Forest Grey	Charcoal	Copper	Bronze
				
Midnight	Pewter	Shadow	White	

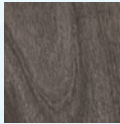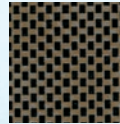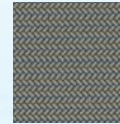CARBON FIBER 3 PATTERNS

		
Greywood	Laser Leather	Lemans

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

DESIGN YOUR 3900 PATIENT CHAIR

STEP 1: SELECT YOUR CHAIR - Select 1 chair and 1 powder coat paint color

STEP
1

Description	Part #	Powder Coat Pain Color	
<input type="checkbox"/> 3900 Patient Chair	3900	<input type="checkbox"/> Forest Grey*	<input type="checkbox"/> Shadow*
<input type="checkbox"/> 3900 Patient Chair, Designer Friendly	3900C	<input type="checkbox"/> Arctic Silver	<input type="checkbox"/> Bronze
*White, Shadow, and Forest Grey chairs come standard with the baseplate painted Shadow, all other colors come standard with the baseplate painted the same color as the chair.		<input type="checkbox"/> Charcoal	<input type="checkbox"/> Copper
		<input type="checkbox"/> Midnight	<input type="checkbox"/> Pewter
		<input type="checkbox"/> *White	

STEP 2: SELECT YOUR UPHOLSTERY - Select 1 option and 1 upholstery color

STEP
2

Description	Part #	Upholstery Color (used as inner color if dual color upholstery is selected)	
<input type="checkbox"/> Wide Back, Plush	3907-PL-WB	Standard Ultraleather® Colors	
<input type="checkbox"/> Wide Back, Seamless	3907-AS-WB	<input type="checkbox"/> Autumn Leaf	<input type="checkbox"/> Baltic
<input type="checkbox"/> Wide Back, Memory Comfort	3907-MC-WB	<input type="checkbox"/> Boysenberry	<input type="checkbox"/> Charcoal
<input type="checkbox"/> Narrow Back, Plush	3907-PL-NB	<input type="checkbox"/> Doe	<input type="checkbox"/> Diplomat Blue
<input type="checkbox"/> Narrow Back, Seamless	3907-AS-NB	<input type="checkbox"/> Dove Grey	<input type="checkbox"/> Papyrus
<input type="checkbox"/> Narrow Back with Slings, Plush	3907-PL-NS	<input type="checkbox"/> Parrot	<input type="checkbox"/> Raven Wing
<input type="checkbox"/> Narrow Back with Slings, Seamless	3907-AS-NS	<input type="checkbox"/> St. John	<input type="checkbox"/> Stone
		<input type="checkbox"/> Midnight*	<input type="checkbox"/> Pewter*
		<input type="checkbox"/> Tin Man*	<input type="checkbox"/> Pharoah*
		<input type="checkbox"/> Oz*	
		<input type="checkbox"/> *Pearlized	
		<input type="checkbox"/> Special color, select PN 3907-999-SP from Step 3 and write color choice here: _____	

STEP 3: ADDITIONAL UPHOLSTERY OPTIONS - Select at least 1 of the below options

STEP
3

<input type="checkbox"/> No Additional Options	N/A	
<input type="checkbox"/> Special Color Ultraleather Upholstery	3907-999-SP	Color: _____
<input type="checkbox"/> Dual Color, Plush - specify outer color (Not available with seamless or memory comfort)	3907-999-DC	Outer Color: (inner color selected in Step 2 above)
<input type="checkbox"/> Dual Color for Memory Comfort specify outer color	3907-MC-DC	Outer Color: (inner color selected in Step 2 above)
<input type="checkbox"/> Contrast/Deviated Stitching, Plush (Not available with seamless or memory comfort)	3907-999-DS	Thread Color: _____
<input type="checkbox"/> Neck Pillow	3907-318	Color: _____
<input type="checkbox"/> Knee Pillow	3907-317	Color (inner if dual color): _____
<input type="checkbox"/> Dual Color Upholstery, Knee Pillow	3907-317-DC	Outer Color: _____
<input type="checkbox"/> Contrast Stitch Upholstery, Knee Pillow	3907-317-DS	Thread Color: _____
<input type="checkbox"/> Pedo Booster Set	3907-180	
<input type="checkbox"/> 29" Toe Guard in Lieu of 21"	3907-028-CS	

STEP 4: MASSAGE OPTION - Select at least 1 of the below options

STEP
4

<input type="checkbox"/> No Massage Options	N/A	
<input type="checkbox"/> Soft Touch Massage*	3900-100	Not available with Memory Comfort Upholstery

*If a chair-mounted delivery unit is being ordered then a utility center housing is included with purchase. Otherwise we recommend adding a utility center for the massage air compressor (Part numbers below).

STEP 5: UTILITY CENTER SELECTION - Select at least 1 of the below options

STEP
5

<input type="checkbox"/> No Additional Options	N/A	Select this option if you are buying a chair-mount delivery unit and no massager
Contoured Utility Center Housing and Cover		
<input type="checkbox"/> Forest Grey	3900-101-G	
<input type="checkbox"/> Shadow	3900-101-SH	
<input type="checkbox"/> White	3900-101	
<input type="checkbox"/> Designer Friendly®	3900-101-C	

STEP 6: CHAIR CONTROL LOCATIONIN - Select at least 1 chair control

Maximum of 1 control can be mounted on the delivery unit, maximum of 2 on the chair (left & right), 4 total maximum number of controls per operatory

STEP
6

<input type="checkbox"/> First Chair Control (Required) First control included with chair purchase	3914-998-CS	<input type="checkbox"/> Delivery Unit	<input type="checkbox"/> Footswitch
		<input type="checkbox"/> Chair, right side	<input type="checkbox"/> Chair, left side
<input type="checkbox"/> Second Chair Control (Optional)	3914-999-CS	<input type="checkbox"/> Delivery Unit	<input type="checkbox"/> Footswitch
		<input type="checkbox"/> Chair, right side	<input type="checkbox"/> Chair, left side

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

PATIENT CHAIR

STEP

1

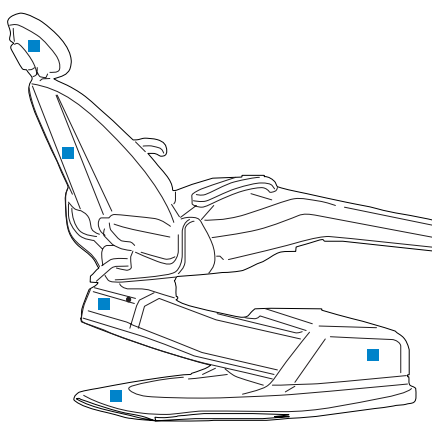
Chair paint color



STEP

2


Upholstery color (inner color if dual color selected)



STEP

3

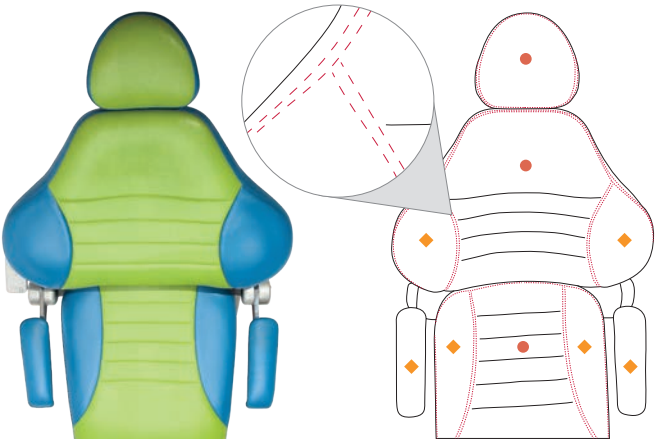Outer upholstery color (if dual color selected, PN 3907-99-DC)
Contrast stitch color (if color stitch selected, PN 3907-999-DS) -----



CHAIR IMAGES



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information



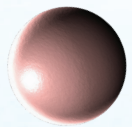
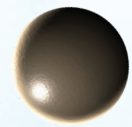
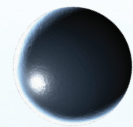
Warranty

DESIGN YOUR STANDARD IC DELIVERY SYSTEM






SURFACE FINISH OPTIONS

POWDER COAT PAINT 9 COLORS

		
Arctic Silver	Forest Grey	Charcoal
		
Copper	Bronze	Midnight
		
Pewter	Shadow	White

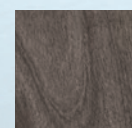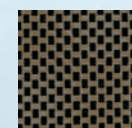CARBON FIBER 3 PATTERNS

		
Greywood	Laser Leather	Lemans

[TABLE OF CONTENTS](#)

- [Pivot & Fixed Chair Mount Packages](#)
- [Pivot & Fixed Chair Mount Package Ordering Guide](#)
- [Fusion Packages](#)
- [Fusion Package Ordering Guide](#)
- [Designer Friendly Colors & Surfaces](#)
- [Chairs & Stools](#)
- [Lights and Monitor Mounts](#)
- [Rear Mount Side Delivery Systems](#)
- [Ancillary Ordering & Options](#)
- [Control Head Accessories](#)
- [Simplified Cuspidors & Assistant's Instrumentation](#)
- [Simplified Rear & Side Delivery Systems](#)
- [Factory Options & Parts](#)
- [General Information](#)
- [Warranty](#)

DESIGN YOUR STANDARD IC DELIVERY SYSTEM

STEP
7

STEP 7: SELECT YOUR DELIVERY SYSTEM - Select 1 delivery system

Description	Part #
<input type="checkbox"/> Standard IC Delivery Head, Pivot Chair Mount	4486PI
<input type="checkbox"/> Standard IC Delivery Head, Fixed Chair Mount with Sidebox	4485SI

STEP
8

STEP 8: SURFACE TREATMENTS - Select 1 powder coat and 1 plastic finish

Powder Coat Paint:	Plastic Finish:
<input type="checkbox"/> Forest Grey <input type="checkbox"/> Shadow <input type="checkbox"/> Arctic Silver <input type="checkbox"/> Bronze <input type="checkbox"/> Charcoal <input type="checkbox"/> Copper <input type="checkbox"/> Midnight <input type="checkbox"/> Pewter <input type="checkbox"/> White	<input type="checkbox"/> Forest Grey <input type="checkbox"/> Greywood <input type="checkbox"/> Laser Leather <input type="checkbox"/> Lemans <input type="checkbox"/> White

STEP
9

STEP 9: CONFIGURE YOUR DELIVERY SYSTEM

Note: Please refer to ancillary and options graph in the catalog to determine number of ancillaries and options possible

Description	Part #	Handpiece Positions / Comments						
		0	1	2	3	4	5	6
<input type="checkbox"/> Air/Water Syringe (included standard)								
<input type="checkbox"/> 4-hole handpiece tubing (included standard)								
<input type="checkbox"/> 5-hole system (1 tubing included, system accomodates a 2nd tubing)	1700-021-CS							
<input type="checkbox"/> - Second 5-hole tubing	0017-615-CS							
<input type="checkbox"/> 6-pin system (1 tubing included, system accomodates a 2nd tubing)	1700-019-CS							
<input type="checkbox"/> - Second 6-pin tubing	0017-603-CS							
<input type="checkbox"/> Star® ETorque Electric System, Factory Installed*	STAR-ETORQUE							
<input type="checkbox"/> Star Blis-sonic™ SW Scaler (with 4 line swivel)*	STAR-BLIS-SW							
<input type="checkbox"/> Star Blis-sonic K Scaler (with fixed back-end)*	STAR-BLIS-K							
<input type="checkbox"/> Satelec SP Scaler*	SATELEC-S							
<input type="checkbox"/> Satelec SP LED Scaler*	SATELEC-S-LED							
<input type="checkbox"/> Cavitron Scaler*	CAVITRON							
<input type="checkbox"/> Additional Handpiece Position, Wet	1100-189-CS							
<input type="checkbox"/> Additional Handpiece Position, Air Drive Only	1100-190-CS							
<input type="checkbox"/> Assistant's hygiene instrumentation, control head-mounted Select 1 position for the SE and 1 position for the HVE	1408-CS	SE						
		HVE						
<input type="checkbox"/> Ancillary Shelf	0010-310							
<input type="checkbox"/> Star Handpiece Bundle (Write-in PN and price)								
<input type="checkbox"/> Quick Disconnect w/Flow Control, 1/4" (for water)	0016-454-CS							
<input type="checkbox"/> Quick Disconnect w/Flow Control, 3/8" (for air)	0016-456-CS							
<input type="checkbox"/> Swivel Tray	1110-468							
<input type="checkbox"/> Less sidebox telescoping assistant's instrumentation	1400-000-CS							
<input type="checkbox"/> Chrome Accents	9000-998-CS							
<input type="checkbox"/> Additional Umbilical, Per Foot (max of 9 feet)	1115-999-CS							
<input type="checkbox"/> Routing Valve Kit w/Water shutoff & Regulator, Factory inst.	1115-315-CS							
<input type="checkbox"/> Routing Valve Kit for Cuspidor, Factory Installed	1115-316-CS							

* Ancillaries utilize an existing handpiece position - some integrated options may require purchase of Ancillary Shelf PN 0010-310

STEP
10

STEP 10: ASSISTANT'S INSTRUMENTATION AND CUSPIDORS - Must select 1 option for order to be valid

<input type="checkbox"/> No Assistant's Instrumentation or Cuspidor*	N/A							
<input type="checkbox"/> Rear Pivot Instrumentation, Quad Holder Bar	1095							
<input type="checkbox"/> Rear Pivot Instrumentation with Integrated Touchpad	1096							
<input type="checkbox"/> Rear Pivot Instrumentation with Deluxe Cuspidor**	1590	Bowl Color	<input type="checkbox"/> Black	<input type="checkbox"/> White				
<input type="checkbox"/> Deluxe Cuspidor only, Rear Pivot**	1580	Bowl Color	<input type="checkbox"/> Black	<input type="checkbox"/> White				
<input type="checkbox"/> Basic Cuspidor, Sidebox Mounted (4485SI only)	1679							
<input type="checkbox"/> Deluxe Cuspidor, Sidebox Mounted (4485SI only)**	1577	Bowl Color	<input type="checkbox"/> Black	<input type="checkbox"/> White				

*Sidebox mounted telescoping assistant's instrumentation included standard with 4485SI

**Not available with 4485SI

STEP
11

STEP 11: ASSISTANT'S INSTRUMENTATION AND CUSPIDOR OPTIONS

<input type="checkbox"/> Quick Disconnect w/Flow Control, 1/4" (For 1095 instrumentation)	0016-455-CS	Specify One	<input type="checkbox"/> Air	<input type="checkbox"/> Water
<input type="checkbox"/> Additional Autoclavable HVE with 5' Tubing & Dry Oral Cup	1400-113-CS			
<input type="checkbox"/> Hands-Free HVE Holder with HVE and 5' Tubing (Select PN) Colors (holder/tubing): GRY (white/grey), STG (white/sterling), BLK (black/black), SHW (black/shadow)		<input type="checkbox"/> 1400-279-STG	<input type="checkbox"/> 1400-279-BLK	
		<input type="checkbox"/> 1400-279-SHW	<input type="checkbox"/> 1400-279-GRY	
<input type="checkbox"/> Self-Contained Clean Water System w/Quick Disconnect, 1.5L bottle (For 1095 Instrumentation only)	1700-319			
<input type="checkbox"/> Tray Assembly (For 1095 Instrumentation only)	1110-680			
<input type="checkbox"/> Deluxe Cuspidor Unit-Mounted Activation (For 1577, 1580, 1590)	1100-480-CS			
<input type="checkbox"/> Chair Control Touchpad - Mounted to 1095 Instrumentation	1700-088-CS			

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

STANDARD IC DELIVERY UNIT

STEP

8

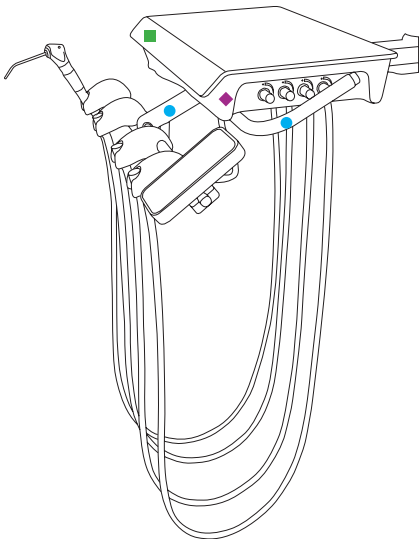
Powder coat paint color



STEP

8


Powder coat paint with upgrade option to chrome/polished finish (PN 9000-998-CS, Step 9)



STEP

8

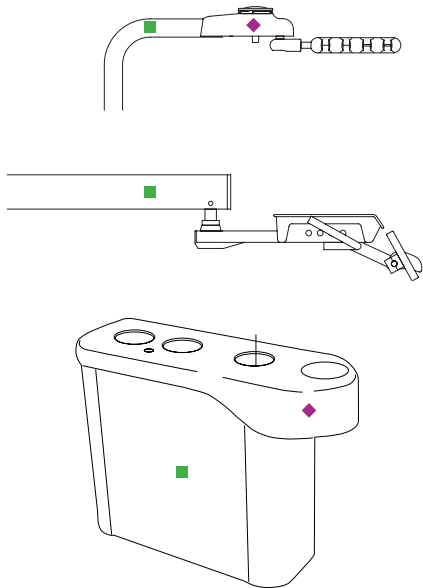Plastic finish



UNIT IMAGES



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

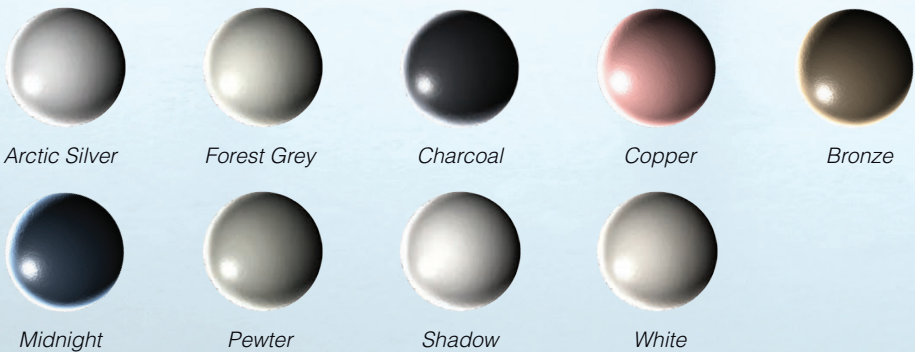
Warranty

DESIGN YOUR EURO DELIVERY SYSTEM




SURFACE FINISH OPTIONS

POWDER COAT PAINT 9 COLORS



CARBON FIBER 3 PATTERNS

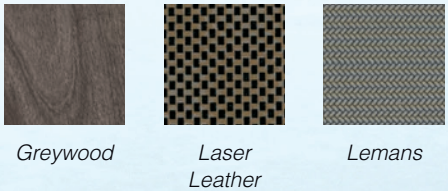

- Pivot & Fixed Chair Mount Packages
- Pivot & Fixed Chair Mount Package Ordering Guide
- Fusion Packages
- Fusion Package Ordering Guide
- Designer Friendly Colors & Surfaces
- Chairs & Stools
- Lights and Monitor Mounts
- Rear Mount Side Delivery Systems
- Ancillary Ordering & Options
- Control Head Accessories
- Simplified Cuspidors & Assistant's Instrumentation
- Simplified Rear & Side Delivery Systems
- Factory Options & Parts
- General Information
- Warranty

DESIGN YOUR EURO DELIVERY SYSTEM

STEP
7

STEP 7: SELECT YOUR DELIVERY SYSTEM - Select 1 delivery system

Description	Part #	Comments
<input type="checkbox"/> Euro Delivery Head, Pivot Chair Mount	4988PI	
<input type="checkbox"/> Euro Delivery Head, Fixed with Sidebox	4987SI	

STEP
8

STEP 8: SURFACE TREATMENTS - Select 1 powder coat and 1 plastic finish

Powder Coat Paint:	Plastic Finish:
<input type="checkbox"/> Forest Grey <input type="checkbox"/> Shadow <input type="checkbox"/> Arctic Silver <input type="checkbox"/> Bronze <input type="checkbox"/> Charcoal	<input type="checkbox"/> Forest Grey <input type="checkbox"/> Greywood <input type="checkbox"/> Laser Leather
<input type="checkbox"/> Copper <input type="checkbox"/> Midnight <input type="checkbox"/> Pewter <input type="checkbox"/> White	<input type="checkbox"/> Lemans <input type="checkbox"/> White

STEP
9

STEP 9: CONFIGURE YOUR DELIVERY SYSTEM

Note: Please refer to ancillary and options graph in the catalog to determine number of ancillaries and options possible

Description	Part #	Handpiece Positions / Comments				
		0	1	2	3	4
<input type="checkbox"/> Air/Water Syringe (included standard)		X				
<input type="checkbox"/> 4-hole handpiece tubing (included standard)						
<input type="checkbox"/> 5-hole system (1 tubing included, system accomodates a 2nd tubing)	1700-021-CS					
<input type="checkbox"/> - Second 5-hole tubing	0017-615-CS					
<input type="checkbox"/> 6-pin system (1 tubing included, system accomodates a 2nd tubing)	1700-019-CS					
<input type="checkbox"/> - Second 6-pin tubing	0017-603-CS					
<input type="checkbox"/> Star® ETorque Electric System, Factory Installed	STAR-ETORQUE					
<input type="checkbox"/> Star Blis-sonic™ SW Scaler (with 4 line swivel)	STAR-BLIS-SW					
<input type="checkbox"/> Star Blis-sonic K Scaler (with fixed back-end)	STAR-BLIS-K					
<input type="checkbox"/> Satelec SP Scaler	SATELEC-S					
<input type="checkbox"/> Satelec SP LED Scaler	SATELEC-S-LED					
<input type="checkbox"/> Cavitron Scaler	CAVITRON					
<input type="checkbox"/> Star Handpiece Bundle (Write-in PN and price)						
<input type="checkbox"/> Quick Disconnect w/Flow Control, 1/4" (for water)	0016-454-CS	Included with 4987SI				
<input type="checkbox"/> Quick Disconnect w/Flow Control, 3/8" (for air)	0016-456-CS	Included with 4987SI				
<input type="checkbox"/> Less sidebox telescoping assistant's instrumentation	1400-000-CS	4987SI only				
<input type="checkbox"/> Chrome Accents	9000-998-CS					
<input type="checkbox"/> Additional Umbilical, Per Foot (max of 9 feet)	1115-999-CS					
<input type="checkbox"/> Routing Valve Kit w/Water shutoff & Regulator, Factory inst.	1115-315-CS	Included with 4987SI				
<input type="checkbox"/> Routing Valve Kit for Cuspidor, Factory Installed	1115-316-CS					

STEP
10

STEP 10: ASSISTANT'S INSTRUMENTATION AND CUSPIDORS - Must select 1 option for order to be valid

<input type="checkbox"/> No Assistant's Instrumentation or Cuspidor*	N/A			
<input type="checkbox"/> Rear Pivot Instrumentation, Quad Holder Bar	1095			
<input type="checkbox"/> Rear Pivot Instrumentation with Integrated Touchpad	1096			
<input type="checkbox"/> Rear Pivot Instrumentation with Deluxe Cuspidor**	1590	Bowl Color	<input type="checkbox"/> Black	<input type="checkbox"/> White
<input type="checkbox"/> Deluxe Cuspidor only, Rear Pivot**	1580	Bowl Color	<input type="checkbox"/> Black	<input type="checkbox"/> White
<input type="checkbox"/> Basic Cuspidor, Sidebox Mounted (4987SI only)	1679			
<input type="checkbox"/> Deluxe Cuspidor, Sidebox Mounted (4987SI only)**	1577	Bowl Color	<input type="checkbox"/> Black	<input type="checkbox"/> White
*Sidebox mounted telescoping assistant's instrumentation included standard with 4987SI		**Not available with 4987SI		

STEP
11

STEP 11: ASSISTANT'S INSTRUMENTATION OPTIONS

<input type="checkbox"/> Quick Disconnect w/Flow Control, 1/4" (For 1095 instrumentation)	0016-455-CS	Specify	<input type="checkbox"/> Air	<input type="checkbox"/> Water
<input type="checkbox"/> Additional Autoclavable HVE with 5' Tubing & Dry Oral Cup	1400-113-CS			
<input type="checkbox"/> Hands-Free HVE Holder with HVE and 5' Tubing (Select PN)		<input type="checkbox"/> 1400-279-STG	<input type="checkbox"/> 1400-279-BLK	
Colors (holder/tubing): GRY (white/grey), STG (white/sterling), BLK (black/black), SHW (black/shadow)		<input type="checkbox"/> 1400-279-SHW	<input type="checkbox"/> 1400-279-GRY	
<input type="checkbox"/> Self-Contained Clean Water System w/Quick Disconnect, 1.5L bottle (For 1095 Instrumentation only)	1700-319			
<input type="checkbox"/> Tray Assembly (For 1095 Instrumentation only)	1110-680			
<input type="checkbox"/> Deluxe Cuspidor Unit-Mounted Activation (For 1577, 1580, 1590)	1100-480-CS			
<input type="checkbox"/> Chair Control Touchpad - Mounted to 1095 Instrumentation	1700-088-CS			

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

EURO DELIVERY UNIT

STEP

8

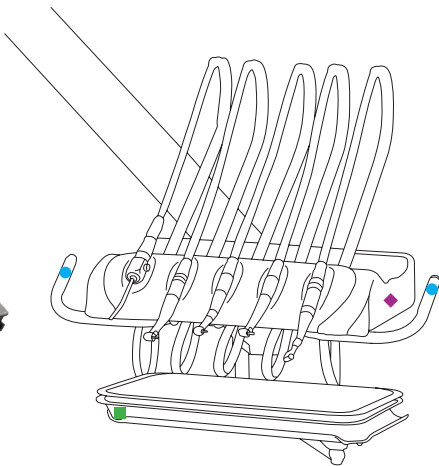
Powder coat paint color



STEP

8


Powder coat paint with upgrade option to chrome/polished finish (PN 9000-998-CS, Step 9)



STEP

8

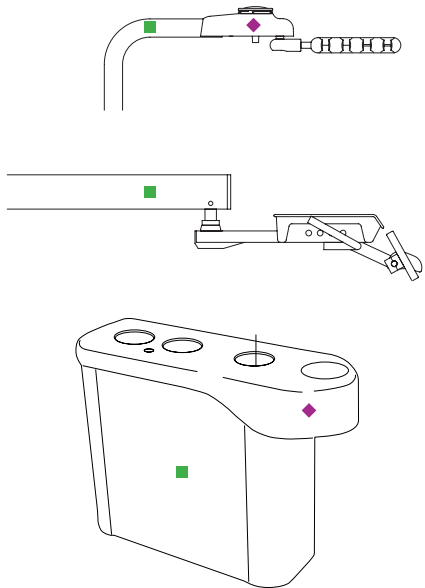Plastic finish



UNIT IMAGES



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

DESIGN YOUR LED OPERATORY LIGHT



SURFACE FINISH OPTIONS

POWDER COAT PAINT 9 COLORS

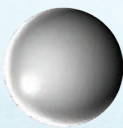
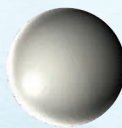
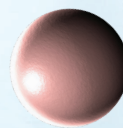
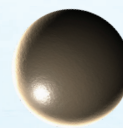
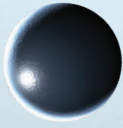
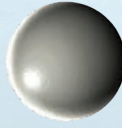
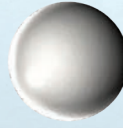
				
Arctic Silver	Forest Grey	Charcoal	Copper	Bronze
				
Midnight	Pewter	Shadow	White	

TABLE OF CONTENTS

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

DESIGN YOUR LED OPERATORY LIGHT

STEP 12: OPERATORY LIGHT - Surface Treatments, 1 powder coat paint selection required

STEP
12

Powder Coat Paint:

- ☐ Forest Grey
- ☐ Shadow
- ☐ Arctic Silver
- ☐ Bronze
- ☐ Charcoal
- ☐ Copper
- ☐ Midnight
- ☐ Pewter
- ☐ White

Select Light Head Mounting Option - 1 light and mounting option required

	9083 Halogen		9078 LED		9073 LED	
	Part #	Qty	Part #	Qty	Part #	Qty
Chair Mount	<input type="checkbox"/> 9083		<input type="checkbox"/> 9078		<input type="checkbox"/> 9073	
Ceiling Mount	<input type="checkbox"/> 9083-CL		<input type="checkbox"/> 9078-CL		<input type="checkbox"/> 9073-CL	
Wall Mount	<input type="checkbox"/> 9093W		<input type="checkbox"/> 9079W		<input type="checkbox"/> 9074W	
Cabinet Mount	<input type="checkbox"/> 9093C		<input type="checkbox"/> 9079C		<input type="checkbox"/> 9074C	
Track Mount	<input type="checkbox"/> 9083T		<input type="checkbox"/> 9078T		<input type="checkbox"/> 9073T	

STEP 13: OPERATORY LIGHT CONFIGURATION AND OPTIONS

Configure based on the light-mounting option selected above

STEP
13

Chair Mounted Light Options (Pivot, Sidebox, and Fixed)

Description	Part #	Qty
<input type="checkbox"/> 20" Rigid Arm in Lieu of Standard 32" Arm	1110-227-CS	
Light Posts (Select 1)		
<input type="checkbox"/> Sidebox, 44" Post	1110-900	
<input type="checkbox"/> Sidebox w/Monitor Cord slot, 44"	1110-170	
<input type="checkbox"/> Pivot without chair adapter (adapter included with unit purchase)	1113-654	
<input type="checkbox"/> Pivot with 3800/3900 chair adapter (when chair-mounted unit not ordered)	1113-654-CA	
Monitor Mount Assemblies		
<input type="checkbox"/> Ring Post Mount Assembly (Must pair with either PN 1110-170 or 1110-045 above) This option allows you to have both a light and monitor mounted to your chair	1300RM	
<input type="checkbox"/> Top Post Mount Assembly (Must pair with one of the post options below) This option allows a monitor instead of a light on your chair	1300TPM	
<input type="checkbox"/> Sidebox: Monitor Post, Dual Slot, Fixed Sidebox Mount, 30"	1331-998	
<input type="checkbox"/> Pivot: Monitor in Lieu of Light Post, Pivot	1331-654	
<input type="checkbox"/> Pivot: Monitor in Lieu of Light Post, Pivot with Chair Adapter	1331-654-CA	
Ceiling Mounted Light Configuration (Select Mounting Height and Post Length)		
Mounting Height (Required) <input type="checkbox"/> 8' <input type="checkbox"/> 8'6" <input type="checkbox"/> 9' <input type="checkbox"/> 9'6" <input type="checkbox"/> 10' <input type="checkbox"/> 10'6" <input type="checkbox"/> 11' <input type="checkbox"/> Custom Height, Specify:		
Post Length (Required) <input type="checkbox"/> 12" <input type="checkbox"/> 18" <input type="checkbox"/> 24" <input type="checkbox"/> 30" <input type="checkbox"/> 36" <input type="checkbox"/> 42" <input type="checkbox"/> 48" <input type="checkbox"/> Custom Length (max 48"), Specify:		
Cabinet and Wall Mounted Light Configuration and Options		
Description	Part #	Qty
Arm Height (Required): <input type="checkbox"/> Standard <input type="checkbox"/> Tall	N/A	
Specify Cabinet or Wall Mounting Height (write-in):	N/A	
<input type="checkbox"/> Horizontal Top Cabinet Mount in Lieu of Internal Mount	1110-881-CS	
<input type="checkbox"/> Vertical Exterior Cabinet Mount in Lieu of Internal Mount	1110-886-CS	
Track Mounted Light Configuration and Options		
<input type="checkbox"/> Monitor Mount Assembly (specify monitor post length in next section)	1305T	
<input type="checkbox"/> S Video Cable, Pre-Installed	0014-282	
<input type="checkbox"/> RCA Cable, Pre-Installed	0014-283	
Select Track Mount Light Post Lengths		
Mounting Height (Required) <input type="checkbox"/> 8' <input type="checkbox"/> 8'6" <input type="checkbox"/> 9' <input type="checkbox"/> 9'6" <input type="checkbox"/> 10' <input type="checkbox"/> 10'6" <input type="checkbox"/> 11' <input type="checkbox"/> Custom Height, Specify:		
Light Post Length (Required) <input type="checkbox"/> 7.25" <input type="checkbox"/> 13.75" <input type="checkbox"/> 19.75" <input type="checkbox"/> 25.75" <input type="checkbox"/> 31.75" <input type="checkbox"/> 37.75" <input type="checkbox"/> 43.75" <input type="checkbox"/> Custom Length (max 43.75"), Specify:		
Monitor Post Length (Required if 1305T selected) <input type="checkbox"/> 20.25" <input type="checkbox"/> 26.25" <input type="checkbox"/> 32.25" <input type="checkbox"/> 38.25" <input type="checkbox"/> 44.25" <input type="checkbox"/> 50.25" <input type="checkbox"/> 56.25" <input type="checkbox"/> Custom Length (max 56.25"), Specify:		

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

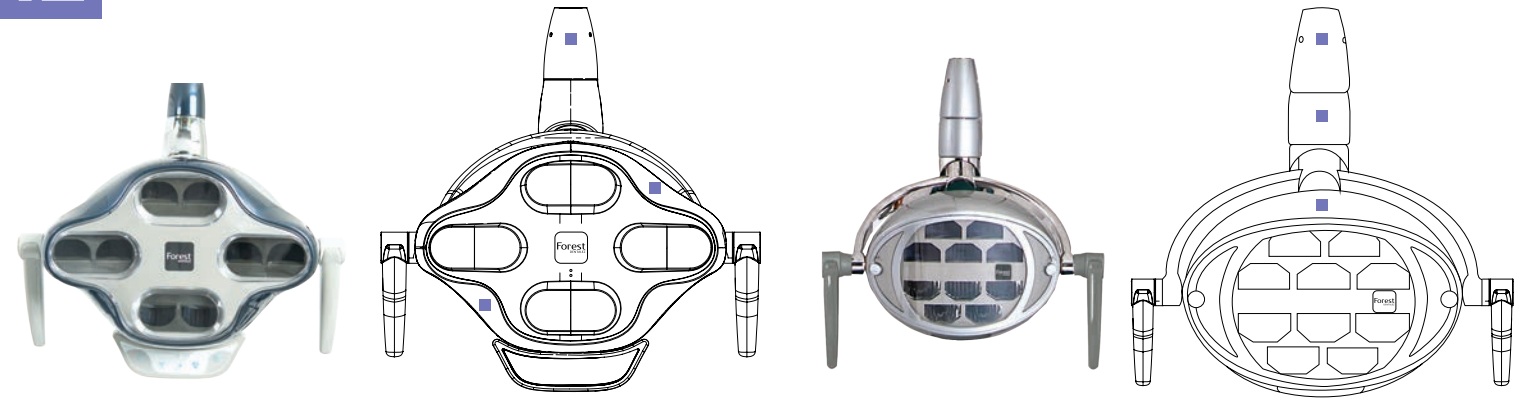
Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Powder coat paint



LIGHT IMAGES



- Pivot & Fixed Chair Mount Packages
- Pivot & Fixed Chair Mount Package Ordering Guide
- Fusion Packages
- Fusion Package Ordering Guide
- Designer Friendly Colors & Surfaces
- Chairs & Stools
- Lights and Monitor Mounts
- Rear Mount Side Delivery Systems
- Ancillary Ordering & Options
- Control Head Accessories
- Simplified Cuspidors & Assistant's Instrumentation
- Simplified Rear & Side Delivery Systems
- Factory Options & Parts
- General Information
- Warranty

THE FUSION PACKAGE

Forest

DENTALEZ



LUXURY MEETS AFFORDABILITY

The Fusion Package integrates thoughtful design with premium materials for years of beauty at an exceptional value.

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Base Package



Popular Options



Dr's Economy Stool,
Naugahyde Upholstery

Part number: 6156V



Assistants Economy Stool,
Naugahyde Upholstery

Part number: 6157V



E- Torque Electric
Handpiece System

Part number:
STAR-ETORQUE



Assistants Instrumentation
Rear Pivot with Touchpad

Part number: 1296

Base Package As Shown

Description	Part Number
3800 Patient Chair, Pretubed	3800PT
Delivery System, w/Fusion, Pivot	4296PI
Assistants Instrumentation, Rear Pivot	1295
LED Operatory Light, 9078	9078
Light Post, Pivot	1113-654
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Popular Options

Description	Part Number
Dr. Economy Stool, Naugahyde Upholstery	6156V
Assistants Economy Stool, Naugahyde Upholstery	6157V
E-Torque Electric HP System	STAR-ETORQUE
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Quick Disconnect for Water	0016-454-CS
Quick Disconnect for Air	0016-456-CS
Routing Valve Kit	1115-315-CS
Additional HVE w/Dry Oral Cup	1400-113-CS
Ultraleather Upgrade	3807-999-UL
Assistants Instrumentation, Rear Pivot with Touchpad	1296

Equipment dimensions, see page 95
Shipping weight estimate:
Delivery Unit & Chair 445 lb / 203 kg
Base package 445 lb / 203 kg
As shown package 585 lb / 266 kg




Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Base Package




Popular Options




Dr's Economy Stool,
Naugahyde Upholstery

Part number: 6156V



Assistants Economy Stool,
Naugahyde Upholstery

Part number: 6157V



E- Torque Electric
Handpiece System

Part number:
STAR-ETORQUE

Base Package As Shown

Description	Part Number
3800 Patient Chair, Pretubed	3800PT
Delivery System, Sidebox w/Fusion	4295SI
Routing Valve	Included
Quick Disconnect for Water	Included
Quick Disconnect for Air	Included
Deluxe Cuspidor	1577
LED Operatory Light, 9078	9078
Light Post, Sidebox	1110-900
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Equipment dimensions, see page 96

Shipping weight estimate:

Delivery unit & chair	440 lb / 202 kg
Base package	440 lb / 202 kg
As shown package	561 lb / 255 kg



Popular Options

Description	Part Number
Dr. Economy Stool, Naugahyde Upholstery	6156V
Assistants Economy Stool, Naugahyde Upholstery	6157V
E-Torque Electric HP System	STAR-ETORQUE
Satelec SP	SATELEC-S
Satelec SP, LED lit	SATELEC-S-LED
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Swivel Tray	1110-468
Additional Handpiece Position, Wet	1100-189-CS
Deluxe Cuspidor Remote Activation	1100-480-CS

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Fusion Packages

Fusion, 2186BC Rear Delivery

Base Package



Popular Options



Dr's Economy Stool,
Naugahyde Upholstery
Part number: 6156V



Assistants Economy Stool,
Naugahyde Upholstery
Part number: 6157V

Base Package As Shown

Description	Part Number
3800 Patient Chair	3800
Delivery System, Beneath Counter Duo	2186BC
LED Operator Light, 9078	9078
Light Post, Pivot w/Chair Adapter	1113-654-CA
Chair Control (select from unit-mounted, chair-mounted, or foot control)	3914-998-CS

Popular Options

Description	Part Number
Dr. Economy Stool, Naugahyde Upholstery	6156V
Assistants Economy Stool, Naugahyde Upholstery	6157V
6-PIN Handpiece Illumination w/tubing	1700-019-CS
Quick Disconnect for Water	0016-454-CS
Quick Disconnect for Air	0016-456-CS
Routing Valve Kit	1115-315-CS
Additional HVE w/Dry Oral Cup	1400-113-CS
Ultraleather Upgrade	3807-999-UL
LED Operator Light, Ceiling Mounted	9078-CL
LED Operator Light, Wall Mounted	9079W
Deluxe Utility Center Housing and Cover	1115-007-G

Shipping weight estimate:
Delivery Unit & Chair 445 lb / 203 kg
Base package 445 lb / 203 kg
As shown package 585 lb / 266 kg



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

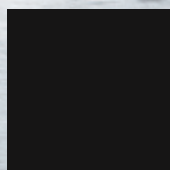
Forest

DENTALEZ

FUSION PACKAGE ORDERING GUIDE




NAUGAHYDE® UPHOLSTERY OPTIONS



Black



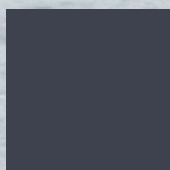Cocoa



Deep
Sapphire



Saddle



Sand



Taupe

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

DESIGN YOUR FUSION OPERATORY CHAIR

STEP 1

STEP 1: SELECT YOUR CHAIR - Select 1 chair

Description	Part #	Qty.	Comments, color
<input type="checkbox"/> 3800 Patient Chair, non-tubed	3800		Forest Grey w/Shadow baseplate only
<input type="checkbox"/> 3800 Patient Chair, pre-tubed	3800PT		Forest Grey w/Shadow baseplate only

STEP 2

STEP 2: SELECT YOUR UPHOLSTERY - Select 1 option and 1 upholstery color

<input type="checkbox"/> Standard Naugahyde® Upholstery, Seamless (Included)	N/A		<input type="checkbox"/> Black <input type="checkbox"/> Deep Sapphire <input type="checkbox"/> Sand <input type="checkbox"/> Cocoa <input type="checkbox"/> Saddle <input type="checkbox"/> Taupe
<input type="checkbox"/> Upgrade to standard Ultraleather® color, Seamless	3807-999-UL		Standard Ultraleather Colors <input type="checkbox"/> Autumn Leaf <input type="checkbox"/> Boysenberry <input type="checkbox"/> Diplomat Blue <input type="checkbox"/> Dove Grey <input type="checkbox"/> Parrot <input type="checkbox"/> St. John <input type="checkbox"/> Baltic <input type="checkbox"/> Charcoal <input type="checkbox"/> Doe <input type="checkbox"/> Papyrus <input type="checkbox"/> Raven Wing <input type="checkbox"/> Stone
<input type="checkbox"/> Upgrade, special color Ultraleather, Seamless	3807-999-UL-SP		Standard Pearlized Ultraleather Colors <input type="checkbox"/> Midnight <input type="checkbox"/> Oz <input type="checkbox"/> Tin Man <input type="checkbox"/> Pewter <input type="checkbox"/> Pharaoh
<input type="checkbox"/> Upgrade, special color Ultraleather, Seamless	3807-999-UL-SP		Specify Color:

STEP 3

STEP 3: ADDITIONAL UPHOLSTERY OPTIONS - Select at least 1 of the below options

Description	Part #	Qty.	Comments, color
<input type="checkbox"/> No Additional Options	N/A		
<input type="checkbox"/> Neck Pillow (only available in Ultraleather)	3907-318		Color:
<input type="checkbox"/> Knee Pillow (only available in Ultraleather)	3907-317		Color (inner if dual color):
<input type="checkbox"/> Dual Color Upholstery, Knee Pillow	3907-317-DC		Outer Color:
<input type="checkbox"/> Contrast Stitch Upholstery, Knee Pillow	3907-317-DS		Thread Color:
<input type="checkbox"/> Pedo Booster Seat (only available in Ultraleather)	3907-180		Color:

STEP 4

STEP 4: UTILITY CENTER OPTIONS

Description	Part #	Qty.	Comments, color
<input type="checkbox"/> No Utility Center*	N/A		Select this option if you are buying a chair-mounted delivery unit
<input type="checkbox"/> Deluxe Utility Center Housing & Cover, Rear Umbilical Access, Forest Grey	1115-007-G		
<input type="checkbox"/> Deluxe Utility Center Housing & Cover, Right Umbilical Access, Forest Grey	1115-009-G		

*If you are ordering a chair-mounted delivery unit then select "No Utility Center" because one will be included with your delivery unit

STEP 5

STEP 5: CHAIR CONTROL LOCATION - Select at least 1 chair control

Maximum of 1 control can be mounted to the delivery unit, maximum of 2 on the chair (left and right), 4 total maximum number of controls per operatory			
<input type="checkbox"/> First Chair Control (Required) First control included with chair purchase	3914-998-CS	<input type="checkbox"/> Delivery Unit <input type="checkbox"/> Chair, right side	<input type="checkbox"/> Footswitch <input type="checkbox"/> Chair, left side
<input type="checkbox"/> Second Chair Control (Optional)	3914-999-CS	<input type="checkbox"/> Delivery Unit <input type="checkbox"/> Chair, right side	<input type="checkbox"/> Footswitch <input type="checkbox"/> Chair, left side

Pivot & Fixed Chair Mount Packages

Pivot & Fixed Chair Mount Package Ordering Guide

Fusion Packages

Fusion Package Ordering Guide

Designer Friendly Colors & Surfaces

Chairs & Stools

Lights and Monitor Mounts

Rear Mount Side Delivery Systems

Ancillary Ordering & Options

Control Head Accessories

Simplified Cuspidors & Assistant's Instrumentation

Simplified Rear & Side Delivery Systems

Factory Options & Parts

General Information

Warranty

DESIGN YOUR FUSION OPERATORY DELIVERY UNIT

STEP 6

STEP 6: SELECT YOUR DELIVERY SYSTEM - Select 1 delivery system

Description	Part #	Qty.	Price	Comments
<input type="checkbox"/> Fusion Delivery Head, Pivot Chair Mount	4296PI			Forest Grey paint & plastics only
<input type="checkbox"/> Fusion Delivery Head, Fixed Chair Mount with Sidebox	4295SI			Forest Grey paint & plastics only

STEP 7

STEP 7: CONFIGURE YOUR DELIVERY SYSTEM

Note: Please refer to ancillary and options graph in the catalog to determine number of ancillaries and options possible

Description	Part #	Qty.	0	1	2	3	4	5	6
<input type="checkbox"/> Air/Water Syringe (included standard)									
<input type="checkbox"/> 4-hole handpiece tubing (included standard)									
<input type="checkbox"/> 5-hole system (1 tubing included, system accomodates a 2nd tubing)	1700-021-CS								
<input type="checkbox"/> - Second 5-hole tubing	0017-615-CS								
<input type="checkbox"/> 6-pin system (1 tubing included, system accomodates a 2nd tubing)	1700-019-CS								
<input type="checkbox"/> - Second 6-pin tubing	0017-603-CS								
<input type="checkbox"/> Star® ETorque Electric System, Factory Installed*	STAR-ETORQUE								
<input type="checkbox"/> Star Blis-sonic™ SW Scaler (with 4 line swivel)*	STAR-BLIS-SW								
<input type="checkbox"/> Star Blis-sonic K Scaler (with fixed back-end)*	STAR-BLIS-K								
<input type="checkbox"/> Satelec SP Scaler*	SATELEC-S								
<input type="checkbox"/> Satelec SP LED Scaler*	SATELEC-S-LED								
<input type="checkbox"/> Cavitron Scaler*	CAVITRON								
<input type="checkbox"/> Additional Handpiece Position, Wet	1100-189-CS								
<input type="checkbox"/> Additional Handpiece Position, Air Drive Only	1100-190-CS								
<input type="checkbox"/> Assistant's hygiene instrumentation, control head-mounted Select 1 position for the SE and 1 position for the HVE	1408-CS		SE						
			HVE						
<input type="checkbox"/> Star Handpiece Bundle (Write-in PN and price)									
<input type="checkbox"/> Quick Disconnect w/Flow Control, 1/4" (for water)	0016-454-CS		Included with 4295SI						
<input type="checkbox"/> Quick Disconnect w/Flow Control, 3/8" (for air)	0016-456-CS		Included with 4295SI						
<input type="checkbox"/> Swivel Tray	1110-468								
<input type="checkbox"/> Less sidebox telescoping assistant's instrumentation	1400-000-CS		4295SI only						
<input type="checkbox"/> Chrome Accents	9000-998-CS								
<input type="checkbox"/> Additional Umbilical, Per Foot (max of 9 feet)	1115-999-CS								
<input type="checkbox"/> Routing Valve Kit w/Water shutoff & Regulator, Factory inst.	1115-315-CS		Included with 4295SI						
<input type="checkbox"/> Routing Valve Kit for Cuspidor, Factory Installed	1115-316-CS								

*Ancillaries utilize an existing handpiece position

STEP 8

STEP 8: ASSISTANT'S INSTRUMENTATION AND CUSPIDORS - Must select 1 option for order to be valid

<input type="checkbox"/> No Assistant's Instrumentation or Cuspidor*	N/A		
<input type="checkbox"/> Rear Pivot Instrumentation, Quad Holder Bar	1295		Not available with 4295SI
<input type="checkbox"/> Rear Pivot Instrumentation with Integrated Touchpad	1296		Not available with 4295SI
<input type="checkbox"/> Rear Pivot Instrumentation with Deluxe Cuspidor	1592		Not available with 4295SI
<input type="checkbox"/> Deluxe Cuspidor only, Rear Pivot	1582		Not available with 4295SI
<input type="checkbox"/> Basic Cuspidor, Sidebox Mounted	1679		Only available with 4295SI
<input type="checkbox"/> Deluxe Cuspidor, Sidebox Mounted	1577		Only available with 4295SI

*Sidebox mounted telescoping assistant's instrumentation included standard with 4295SI

STEP 9

STEP 9: ASSISTANT'S INSTRUMENTATION AND CUSPIDOR OPTIONS

<input type="checkbox"/> Quick Disconnect w/Flow Control, 1/4" (For 1095 instrumentation)	0016-455-CS	Specify One	<input type="checkbox"/> Air	<input type="checkbox"/> Water
<input type="checkbox"/> Additional Autoclavable HVE with 5' Tubing & Dry Oral Cup	1400-113-CS			
<input type="checkbox"/> Hands-Free HVE Holder with HVE and 5' Tubing	1400-279-STG			
<input type="checkbox"/> Self-Contained Clean Water System w/Quick Disconnect, 1.5L bottle (For 1095 Instrumentation only)	1700-319			
<input type="checkbox"/> Tray Assembly (For 1095 Instrumentation only)	1110-680			
<input type="checkbox"/> Deluxe Cuspidor Unit-Mounted Activation (For 1577, 1582, 1592)	1100-480-CS			
<input type="checkbox"/> Chair Control Touchpad - Mounted to 1095 Instrumentation	1700-088-CS			

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

DESIGN YOUR OPERATORY LIGHT SELECTION

STEP 10: OPERATORY LIGHT

STEP
10

Select Light Head Mounting Option - 1 light and mounting option required

	9083 Halogen		9078 LED	
	Part #	Qty	Part #	Qty
Lights with Fusion packages come in Forest Grey Only				
Chair Mount	<input type="checkbox"/> 9083		<input type="checkbox"/> 9078	
Ceiling Mount	<input type="checkbox"/> 9083-CL		<input type="checkbox"/> 9078-CL	
Wall Mount	<input type="checkbox"/> 9093W		<input type="checkbox"/> 9079W	
Cabinet Mount	<input type="checkbox"/> 9093C		<input type="checkbox"/> 9079C	
Track Mount	<input type="checkbox"/> 9083T		<input type="checkbox"/> 9078T	

STEP 11: OPERATORY LIGHT CONFIGURATION AND OPTIONS

STEP
11

Chair Mounted Light Options (Pivot, Sidebox, and Fixed)

Description

Part #

Qty

☐ 20" Rigid Arm in Lieu of Standard 32" Arm

1110-227-CS

Light Posts (Select 1)

☐ Sidebox, 44" Post

1110-900

☐ Sidebox w/Monitor Cord slot, 44"

1110-170

☐ Pivot without chair adapter (adapter included with unit purchase)

1113-654

☐ Pivot with 3800/3900 chair adapter (when chair-mounted unit not ordered)

1113-654-CA

Monitor Mount Assemblies

☐ Ring Post Mount Assembly (Must pair with either PN 1110-170 or 1110-045 above)
This option allows you to have both a light and monitor mounted to your chair

1300RM

☐ Top Post Mount Assembly (Must pair with one of the post options below)
This option allows a monitor instead of a light on your chair

1300TPM

☐ Sidebox: Monitor Post, Dual Slot, Fixed Sidebox Mount, 30"

1331-998

☐ Pivot: Monitor in Lieu of Light Post, Pivot

1331-654

☐ Pivot: Monitor in Lieu of Light Post, Pivot with Chair Adapter

1331-654-CA

Ceiling Mounted Light Configuration (Select Mounting Height and Post Length)

Mounting Height (Required)

☐ 8' ☐ 8'6" ☐ 9' ☐ 9'6" ☐ 10' ☐ 10'6" ☐ 11' ☐ Custom Height, Specify:

Post Length (Required)

☐ 12" ☐ 18" ☐ 24" ☐ 30" ☐ 36" ☐ 42" ☐ 48" ☐ Custom Length (max 48"), Specify:

Cabinet and Wall Mounted Light Configuration and Options

Description

Part #

Qty

Arm Height (Required): ☐ Standard ☐ Tall

N/A

Specify Cabinet or Wall Mounting Height (write-in):

N/A

☐ Horizontal Top Cabinet Mount in Lieu of Internal Mount

1110-881-CS

☐ Vertical Exterior Cabinet Mount in Lieu of Internal Mount

1110-886-CS

Track Mounted Light Configuration and Options

☐ Monitor Mount Assembly (specify monitor post length in next section)

1305T

☐ S Video Cable, Pre-Installed

0014-282

☐ RCA Cable, Pre-Installed

0014-283

Select Track Mount Light Post Lengths

Mounting Height (Required)

☐ 8' ☐ 8'6" ☐ 9' ☐ 9'6" ☐ 10' ☐ 10'6" ☐ 11' ☐ Custom Height, Specify:

Light Post Length (Required)

☐ 8.25" ☐ 14.25" ☐ 20.25" ☐ 26.25" ☐ 32.25" ☐ 38.25" ☐ 44.25" ☐ Custom Length (max 44.25"), Specify:

Monitor Post Length
(Required if 1305T selected)

☐ 20.25" ☐ 26.25" ☐ 32.25" ☐ 38.25" ☐ 44.25" ☐ 50.25" ☐ 56.25" ☐ Custom Length (max 56.25"), Specify:

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Surface Treatments

Forest is the only dental equipment with a WOW factor. It really is a kick watching a dentist react to Forest equipment when they first see it. A common reaction is “wow, that is beautiful.”

Carbon Fiber

- Control Heads
- Standard IC - Chassis
 - Euro - Lid

- Sidebox - Cover
- Basic Cuspidor - Cover
- Assistant’s Instrumentation
- 1095 - Arm Cover

Powder Coat Paint

Dentists and interior designers love being able to put something other than white into their beautiful patient friendly decors.

- Control Heads
- Standard IC - Lid, Brake Handles and Holder Bar
 - Euro - Chassis, Whip Arms, Brake Handles and Swivel Tray
 - 7000 Series Rear Mount, Side Delivery and Mobile Cart
 - 7020 PRO Series Rear Mount, Side Delivery and Mobile Cart

- Sidebox - Housing
- Basic Cuspidor - Housing and Holder Bar
- Assistant’s Instrumentation
- 1095 - Arm
 - 1096 - Arm and Integrated Holder

Ultraleather® (standard)

Once you see a Forest chair, delivery system and light package with color coordinated paint and upholstery in a well designed treatment room, you wonder why it’s taken so long for a manufacturer to offer dental equipment with a full range of surface treatments.

- 3900 Chair
- Plush Upholstery
 - Seamless Upholstery
 - Knee Pillow, Neck Pillow and Pedo Booster Seat
 - Deluxe and Economy Stools

Solid Surface (standard)

- Control Heads
- 7010 Series Rear Mount Side Delivery and Mobile Cart
 - 7020 PRO Series Rear Mount Side Delivery and Mobile Cart
 - 3486 Rear Mount Side Delivery and Mobile Cart
 - 2184BC/2186BC Rear and Side Delivery

Laminate (optional)

- Control Heads
- 3486 Rear Mount Side Delivery
 - 2184BC/2186BC Rear and Side Delivery
 - 7010 Series Rear Mount Side Delivery and Mobile Cart

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant’s Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Pivot & Fixed Chair Mount Packages

Pivot & Fixed Chair Mount Package Ordering Guide

Fusion Packages

Fusion Package Ordering Guide

Designer Friendly Colors & Surfaces

Chairs & Stools

Lights and Monitor Mounts

Rear Mount Side Delivery Systems

Ancillary Ordering & Options

Control Head Accessories

Simplified Cuspidors & Assistant's Instrumentation

Simplified Rear & Side Delivery Systems

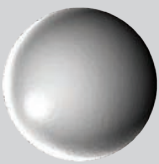
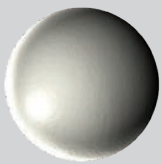
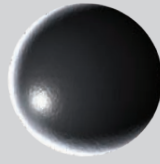
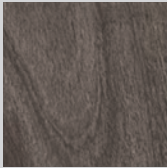
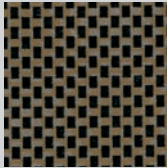
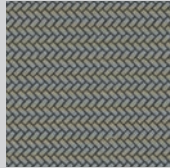
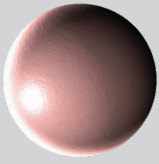
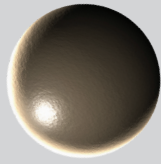
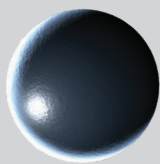
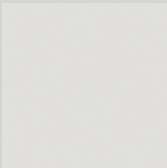
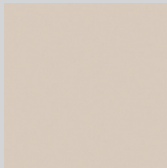
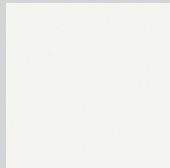
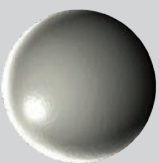
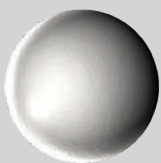
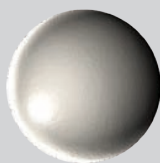

Factory Options & Parts

General Information







Warranty

POWDER COAT PAINT (9)

CARBON FIBER (3) FORMICA® (4)


 <p>ARCTIC SILVER Powdercoat Paint</p>	 <p>FOREST GREY Powdercoat Paint</p>	 <p>CHARCOAL Powdercoat Paint</p>	 <p>GREYWOOD Carbon Fiber</p>	 <p>LASER LEATHER Carbon Fiber</p>	 <p>LEMANS Carbon Fiber</p>
 <p>COPPER Powdercoat Paint</p>	 <p>BRONZE Powdercoat Paint</p>	 <p>MIDNIGHT Powdercoat Paint</p>	 <p>FOLKSTONE Formica®</p>	 <p>OYSTER GREY Formica®</p>	 <p>SURF Formica®</p>
 <p>PEWTER Powdercoat Paint</p>	 <p>SHADOW Powdercoat Paint</p>	 <p>WHITE Powdercoat Paint</p>	 <p>WHITE BLANC Formica®</p>		

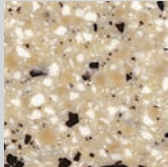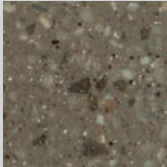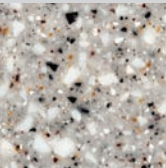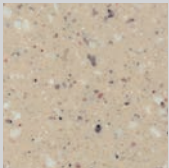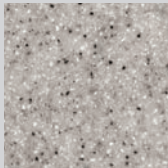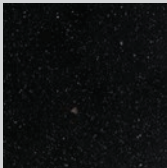SOLID SURFACES (10)

 <p>ANTARCTICA Corian®</p>	 <p>ARCTIC MELANGE Wilsonart®</p>	 <p>BAJA MELANGE Wilsonart®</p>	 <p>CINDER GRANITE Meganite®</p>	 <p>GOLDRUSH GEMSTONE Meganite®</p>	 <p>MOONSTONE Meganite®</p>
 <p>PLATINUM GRANITE Meganite®</p>	 <p>ROCKY ROAD GRANITE Meganite®</p>	 <p>SILVER MIST Meganite®</p>	 <p>STARRY STARRY NIGHT Meganite®</p>		

NAUGAHYDE®* (6)

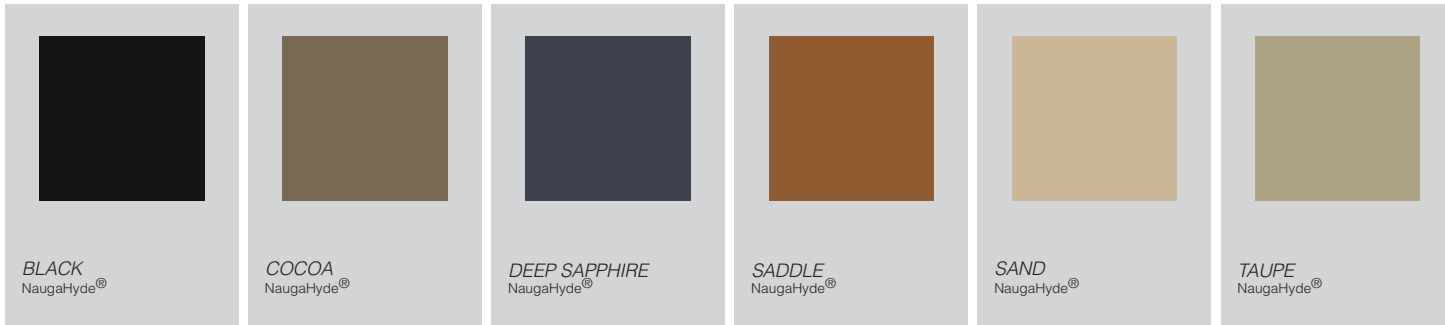*For The Fusion Packages only. Pages 22-29



ULTRALEATHER®* (12 CLASSIC) (5 PEARLIZED)

*Pearlized



CUSTOMIZE YOUR TREATMENT ROOMS WITH DENTALEZ

Powder coat paint your drawer and door pulls, grab handles, trash receptacles, or other metal surfaces in your treatment room to match your delivery system paint color

Source custom powder coat paint colors to match your provided sample

Special order any UltraLeather® upholstery color shown on the Ultrafabrics® website (from original or pearlized palette)

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Forest

CHAIRS AND STOOLS



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty



3800 CHAIR

- Material: hydraulic cantilever mechanism
- Benefit: smooth and steady heavy lifting
- Material: cast aluminum chair back, base, armrest and headrest
- Benefit: strong, impervious to fading or cracking
- Material: flip down armrests
- Benefit: easy patient entry and exit
- Material: optional soft-touch back massage with integrated compressor
- Benefit: saves central compressor capacity for delivery systems

3900 CHAIR

- Material: cast aluminum chair back, chair base and utility center base
- Benefit: strong, impervious to fading or cracking
- Material: pressure-formed plastic chair and utility center covers
- Benefit: warp, impact and stain resistant
- Material: aluminum hydraulic manifold and speed adjustment screws
- Benefit: independent speed adjustment of base and back both up and down
- Material: optional soft-touch back massage with integrated compressor
- Benefit: saves central compressor capacity for delivery systems
- Material: Ultraleather® upholstery standard
- Benefit: withstands heavy use and cleaning, with glove-soft leather feel
- Material: Forest Grey or Shadow powder coat standard on 3900 chair
- Benefit: complements treatment room interior design color palette

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

3900 Hydraulic Patient Chair




Plush

Model 3900
Part Number 3900
Retail Price \$10,463



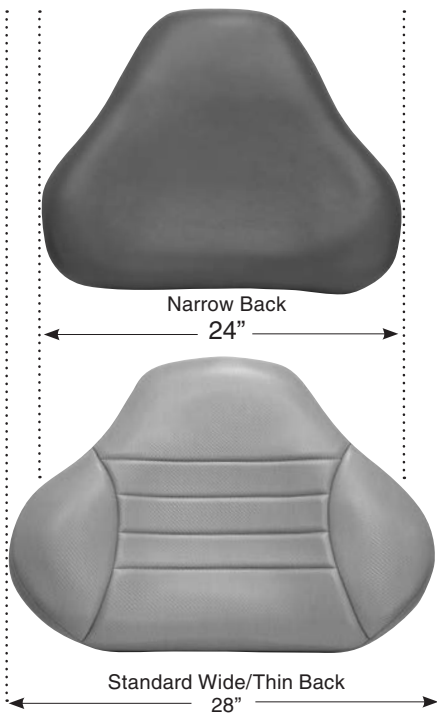Seamless



Standard Features	Factory Options	Part No
Hydraulic Chair 3900 Ultraleather® upholstery with wide/thin back for elbow support Plastic toe guard 21" Double articulating pinch-style headrest Thin back w/full lumbar support Right/left patient access Swing-down locking armrests Hydraulic cantilever mechanism Emergency shut-off 30 degree base rotation w/lock Multifunction touchpad w/4 programmable preset positions Seat height: minimum 17" maximum 29"	CHAIR (Forest Grey or Shadow) CHAIR Designer Friendly® CHAIR PRO Designer Friendly Includes: Plush dual color and deviated/ contrast stitching upholstery, massage, contoured utility center housing and cover, knee pillow and neck pillow CHAIR PRO Designer Friendly less utility center housing and cover One chair control included (specify type below) ○ Footswitch control ○ Touchpad, chair-mounted ! Additional chair control (limit 2 on chair, limit 1 on unit) ○ Footswitch control ○ Touchpad, chair-mounted Chair adapter, fixed (for external umbilical) Contoured utility center housing & cover Contoured utility center housing & cover Designer Friendly ! UPHOLSTERY SELECTIONS: Pedo booster seat (see page 40) Knee pillow (see page 40) Neck pillow (see page 40) Soft-touch back massage option (not available with Memory Comfort) Plastic toe guard 29" in lieu of 21" ! LIGHT POSTS (for unit mount lights/post mount monitor assemblies see pages 46) Fixed (chair adaptor sold separately) 44" light post without slot, chair mount 44" light post with slot for monitor cord Pivot (includes pivot cast chair adaptor, bushing, bearing kit - mounts on Forest chairs only) Pivot monitor post, (includes pivot cast chair adaptor, bushing, bearing kit - mounts on Forest chairs only) (monitor mount assembly not included) Monitor mount assembly, top post mount (see page 49) Monitor mount assembly, ring mount (see page 49) ! Often overlooked Δ Lead time varies	3900 3900C 3900C-PRO 3900C-PRO-LUC 3914-998-CS* 3914-999-CS 1115-195 3900-101 3900-101-C 3907-180 3907-317 3907-318 3900-100 3907-028-CS 1110-794 1110-050 1113-654-CA 1331-654-CA 1300TPM 1300RM

! Designer Friendly® Surface Treatments

STEP 1 SELECT STANDARD POWDER COAT PAINT (CHAIR):

☐ Forest Grey

☐ Shadow

UPGRADE: Chair Base/Back Powder Coat Paint
UPGRADE: PRO Chair
(Less Utility Center Housing and Cover)

☐ Arctic Silver

☐ Bronze

☐ Charcoal

☐ Copper

☐ Midnight

☐ Pewter

☐ White

STEP 2 SELECT ULTRALEATHER® UPHOLSTERY COLOR CHOICES FOR CHAIR AND DELUXE /ECONOMY STOOLS:

☐ Autumn Leaf☐ Baltic☐ Charcoal

☐ Diplomat Blue☐ Doe☐ Dove Grey

☐ Boysenberry☐ Papyrus☐ Parrot

☐ Raven Wing☐ Stone☐ St. John

Pearlized Colors: ☐ Midnight ☐ Pewter ☐ Oz
☐ Pharaoh ☐ Tin Man

Δ OR SELECT Special Color Ultraleather (Call for leadtime)



☐ Chair

3907-999-SP

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

37

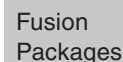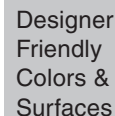TABLE OF CONTENTS

Fusion
Package
Ordering GuideChairs
& Stools

Lights and Monitor Mounts

Rear Mount Side Delivery Systems

Ancillary Ordering & Options

Control Head
Accessories

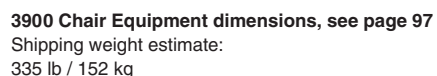
Simplified Cuspidors & Assistant's Instrumentation

Simplified Rear & Side Delivery Systems

Factory
Options & Parts

General Information

Warranty



3800 Hydraulic Patient Chair

Model 3800

Part Number 3800
Retail Price \$7,613



Standard Features

Hydraulic Chair 3800

Naugahyde® upholstery
w/wide back for elbow support
Plastic toe guard 21"
Double articulating pinch-style headrest
Right/left patient access
Swing-down locking armrests
Hydraulic cantilever mechanism
Aluminum casting chair back, base,
armrest and headrest
Emergency shut-off
30-degree swivel rotation w/lock

Factory Options

Part No

CHAIR (pre-tubed) 3800PT
CHAIR (non-tubed) 3800
One chair control included (specify type below) 3914-998-CS
☐ Footswitch control
☐ Touchpad, chair-mounted
! Additional chair control (limit 2 on chair, limit 1 on unit) 3914-999-CS
☐ Footswitch control
☐ Touchpad, chair-mounted
Chair adapter, fixed (for external umbilical) 1115-195
Deluxe utility center housing and cover, rear umbilical access 1115-007-G
Deluxe utility center housing and cover, right umbilical access 1115-009-G
! **UPHOLSTERY SELECTIONS (only available in ultraleather):**
Pedo booster seat (see page 40) 3907-180
Knee pillow (see page 40) 3907-317
Neck pillow (see page 40) 3907-318
! **LIGHT POSTS** (for unit mount lights/post mount monitor assemblies see pages 46)
Fixed (chair adaptor sold separately) 1110-794
44" light post without slot, chair mount 1110-050
44" light post with slot for monitor cord 1110-050
Pivot (includes pivot cast chair adaptor, bushing, bearing kit - mounts on Forest chairs only) 1113-654-CA
Pivot monitor post, (includes pivot cast chair adaptor, bushing, bearing kit - mounts on Forest chairs only) (monitor mount assembly not included) 1331-654-CA
Monitor mount assembly, top post mount (see page 49) 1300TPM
Monitor mount assembly, ring mount (see page 49) 1300RM

! Designer Friendly® Surface Treatments

STEP 1 SELECT STANDARD POWDER COAT PAINT (CHAIR):

☐ Forest Grey

STEP 2 SELECT NAUGAHYDE® UPHOLSTERY COLOR FOR YOUR CHAIR:

☐ Black ☐ Cocoa ☐ Deep Sapphire
☐ Saddle ☐ Sand ☐ Taupe

UPGRADE: Standard Ultraleather® upholstery (select from Forest standard options) 3807-999-UL

Δ **OR SELECT** Special Color upholstery (Call for leadtime)
Ultraleather 3807-999-UL-SP
Naugahyde 3807-999-SP

! Often overlooked
Δ Lead time varies



3800 Chair Equipment dimensions, see page 96
305 lb / 152 kg

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Chair Accessories

Neck Pillow



3907-318

Provides ergonomic support for patients with neck discomfort and/or during long procedures.

Neck support offers patients a comfortable relaxed experience.

Specify Ultraleather® upholstery color

Knee Pillow



3907-317

Knee support ideal for patients with joint discomfort and/or long procedures.

Zero gravity posture, eliminates pressure on knees.

Specify Ultraleather upholstery color

3900 dual color upholstery option (upcharge) 3907-317-DC

3900 contrast/deviated stitching option (upcharge) 3907-317-DS

Pedo Booster Seat (to accommodate smaller patients)



3907-180

A cushion that takes the lumbar out of the chair, essentially turning it into a bed, ideal for children.

Specify Ultraleather upholstery color

3900/3800 Chair Control, Footswitch

Additional Chair Control, Footswitch



Factory installed 3914-998-CS

Factory installed 3914-999-CS

After Market 3914-010-SH

4 programmable preset positions

3900/3800 Remote Mount Touchpad with 25' Cable



FIXED HORIZONTAL MOUNT

White 1700-079

Forest Grey 1700-079-G

Designer Friendly® (specify color) 1700-079-C

SWING HORIZONTAL MOUNT

White 1700-089-HM

Forest Grey 1700-089-HM-G

Designer Friendly (specify color) 1700-089-HM-C

SWING VERTICAL MOUNT

White 1700-089-VM

Forest Grey 1700-089-VM-G

Designer Friendly (specify color) 1700-089-VM-C

• 4 programmable preset positions

3900 Contoured Utility Center Housing & Cover



Factory option: 3900-101

Forest Grey 3900-101-G

Shadow 3900-101-SH

Designer Friendly (specify color) 3900-101-C

Included w/4486PI, 4485SI, 4987SI and 4988PI chair mounted packages

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Economy and Deluxe Stools

Nomenclature Key Example

C = Ultraleather	6157C
V = Naugahyde	6156V



Model 6156	Model 6157	Model 6159	Model 6160	Model 6161
Part Number 6156C	Part Number 6157C	Part Number 6159C	Part Number 6160C	Part Number 6161C
Part Number 6156V	Part Number 6157V	Part Number 6159V	Part Number 6160V	Part Number 6161V

Standard Features

Economy and Deluxe Stools

Ultraleather® or Naugahyde® Upholstery
Injection molded foam
5-leg fully polished base
Twin wheel hard floor casters
Accommodates maximum of 300 lbs.

Doctor's Economy 6156C/6156V

Backrest, non-adjusting
Round seat with single lever height adjustment
Pneumatic cylinder height adjustment, 18" - 23"

Assistant's Economy 6157C/6157V

Round seat with single lever height adjustment
Torso support, non-ratcheting, with height adjustment
Pneumatic cylinder height adjustment, 21 1/2" - 28"
Foot ring affixed to base

Doctor's Deluxe 6159C/6159V

Backrest with height adjustment
Contoured seat with triple lever adjustment
Pneumatic cylinder height adjustment, 18 1/2" - 23"

Assistant's Deluxe 6160C/6160V

Contoured seat with single lever adjustment
Torso support with ratcheting and height adjustment
Pneumatic cylinder height adjustment, 21 1/2" - 28"
Foot ring with adjustment

Assistant's Deluxe with Back 6161C/6161V

Backrest with height adjustment
Contoured seat with triple lever adjustment
Torso support with ratcheting and height adjustment
Pneumatic cylinder height adjustment, 21 1/2" - 28"
Foot ring with adjustment

Equipment dimensions, see page 97

Shipping weight estimate:	
6156/6159	35 lb / 15 kg
6161	45 lb / 21 kg
6157/6160	40 lb / 18 kg

Factory Options

STOOL	MODEL CHOICE	Part No
Δ Special color upholstery (upcharge)		6007-999-SP
! CYLINDER CHOICES:		
Doctor's short cylinder, 16" – 21" (6156C/6159C)		6000-628-CS
Doctor's tall cylinder, 21 1/2" – 28" (6156C/6159C)		6000-604-CS
foot ring with support highly recommended		
Assistant's short cylinder, 18" – 23" (6157C/6160C/6161C)		6000-605-CS
Foot ring with support (6156C/6159C)		6000-602
! Often overlooked		
* Specify vinyl manufacturer and color (not all colors may be available.) (Call Forest for ordering details 1-800-423-3555.)		
Δ Lead time varies		

!

Designer Friendly® Surface Treatments

STEP 1 SELECT ULTRALEATHER® UPHOLSTERY COLOR CHOICES FOR CHAIR AND DELUXE /ECONOMY STOOLS:

<input type="radio"/> Autumn Leaf	<input type="radio"/> Baltic	<input type="radio"/> Charcoal	<input type="radio"/> Diplomat Blue	<input type="radio"/> Doe	<input type="radio"/> Dove Grey
<input type="radio"/> Boysenberry	<input type="radio"/> Papyrus	<input type="radio"/> Parrot	<input type="radio"/> Raven Wing	<input type="radio"/> Stone	<input type="radio"/> St. John
Pearlized Colors:	<input type="radio"/> Midnight	<input type="radio"/> Pewter	<input type="radio"/> Pharaoh	<input type="radio"/> Oz	<input type="radio"/> Tin Man

OR SELECT NAUGAHYDE® UPHOLSTERY COLOR:

<input type="radio"/> Black	<input type="radio"/> Cocoa	<input type="radio"/> Deep Sapphire	<input type="radio"/> Saddle	<input type="radio"/> Sand	<input type="radio"/> Taupe
-----------------------------	-----------------------------	-------------------------------------	------------------------------	----------------------------	-----------------------------

OR SELECT Special Color Ultraleather® (Call for leadtime)

<input type="radio"/> Deluxe/Economy Stool	6007-999-SP
--	-------------

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cupidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Comfort Stools



Model 6176

Part Number 6176



Model 6187

Part Number 6187

Standard Features

Comfort Stools

- Air circulating weave construction
- Active tilt® seat (actively moves with user)
- Tilt seat lever function
- Contoured seat
- 5-leg base
- Twin wheel hard floor casters
- Accommodates maximum of 300 lbs.

Doctor’s Comfort 6176

- Triple lever adjustment - independent back and seat
- Backrest with angle & height adjustment
- Pneumatic cylinder height adjustment, 19 1/4”– 24 1/2”

Assistant’s Comfort 6187

- Triple lever adjustment - independent back and seat
- Backrest with angle & height adjustment
- Fully adjustable ratcheting assistant arm
- Foot ring with height adjustment
- Pneumatic cylinder height adjustment, 23”– 30 1/2”

Factory Options

STOOL

! Cylinder choices:

Doctor’s/Office short cylinder, 16 ” – 21” (6176)

Doctor’s/Office tall cylinder, 21 1/2” – 28” (6176)
foot ring with support highly recommended

Assistant’s short cylinder, 18” – 23” (6187)

Foot ring with support (6176)

MODEL CHOICE

Part No

6000-628-CS

6000-604-CS

6000-605-CS

6000-602

! Often overlooked

Equipment dimensions, see page 97

Shipping weight estimate:

6176 48 lb / 23 kg

6187 60 lb / 28 kg

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant’s
Instrumentation

Simplified
Rear & Side
Delivery
Systems

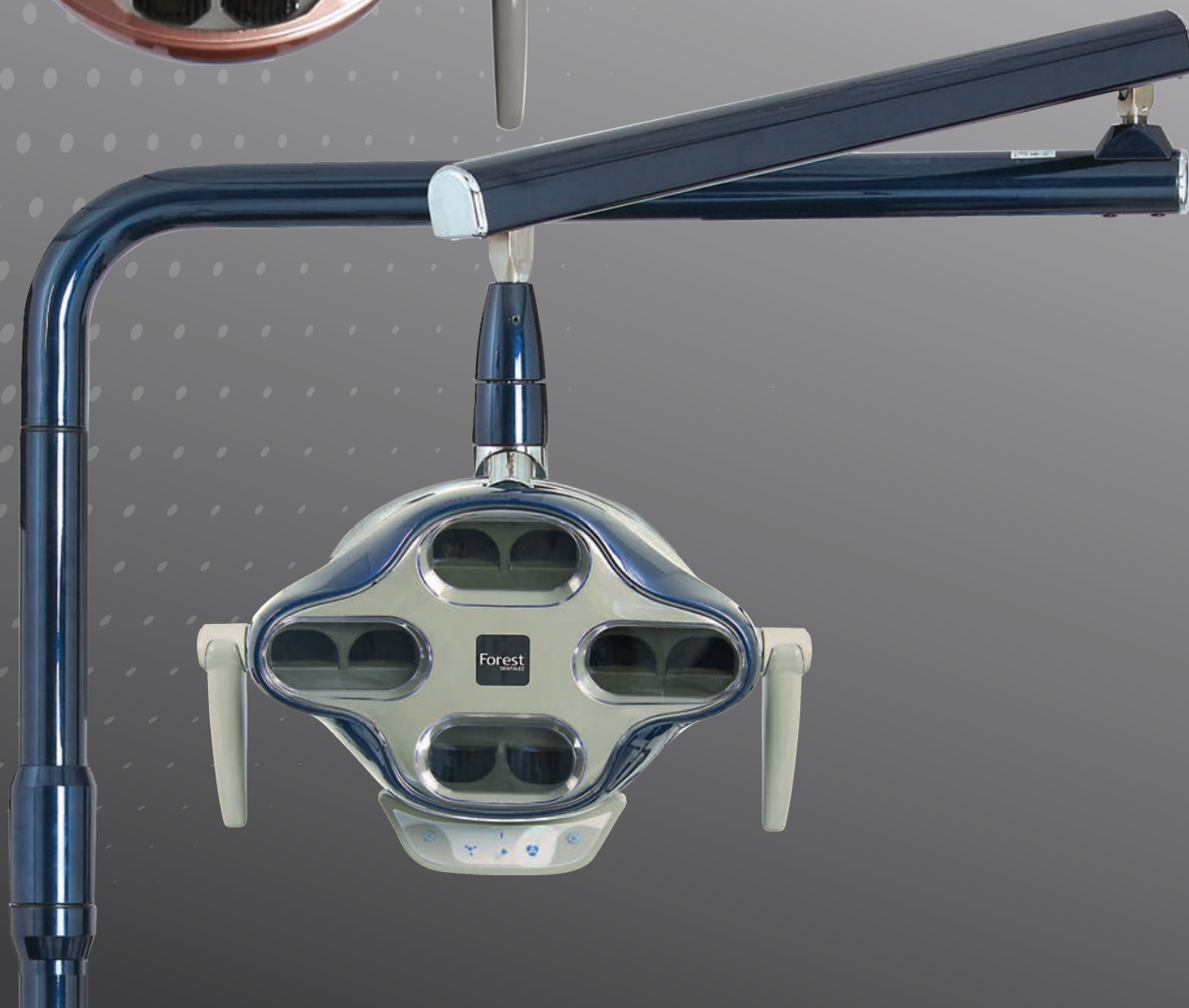
Factory
Options & Parts

General
Information

Warranty

Forest

LIGHTS AND MONITOR MOUNTS



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty



COPPER
Powdercoat Paint

ARCTIC SILVER
Powdercoat Paint

MIDNIGHT
Powdercoat Paint

HALOGEN

- Material: reflector material embedded into glass | Lexan light shield
Benefit: reflector will not peel or crack | shield is scratch resistant
- Material: oil impregnated bronze bushing | 3-axis light head joint
Benefit: no light drift | 3-axis light keeps light pattern parallel w/the long axis of mouth while freeing up space above patient's head

LED

- Material: energy efficient LED lamps
Benefit: creates a light pattern with feathered edges reducing eye fatigue
consumes approximately 75% less power than halogen lights | no heat
22 -year projected LED life span
- Material: cast aluminum light head
Benefit: allows for nine Designer Friendly[®] powder coat paint colors

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Post Mount Lights



Model 9083

Halogen Light

Part Number 9083



Model 9078

LED Light

Part Number 9078



Model 9073

LED Light

Part Number 9073

Standard Features

Halogen Light 9083

Automatic on/off activation w/override
3-Axis light head movement
Additional bulb
Color shade matching temperature setting
4200 to 5500 Kelvin
Designer Friendly® powder coat paint light head

LED Light 9078

Automatic on/off activation w/override
3-Axis light head movement
Easy release autoclavable handles
Light source 10 LED lamps
Color temperature, 4,200 to 6,000 Kelvin
Designer Friendly powder coat paint light head

LED Light 9073

Automatic on/off activation w/override
3-Axis light head movement
Six-mode, polycarbonate button keyboard
Easy release autoclavable handles
Light source 8 1-watt LED lamps
Color temperature, 4,200 to 6,000 Kelvin
Designer Friendly powder coat paint light head

Factory Options

LIGHT

Replacement bulb (9083/9093)
Short rigid arm 20" in lieu of standard 32"

LIGHT POSTS:

44" with slot for sidebox access
44" without slot, chair mount
35" with post clamp

Pivot (mounts on Forest 3900/3800 chair only)
(does not include cast chair adapter)

Pivot (includes 3900/3800 cast chair adapter)

LIGHT POSTS, WITH MONITOR CORD SLOT:

44" with dual slot for sidebox access
44" with single slot, chair mount
35" with single slot & post clamp

MODEL CHOICE

Part No

0014-300
1110-227-CS
1110-900
1110-794
1110-820
1113-654
1113-654-CA
1110-170
1110-050
1110-045

Designer Friendly® Surface Treatments

STEP 1 SELECT POWDER COAT PAINT FOR UNIT & LIGHT

- ☐ Arctic Silver ☐ Bronze ☐ Charcoal ☐ Copper
☐ Forest Grey ☐ Midnight ☐ Pewter ☐ Shadow ☐ White



Equipment dimensions, see page 98
Shipping weight estimate:
35 lb / 16 kg

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Track Mount Lights



Model 9083T
 Halogen Light/Track
 Shown with Optional 1305T Monitor Mount
 Part Number 9083T



Model 1305T
 Monitor Mount
 Monitor Screen and Track Not Included
 Part Number 1305T



Model 9078T
 LED Light/Track
 Part Number 9078T



Model 9073T
 LED Light/Track
 Part Number 9073T

Standard Features

- Light and Track Assembly**

 - Pre-wired through flex arm
 - Track/Trolley assembly
 - Vertical light tube w/6 inch horizontal side offset for clearance of optional monitor
 - Flex arm & light head assembly
 - Trolley mounting plate and shroud
 - ! Drop post (specify mounting height)
- Monitor Mount Assembly**

 - Pre-wired from track through flex arm
 - Low voltage wiring (2 conductor 18 awg wire)
 - HDMI shielded cable
 - Lower positioning handle
- Monitor Requirements:**

 - 12 volt w/power supply
 - weight: up to 18 lbs max, width of 19"

Factory Options

LIGHT / TRACK	MODEL CHOICE	Part No
Monitor mount assembly		1305T
Replacement bulb (9083/9093)		0014-300
S Video, cable pre-installed		0014-282
RCA cables (3 conductor), pre-installed		0014-283

! Often overlooked

! Standard Mounting Options

STEP 1 CHOOSE YOUR STANDARD MOUNTING OPTION:


MOUNTING HEIGHT:	LIGHT POST LENGTH:	MONITOR POST LENGTH:
<input type="radio"/> 8"	<input type="radio"/> 8.25"	<input type="radio"/> 20.25"
<input type="radio"/> 8' 6"	<input type="radio"/> 14.25"	<input type="radio"/> 26.25"
<input type="radio"/> 9'	<input type="radio"/> 20.25"	<input type="radio"/> 32.25"
<input type="radio"/> 9' 6"	<input type="radio"/> 26.25"	<input type="radio"/> 38.25"
<input type="radio"/> 10'	<input type="radio"/> 32.25"	<input type="radio"/> 44.25"
<input type="radio"/> 10' 6"	<input type="radio"/> 38.25"	<input type="radio"/> 50.25"
<input type="radio"/> 11'	<input type="radio"/> 44.25"	<input type="radio"/> 56.25"
<input type="radio"/> _____ Custom mounting height	<input type="radio"/> _____ * Custom Post Length	<input type="radio"/> _____ * Custom Post Length

* Cannot exceed 44.25" and 56.25"

! Designer Friendly® Surface Treatments

STEP 2 SELECT POWDER COAT PAINT FOR UNIT & LIGHT:

<input type="radio"/> Arctic Silver	<input type="radio"/> Bronze	<input type="radio"/> Charcoal	<input type="radio"/> Copper
<input type="radio"/> Forest Grey	<input type="radio"/> Midnight	<input type="radio"/> Pewter	<input type="radio"/> Shadow
			<input type="radio"/> White



Equipment dimensions, see page 98

Shipping weight estimate:

Light/track	115 lb / 53 kg
1305T	140 lb / 64 kg (includes mount/track light weight)

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cupidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Ceiling Mount Lights



Model 9083-CL
Halogen Light
Part Number 9083-CL



Model 9078-CL
LED Light
Part Number 9078-CL



Model 9073-CL
LED Light
Part Number 9073-CL

Standard Features

- Lights**
Ceiling mounting plate and shroud w/drop post
! Specify mounting height
- Halogen Light 9083-C**
Automatic on/off activation w/override
3-Axis light head movement
Additional bulb
Color temperature, 5,500 Kelvin
Designer Friendly® powder coat paint light head

- LED Light 9078-C**
Automatic on/off activation w/override
3-Axis light head movement
Easy release autoclavable handles
Light source 10 LED lamps
Color temperature, 4,200 to 6,000 Kelvin
Designer Friendly powder coat paint light head

- LED Light 9073-C**
Automatic on/off activation w/override
3-Axis light head movement Six-mode
Polycarbonate button keyboard
Easy release autoclavable handles
Light source 8 1-watt LED lamps
Color temperature, 4,200 to 6,000 Kelvin
Designer Friendly powder coat paint light head

Factory Options

LIGHT	MODEL CHOICE
Replacement bulb (9083/9093)	0014-300
! Often overlooked	

Part No

Standard Mounting Options

STEP 1 CHOOSE YOUR STANDARD MOUNTING OPTION:

MOUNTING HEIGHT:	POST LENGTH:
<input type="radio"/> 8'	<input type="radio"/> 12"
<input type="radio"/> 8' 6"	<input type="radio"/> 18"
<input type="radio"/> 9'	<input type="radio"/> 24"
<input type="radio"/> 9' 6"	<input type="radio"/> 30"
<input type="radio"/> 10'	<input type="radio"/> 36"
<input type="radio"/> 10' 6"	<input type="radio"/> 42"
<input type="radio"/> 11'	<input type="radio"/> 48"
<input type="radio"/> _____ Custom mounting height	<input type="radio"/> _____ * Custom post length
* Cannot exceed 48"	

Designer Friendly® Surface Treatments

STEP 2 SELECT POWDER COAT PAINT FOR UNIT & LIGHT:

<input type="radio"/> Arctic Silver	<input type="radio"/> Bronze	<input type="radio"/> Charcoal	<input type="radio"/> Copper	
<input type="radio"/> Forest Grey	<input type="radio"/> Midnight	<input type="radio"/> Pewter	<input type="radio"/> Shadow	<input type="radio"/> White



Equipment dimensions, see page 98
Shipping weight estimate:
50 lb / 23 kg

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Cabinet and Wall Mount Lights



Model 9079C
 (Standard) 9079C-T (Tall)
 LED Cabinet
 Part Number 9079C



Model 9093C
 (Standard) 9093C-T (Tall)
 Halogen Cabinet
 Part Number 9093C



Model 9079W
 (Standard) 9079W-T (Tall)
 LED Wall
 Part Number 9079W



Model 9074C
 (Standard) 9074C-T (Tall)
 LED Cabinet
 Part Number 9074C



Model 9093W
 (Standard) 9093W-T (Tall)
 Halogen Wall
 Part Number 9093W



Model 9074W
 (Standard) 9074W-T (Tall)
 LED Wall
 Part Number 9074W

Standard Features	Factory Options	Part No
Lights Cabinet mount, internal vertical mount w/post Wall mount, vertical mounting plate Standard rigid arm Tall rigid arm (additional 1' drop)	LIGHT Replacement bulb (9083/9093) Cabinet mounting options: (for dimensions see page 95) Horizontal top cabinet mount in lieu of internal mount Vertical exterior cabinet mount in lieu of internal mount	MODEL CHOICE 0014-300 1110-881-CS 1110-886-CS
Halogen Light 9093 Automatic on/off activation w/override 3-Axis light head movement Additional bulb Color shade matching temperature setting 4200 to 5500 Kelvin Designer Friendly® powder coat paint light head	! ** SPECIFY CABINET/WALL MOUNTING HEIGHT ** ! Often overlooked	
LED Light 9079 Automatic on/off activation w/override 3-Axis light head movement Easy release autoclavable handles Light source 10 LED lamps Color temperature, 4,200 to 6,000 Kelvin Designer Friendly powder coat paint light head		
LED Light 9074 Automatic on/off activation w/override 3-Axis light head movement Six-mode Polycarbonate button keyboard Easy release autoclavable handles Light source 8 1-watt LED lamps Color temperature, 4,200 to 6,000 Kelvin Designer Friendly powder coat paint light head	!	Designer Friendly® Surface Treatments STEP 1 SELECT POWDER COAT PAINT FOR UNIT & LIGHT: <input type="radio"/> Arctic Silver <input type="radio"/> Bronze <input type="radio"/> Charcoal <input type="radio"/> Copper <input type="radio"/> Forest Grey <input type="radio"/> Midnight <input type="radio"/> Pewter <input type="radio"/> Shadow <input type="radio"/> White



Equipment dimensions, see page 99
 Shipping weight estimate:
 Cabinet 60 lb / 28 kg
 Wall 65 lb / 30 kg

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Monitor Mount Assemblies



Model 1300RM Ring Post Mount

Part Number 1300RM



Model 1300TPM Top Post Mount

Part Number 1300TPM

Monitor Screens And Chair-Mounted Posts Not Included

Standard Features

Mount Assembly,
Ring Post Mount 1300RM

ICW Quicklink monitor mount
Rex positioning handles
22" ring mount rigid arm

Monitor Mount Assembly,
Top Post Mount 1300TPM

ICW Quicklink monitor mount
Rex positioning handles
22" top post mount rigid arm

Monitor Requirements:
12 volt with power supply
Weight: 8 1/2 to 13 lbs max, width of 19"

Factory Options

	MONITOR MOUNT ASSEMBLY	MODEL CHOICE	Part No
!	Light Post (1300RM):		
	44" with dual slot, fixed sidebox mount, with Forest light bushing		1110-170
	35" with slot & post clamp, fixed mount, with Forest light bushing		1110-045
	Pivot light posts & requirements		(call to order)
!	MONITOR POST (1300TPM):		
	30" with dual slot, fixed sidebox mount, with Oilite bushing		1331-998
	Pivot monitor posts & requirements		(call to order)
	30" with slot & post clamp, fixed mount, with Oilite bushing		1331-997

! Often overlooked

!

Designer Friendly® Surface Treatments

STEP 1 SELECT POWDER COAT PAINT FOR UNIT & LIGHT:

- ☐ Arctic Silver
- ☐ Bronze
- ☐ Charcoal
- ☐ Copper
- ☐ Forest Grey
- ☐ Midnight
- ☐ Pewter
- ☐ Shadow
- ☐ White

Equipment dimensions, see page 99
Shipping weight estimate:
1300RM, 1300TPM 15 lb / 7 kg

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Forest

7000 SERIES



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

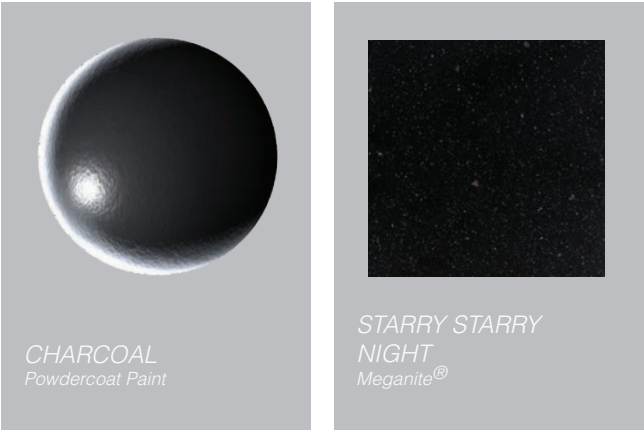
Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty



7000 SERIES REAR MOUNT SIDE DELIVERY

- Material:

Benefit:
- all metal control head chassis
- resistant to UV and cleaners
- Material:

Benefit:
- hinged solid surface lid
- long-lasting surface gives easy access to ancillary products
- Material:

Benefit:
- optional exclusive brass handpiece lockout valves
- allows assistant to remove a handpiece from its holder to change a bur while dentist continues working
- Material:

Benefit:
- four steel mounting configurations
- wall, under cabinet, beneath counter and cart
- Material:

Benefit:
- secondary aluminum work surface
- dedicated work surface for bur block, electric handpiece attachment, mirror, explorer
- Material:

Benefit:
- steel height adjust mechanism
- adjustable work surface and handpiece delivery arm height

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

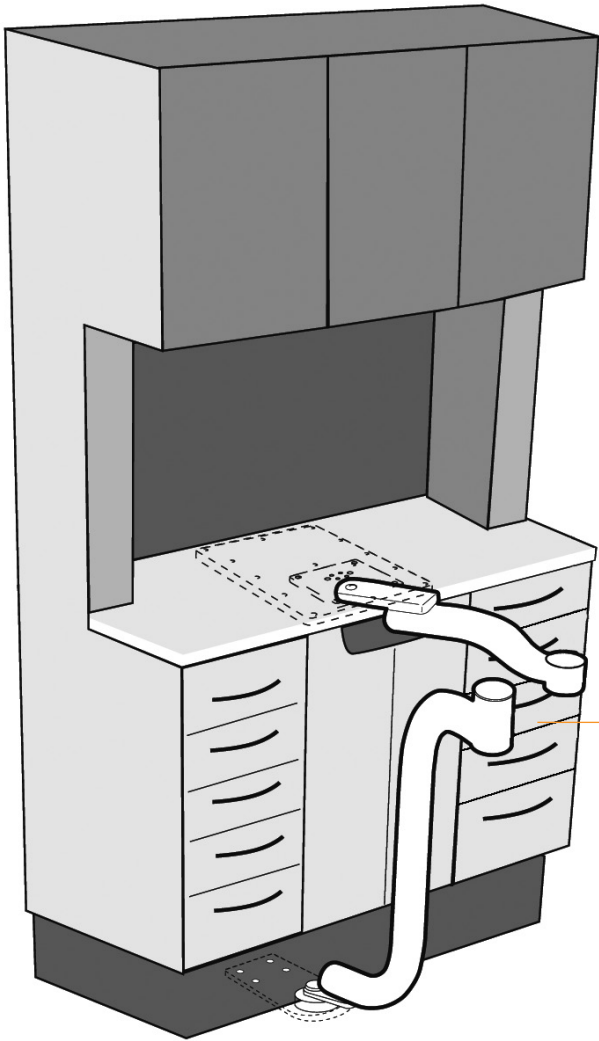
General
Information

Warranty

Cabinetry Mounting Configurations

Made for Each Other

From the ground up, Forest delivery systems are designed for easy cabinet integration.



7000 Series Beneath Counter Mount (BC)

Start with a strong foundation

- Your choice of 3/4" hardened steel 11" x 14-1/2" or optional 10" x 24" mounting plate that bolts to the bottom of your counter top
- Followed by a 6" x 6" leveling plate and a set of needle roller bearings with friction adjustment.
- Durable 2-1/2" steel tube dual arm assembly with needle roller bearings installed between pivoting arm sections.

7000 Series Under Cabinet Floor Mount (UC)

Start with a strong foundation

- Low profile 1/2" steel floor mount plate allows room for your stool casters.
- Needle roller bearings with friction adjustment eliminates drift.
- Durable 2-1/2" steel tube post and arm assembly with needle roller bearings installed between pivoting sections.

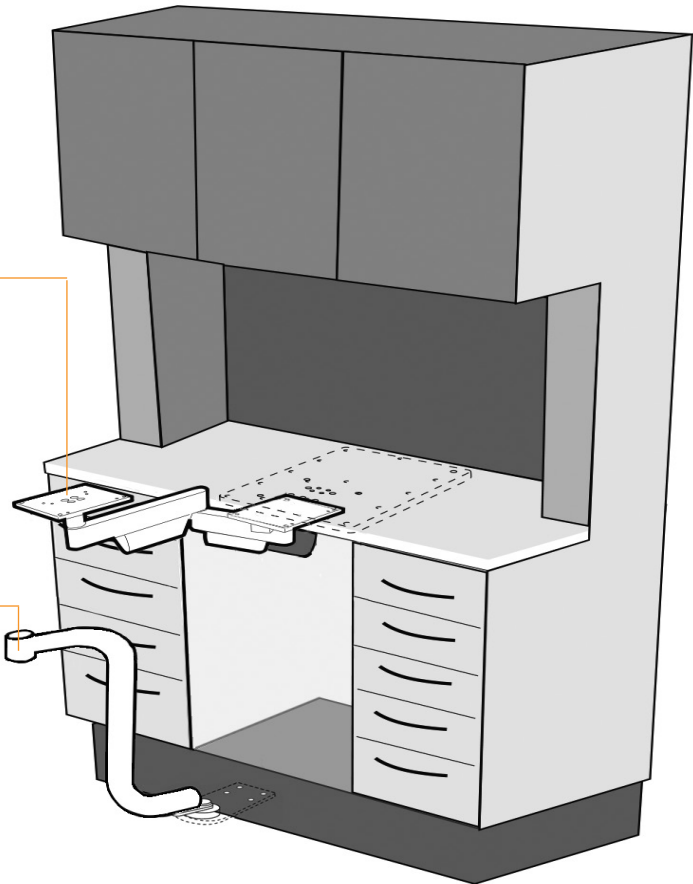
3486 Beneath Counter Mount (BC)

Start with a strong foundation

- Your choice of 3/4" hardened steel 11" x 14-1/2" or optional 10" x 24" mounting plate that bolts to the bottom of your counter top.
- Followed by a 6" x 6" leveling plate with friction adjustment.



3486 Under Cabinet Floor Mount (UC)

- Low profile 1/2" steel floor mount plate allows room for your stool casters.
- Needle roller bearings with friction adjustment eliminates drift.
- Durable 2-1/2" steel tube post and arm assembly with needle roller bearings installed between pivoting sections.



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

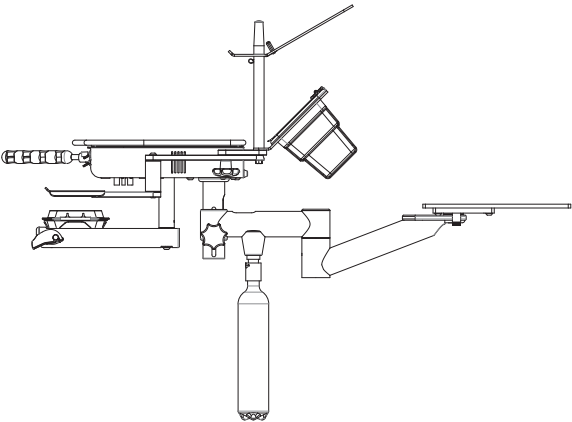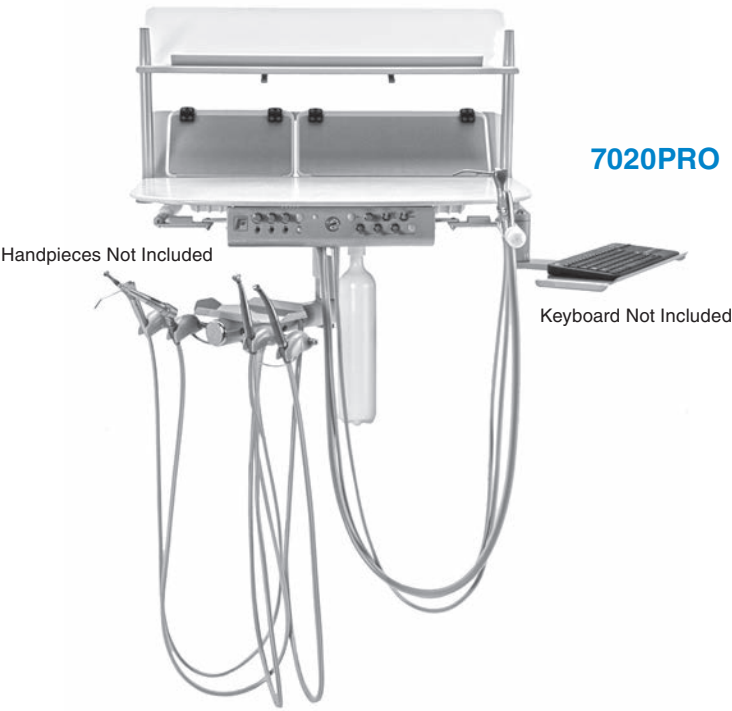7020PRO Duo Swing Delivery System



Model 7020PRO-BC

Beneath Counter Mount

Part Number

7020PRO-BC

Universal heavy-duty mounting plate 11" x 14 1/2" included (mounts to under counter)
Specify Cabinet Manufacturer

Standard Features

Duo Delivery System 7020PRO Base Model

- Duo control
 - automatic control for three handpieces
 - handpiece lockouts
 - 3-way air water syringe w/tubing, 5'
 - 4-hole handpiece tubing, 5'
 - pneumatic master control block w/diaphragm
 - master on/off
 - individual handpiece pressure adjustment
 - individual water coolant adjustment
 - non-retraction water coolant
 - air coolant adjustment
 - handpiece flush
 - handpiece exhaust canister/oil collector
 - wet/dry foot control
- Self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter
- "Burr Buddy" tray holder and covered Zirc tray
- Chrome accents
- Central vacuum accessory package
 - central vacuum canister with disposable solids collector
 - 3-way air/water syringe w/tubing, 5'
 - autoclavable standard HVE lever valve w/tubing, 5'
 - autoclavable saliva ejector lever valve w/tubing, 5'
 - dual assistant's quad holder assembly
- Solid work surface, 17" x 30"
- Pneumatic work surface easy lift/close mechanism
- Height adjustable work surface, 5"
- Large and small organizer bins w/frame, dividers and sharps container
- Procedure shelf
- Keyboard tray
- Stainless trash/CPU tray (MC)
- Universal heavy duty mounting plate 11" x 14 1/2",
 - beneath counter (BC)
 - Wall mount panel 16" x 18" (WM)
- Standard utility center housing and cover (MC)
- Air pressure regulator w/filter
- Master air shut-off valve
- Umbilical, 6'

Mounting plate dimensions, see pages 100

Shipping weight estimate:
7020PRO BC, MC, UC, WM 100 lb / 46 kg



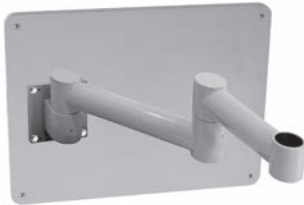
Beneath Counter Mount



Mobile Cart




Under Cabinet Floor Mount




Wall Mount

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

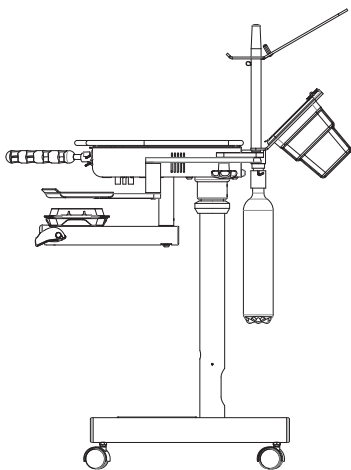
7020PRO Duo Swing Delivery System



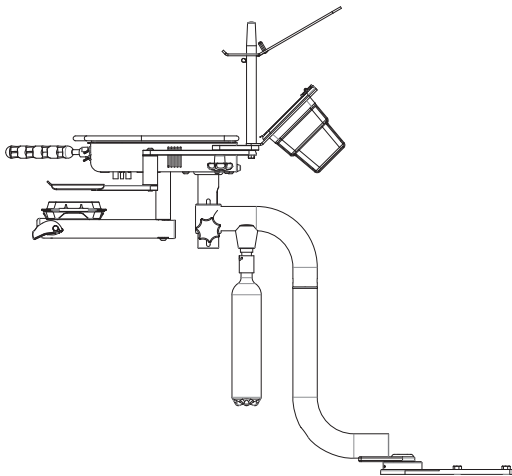Model 7020PRO-MC
Mobile Cart Mount
Part Number 7020PRO-MC



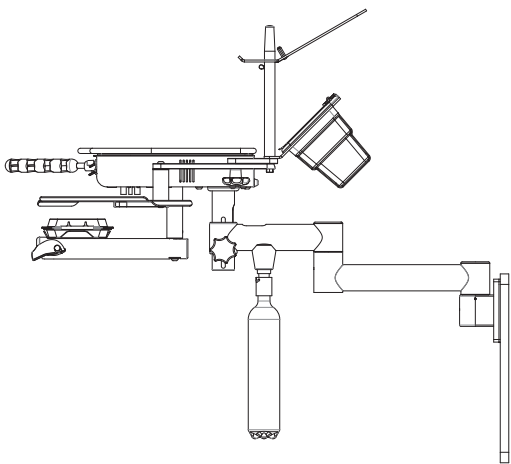Model 7020PRO-UC
Under Cabinet Floor Mount
Part Number 7020PRO-UC



Model 7020PRO-WM
Wall Mount
Part Number 7020PRO-WM

! Designer Friendly® Surface Treatments

STEP 1 SELECT POWDER COAT PAINT:

- ☐ Arctic Silver
- ☐ Bronze
- ☐ Charcoal
- ☐ Copper
- ☐ Forest Grey
- ☐ Midnight
- ☐ Pewter
- ☐ Shadow
- ☐ White

STEP 2 SOLID WORK SURFACE COLORS:

- Corian®

☐ Antarctica
- Wilsonart®:

☐ Arctic Melange☐ Baja Melange
- Meganite®:

☐ Cinder Granite☐ Goldrush Gemstone☐ Moonstone☐ Platinum Granite☐ Rocky Road Granite☐ Silver Mist☐ Starry Starry Night
- Factory Options
- | MONITOR MOUNT ASSEMBLY | MODEL CHOICE | Part No |
|---|-------------------------|---------------|
| Factory Installed STAR ETorque Electric System | | STAR-ETORQUE |
| Handpiece illumination system | | |
| One handpiece (5-hole), accommodates 1 add'l tubing | | 1700-021-CS |
| Add'l handpiece illumination tubing (5-hole) | | 0017-615-CS |
| One handpiece (6-pin), accommodates 1 add'l tubing | | 1700-019-CS |
| Add'l handpiece illumination tubing (6-pin) | | 0017-603-CS |
| ! Ancillaries use existing handpiece position (tubing lengths may vary) (see page 71) | | |
| Write in selection here | | |
| Write in selection here | | |
| ! Handpiece Bundle option (see pages 68-69) | | |
| Write in selection here | | |
| Additional handpiece position, wet | | 1100-189-CS |
| Additional handpiece position, air drive only | | 1100-190-CS |
| Additional standard autoclavable HVE with tubing, 5' & autoclavable oral dry cup | | 1400-113-CS |
| Hands-Free HVE, see page 91, Write in selection here | | |
| Quick disconnect with flow control 1/4" (for water) | | 0016-454-CS |
| Quick disconnect with flow control 3/8" (for air) | | 0016-456-CS |
| ! Monitor mount assembly (monitor screen not included) | | |
| | | 1300-025-CS |
| ! Swivel tray | | |
| | | 1110-909-CS |
| Large bin divider kit (white or black) | 0010-860 / 0010-860-BLK | |
| Small bin divider kit (white or black) | 0010-861 / 0010-861-BLK | |
| Sharps container kit (white or black) | 0010-870 / 0010-870-BLK | |
| ! Touchpad for 3900 chair (25 ft. Cable) | | |
| | | 1700-081-1-CS |
| Standard utility center housing & cover (BC/UC/WM) | | 1115-020 |
| In floor recessed utility center housing & cover | | 1115-305 |
| Additional umbilical (per foot) | | 1115-999-CS |
| Routing valve kit with water shutoff & regulator/factory installation | | 1115-315-CS |
| Universal mounting plate, 10" x 24" in lieu of 11" x 14 1/2" (BC) | | 0010-084-CS |
| Δ Special color solid work surface (call for price and leadtime) | | |
| | | 0023-856-SSV |
- Cabinet mounting configurations, see page 53
- ! Often overlooked


Δ Lead time varies
- TABLE OF CONTENTS
- | |
|--|
| Pivot & Fixed Chair Mount Packages |
| Pivot & Fixed Chair Mount Package Ordering Guide |
| Fusion Packages |
| Fusion Package Ordering Guide |
| Designer Friendly Colors & Surfaces |
| Chairs & Stools |
| Lights and Monitor Mounts |
| Rear Mount Side Delivery Systems |
| Ancillary Ordering & Options |
| Control Head Accessories |
| Simplified Cuspidors & Assistant's Instrumentation |
| Simplified Rear & Side Delivery Systems |
| Factory Options & Parts |
| General Information |
| Warranty |
- 55

Rear Mount Side Delivery Systems

7021PRO Duo Swing Delivery System



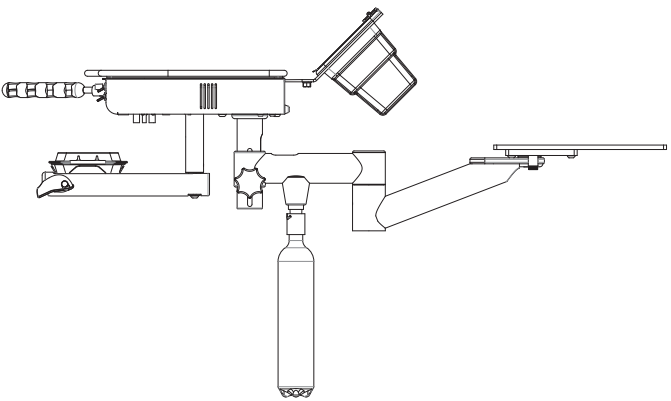7021PRO



Model 7021PRO-BC

Beneath Counter Mount

Part Number

7021PRO-BC

Universal heavy-duty mounting plate 11" x 14 1/2" included (mounts to under counter)

Specify Cabinet Manufacturer

Standard Features

Duo Delivery System 7021PRO

- Duo control
 - automatic control for three handpieces
 - handpiece lockouts
 - 3-way air water syringe w/tubing, 5'
 - 4-hole handpiece tubing, 5'
 - pneumatic master control block w/diaphragm
 - master on/off
 - individual handpiece pressure adjustment
 - individual water coolant adjustment
 - non-retraction water coolant
 - air coolant adjustment
 - handpiece flush
 - handpiece exhaust canister/oil collector
 - wet/dry foot control
- Self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter
- "Burr Buddy" tray holder and covered Zirc tray
- Chrome accents
- Central vacuum accessory package
 - central vacuum canister with disposable solids collector
 - 3-way air/water syringe w/tubing, 5'
 - autoclavable standard HVE lever valve w/tubing, 5'
 - autoclavable saliva ejector lever valve w/tubing, 5'
 - dual assistant's quad holder assembly
- Solid work surface, 17" x 30"
- Pneumatic work surface easy lift/close mechanism
- Height adjustable work surface, 5"
- Large and small organizer bins w/frame, dividers and sharps container
- Universal heavy duty mounting plate 11" x 14 1/2", beneath counter (BC)
- Wall mount panel 16" x 18" (WM)
- Standard utility center housing and cover (MC)
- Air pressure regulator w/filter
- Master air shut-off valve
- Umbilical, 6'



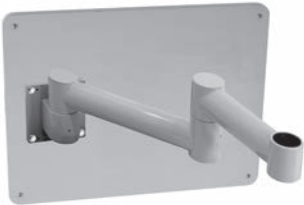
Beneath Counter Mount



Mobile Cart



Under Cabinet Floor Mount




Wall Mount

Mounting plate dimensions, see pages 96-105


Shipping weight estimate:

7021PRO BC, MC, UC, WM 100 lb / 46 kg

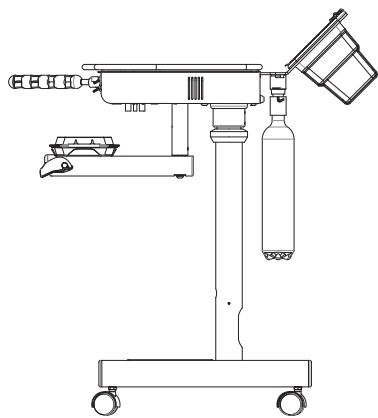
7021PRO Duo Swing Delivery System



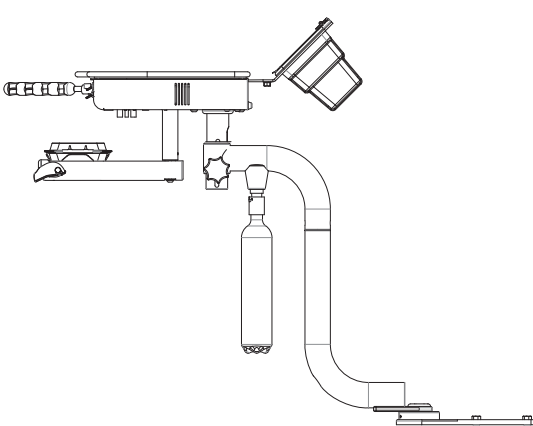Model 7021PRO-MC
Mobile Cart Mount
Part Number 7021PRO-MC



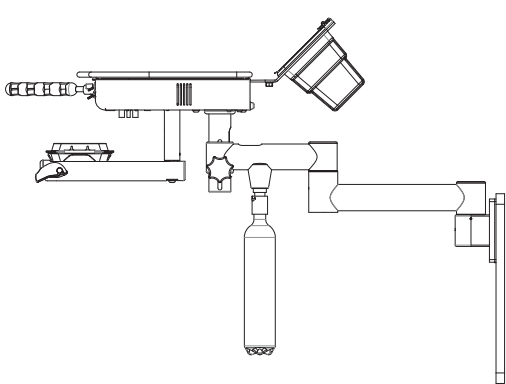Model 7021PRO-UC
Under Cabinet Floor Mount
Part Number 7021PRO-UC



Model 7021PRO-WM
Wall Mount
Part Number 7021PRO-WM

Factory Options		Part No
<div><div>! Designer Friendly® Surface Treatments</div><div>STEP 1 SELECT POWDER COAT PAINT:<div><div><div><input type="radio"/> Arctic Silver</div><div><input type="radio"/> Copper</div><div><input type="radio"/> Pewter</div></div><div><div><input type="radio"/> Bronze</div><div><input type="radio"/> Forest Grey</div><div><input type="radio"/> Shadow</div></div><div><div><input type="radio"/> Charcoal</div><div><input type="radio"/> Midnight</div><div><input type="radio"/> White</div></div></div><div>STEP 2 SOLID WORK SURFACE COLORS:<div><div><div>Corian®</div><div><input type="radio"/> Antarctica</div></div><div><div>Wilsonart®:</div><div><div><input type="radio"/> Arctic Melange</div><div><input type="radio"/> Baja Melange</div></div></div><div><div>Meganite®:</div><div><div><input type="radio"/> Cinder Granite</div><div><input type="radio"/> Goldrush Gemstone</div><div><input type="radio"/> Moonstone</div><div><input type="radio"/> Platinum Granite</div><div><input type="radio"/> Rocky Road Granite</div><div><input type="radio"/> Silver Mist</div><div><input type="radio"/> Starry Starry Night</div></div></div></div></div></div></div>	DELIVERY SYSTEM	MODEL CHOICE
	Factory Installed STAR ETorque Electric System	STAR-ETORQUE
	Handpiece illumination system	
	One handpiece (5-hole), accommodates 1 add'l tubing	1700-021-CS
	Add'l handpiece illumination tubing (5-hole)	0017-615-CS
	One handpiece (6-pin), accommodates 1 add'l tubing	1700-019-CS
	Add'l handpiece illumination tubing (6-pin)	0017-603-CS
	! Ancillaries use existing handpiece position (tubing lengths may vary) (see page 71)	
	Write in selection here	
	Write in selection here	
	! Handpiece Bundle option (see pages 68-69)	
	Write in selection here	
	Additional handpiece position, wet	1100-189-CS
	Additional handpiece position, air drive only	1100-190-CS
	Additional standard autoclavable HVE with tubing, 5' & autoclavable oral dry cup	1400-113-CS
	Hands-Free HVE, see page 91, Write in selection here	
	Quick disconnect with flow control 1/4" (for water)	0016-454-CS
	Quick disconnect with flow control 3/8" (for air)	0016-456-CS
	! Procedure shelf	0010-864-CS
	! Keyboard tray	0010-862-CS
	! Stainless trash/CPU tray, H-frame (MC)	0010-863
	! Monitor mount assembly (monitor screen not included) base model requires purchase of 0010-864-CS	1300-025-CS
	! Swivel tray	1110-909-CS
	Large bin divider kit (white or black)	0010-860 / 0010-860-BLK
	Small bin divider kit (white or black)	0010-861 / 0010-861-BLK
	Sharps container kit (white or black)	0010-870 / 0010-870-BLK
	! Touchpad for 3900 chair (25 ft. Cable)	1700-081-1-CS
	Standard utility center housing & cover (BC/UC/WM)	1115-020
	In floor recessed utility center housing & cover	1115-305
	Additional umbilical (per foot)	1115-999-CS
	Routing valve kit with water shutoff & regulator/factory installation	1115-315-CS
	Universal mounting plate, 10" x 24" in lieu of 11" x 14 1/2" (BC)	0010-084-CS
	Δ Special color solid work surface (call for price and leadtime)	0023-856-SSV

Cabinet mounting configurations, see page 53

! Often overlooked
Δ Lead time varies

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Rear Mount Side Delivery Systems

7010 Duo Swing Delivery System

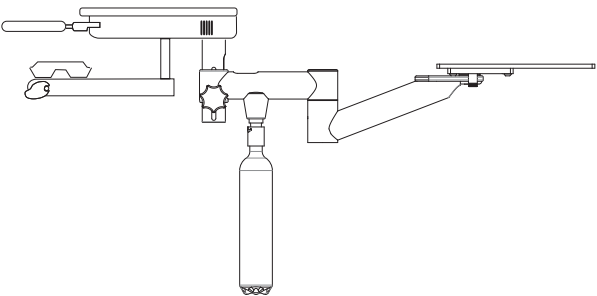


1110-185-CS and 1110-908-CS

Shown with optional storage organizer and swivel tray



Shown with optional handpiece lockouts 0010-748 and additional assistant's quad holder assembly 0024-156



Model 7010BC

Beneath Counter Mount

Part Number

7010BC

Universal heavy-duty mounting plate 11" x 14 1/2" included (mounts to under counter)
Specify Cabinet Manufacturer

Standard Features

Duo Swing Delivery System 7010

- Duo control
- automatic control for three handpieces
 - 3-way air water syringe w/tubing, 4'
 - 4-hole handpiece tubing, 4'
 - pneumatic master control block w/diaphragm
 - master on/off
 - individual handpiece pressure adjustment
 - individual water coolant adjustment
 - non-retraction water coolant
 - air coolant adjustment
 - handpiece flush
 - handpiece exhaust canister/oil collector
 - self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter
 - wet/dry foot control
- Height adjustable work surface, 5" (BC/UC/WM)
- "Burr Buddy" tray holder and covered Zirc tray
- Pneumatic work surface easy lift/close mechanism
- Central vacuum accessory package
- central vacuum canister w/disposable solids collector
 - 3-way air/water syringe w/tubing, 5'
 - autoclavable standard HVE lever valve w/tubing, 5'
 - autoclavable saliva ejector lever valve w/tubing, 5'
- Solid work surface, 19" x 24"
- Universal heavy duty mounting plate 11" x 14 1/2", beneath counter (BC)
- Wall mount panel 16" x 18" (WM)
- Adjustable height frame (MC)
- Twin wheel hard floor casters (MC)
- Standard utility center housing and cover (MC)
- Air pressure regulator w/filter
- Master air shut-off valve
- Umbilical, 7'
- External umbilical, 6' (MC)



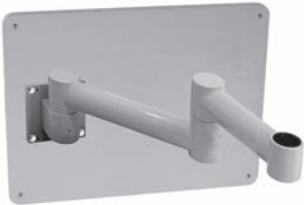
Beneath Counter Mount



Mobile Cart



Under Cabinet Floor Mount




Wall Mount

Mounting plate dimensions, see pages 102-105
Shipping weight estimate:
7010BC, MC, UC, WM 130 lb / 59 kg

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Rear Mount Side Delivery Systems

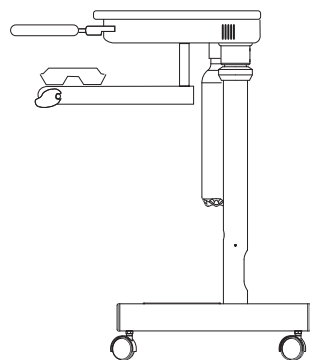
7010 Duo Swing Delivery System



Model 7010MC

Mobile Cart Mount

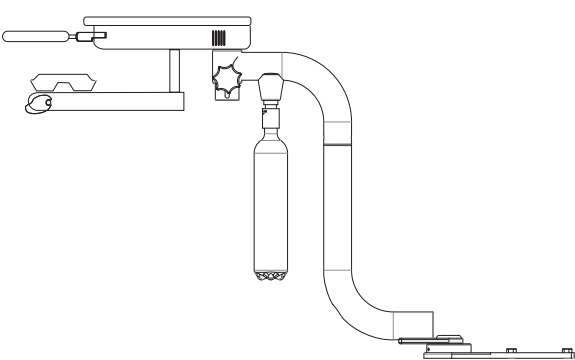
Part Number 7010MC



Model 7010UC

Under Cabinet Floor Mount

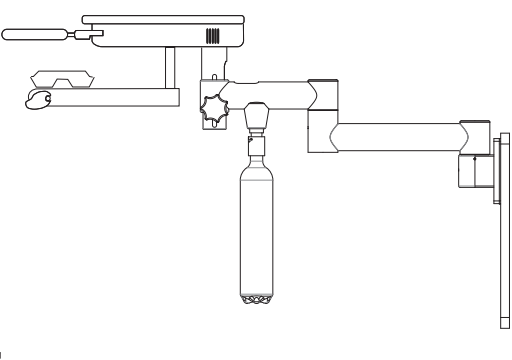Part Number 7010UC



Model 7010WM

Wall Mount

Part Number 7010WM

Factory Options		Part No
<div><div>! Designer Friendly® Surface Treatments</div><div><div>STEP 1 SELECT POWDER COAT PAINT:</div><div><div><div><div><input type="radio"/> Arctic Silver</div><div><input type="radio"/> Bronze</div><div><input type="radio"/> Charcoal</div></div><div><div><input type="radio"/> Copper</div><div><input type="radio"/> Forest Grey</div><div><input type="radio"/> Midnight</div></div><div><div><input type="radio"/> Pewter</div><div><input type="radio"/> Shadow</div><div><input type="radio"/> White</div></div></div><div><div>STEP 2 SOLID WORK SURFACE COLORS:</div><div><div>Corian®<div><input type="radio"/> Antarctica</div></div></div><div><div>Wilsonart®:<div><div><input type="radio"/> Arctic Melange</div><div><input type="radio"/> Baja Melange</div></div></div></div><div><div>Meganite®:<div><div><input type="radio"/> Cinder Granite</div><div><input type="radio"/> Goldrush Gemstone</div><div><input type="radio"/> Moonstone</div><div><input type="radio"/> Platinum Granite</div><div><input type="radio"/> Rocky Road Granite</div><div><input type="radio"/> Silver Mist</div><div><input type="radio"/> Starry Starry Night</div></div></div></div></div></div></div></div>		
DELIVERY SYSTEM	MODEL CHOICE	
Factory Installed STAR ETorque Electric System	STAR-ETORQUE	
Handpiece illumination system		
One handpiece (5-hole), accommodates 1 add'l tubing	1700-021-CS	
Add'l handpiece illumination tubing (5-hole)	0017-615-CS	
One handpiece (6-pin), accommodates 1 add'l tubing	1700-019-CS	
Add'l handpiece illumination tubing (6-pin)	0017-603-CS	
! Ancillaries use existing handpiece position (tubing lengths may vary) (see page 71)		
Write in selection here		
Write in selection here		
! Handpiece Bundle option (see pages 68-69)		
Write in selection here		
Additional handpiece position, wet	1100-189-CS	
Additional handpiece position, air drive only	1100-190-CS	
Lockout feature for 3 handpieces	0010-748	
Additional standard autoclavable HVE with tubing, 5' & autoclavable oral dry cup	1400-113-CS	
Hands-Free HVE, see page 91, Write in selection here		
Additional assistant's quad holder assembly	0024-156	
Quick disconnect with flow control 1/4" (for water)	0016-454-CS	
Quick disconnect with flow control 3/8" (for air)	0016-456-CS	
! Storage organizer with bins & frame	1110-185-CS	
! Swivel tray (requires purchase of 1110-185-CS)	1110-908-CS	
! Swivel tray (mounts without storage organizer)	1110-079-CS	
! Monitor mount assembly (monitor screen not included) base model requires purchase of 0010-864-CS	1300-025-CS	
Swivel tray	1110-909-CS	
Touchpad for 3900 chair (25 ft. Cable)	1700-081-1-CS	
Chrome accents	9000-998-CS	
Standard utility center housing & cover (BC/UC/WM)	1115-020	
In floor recessed utility center housing & cover	1115-305	
Additional umbilical (per foot)	1115-999-CS	
Routing valve kit with water shutoff & regulator/factory installation	1115-315-CS	
Universal mounting plate, 10" x 24" in lieu of 11" x 14 1/2" (BC)	0010-084-CS	
Standard color laminate work surface in lieu of solid surface	0023-776-LAM	
Δ Special color solid or laminate work surface (call for price and leadtime)	0023-776-SSV/0023-776-X	

Cabinet mounting configurations, see page 53

! Often overlooked

Δ Lead time varies

Pivot & Fixed Chair Mount Packages

Pivot & Fixed Chair Mount Package Ordering Guide

Fusion Packages

Fusion Package Ordering Guide

Designer Friendly Colors & Surfaces

Chairs & Stools

Lights and Monitor Mounts

Rear Mount Side Delivery Systems

Ancillary Ordering & Options

Control Head Accessories

Simplified Cuspidors & Assistant's Instrumentation

Simplified Rear & Side Delivery Systems

Factory Options & Parts

General Information

Warranty

3486 Swing Delivery System



1110-179-CS and 1110-908-CS

Shown with optional storage organizer and swivel tray



Model 3486BC

Beneath Counter Mount

Part Number

3486BC

Universal heavy-duty mounting plate 11" x 14 1/2" included (mounts to under counter)
Specify Cabinet Manufacturer

Standard Features

Assistant's Swing Delivery System 3486 Series

- Central vacuum accessory package
 - arm-mounted central vacuum canister w/disposable solids collector
 - 3-way air/water syringe w/tubing, 5'
 - autoclavable standard HVE lever valve w/tubing, 5'
 - autoclavable saliva ejector lever valve w/tubing, 5'
 - self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter
- Assistant's quad holder assembly
- Solid work surface, 19" x 24"
- Beneath counter rigid arm assembly
 - w/heavy duty mounting plate, 11" x 14 1/2" (BC)
- Adjustable height frame (MC)
- Twin wheel hard floor casters (MC)
- Under cabinet floor mount delivery post (UC)
- Wall mount panel, 16" x 18" (WM)
- Standard utility center housing and cover (MC)
- Air pressure regulator w/filter
- Master shut-off valve
- Umbilical, 7'
- electrical outlet/UL (MC)

Note: Arms can be reversed on the beneath counter and wall mount to adjust work surface height.



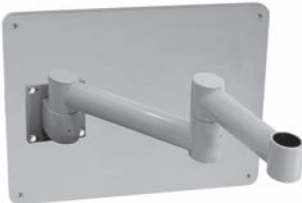
Beneath Counter Mount



Mobile Cart



Under Cabinet Floor Mount




Wall Mount

Mounting plate dimensions, see pages 107-108
Shipping weight estimate:
3486 115 lb / 53 kg

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Storage Organizer, Trays and Monitor Mount

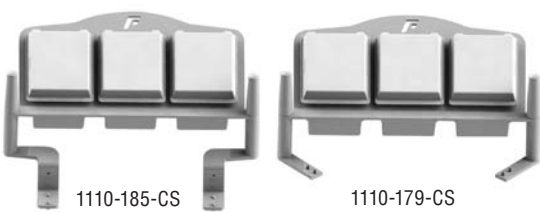Storage Organizer with Bins and Frame



Factory option model: 7010
Factory option models: 2184BC/2186BC/3486

1110-185-CS
1110-179-CS

- Metal powder coated organizer frame
- Three storage bins with removable covers
- Storage dividers

Swivel Tray



Factory option models:
2184BC/2186BC/3486/7010 requires purchase of
1110-179-CS or 1110-185-CS

- Inserts inside storage organizer tower
- Metal powder coated tray holder
- Plastic instrument tray

1110-908-CS

Swivel Tray with Left or Right Mounting Bracket



Factory option model: 7010
Factory option models: 2184BC/2186BC/3486

1110-079-CS
1110-482-CS

- Metal powder coated organizer frame
- Three storage bins with removable covers
- Storage dividers

Swivel Tray with Left and Right Mounting Brackets



Factory option models:

- Metal powder coated tray holder towers and brackets
- Plastic instrument tray

PRO Series and 7010 third tray addition mounts with organizer

1110-909-CS

Monitor Mount



Factory option models:

- Inserts inside storage organizer or PRO Series procedure shelf tower
- Forest's unique locking tilt mechanism with knob 1/2" pin

1300-025-CS

Monitor requirements (not included):
12 volt with power supply
Weight 8-1/4 to 13 lbs. max, width of 19"

PRO Series 7021 requires purchase of 0010-864-CS
7010 requires purchase of 1110-185-CS
2184BC/2186BC/3486 requires purchase of 1110-179-CS

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Forest

DELIVERY SYSTEMS



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

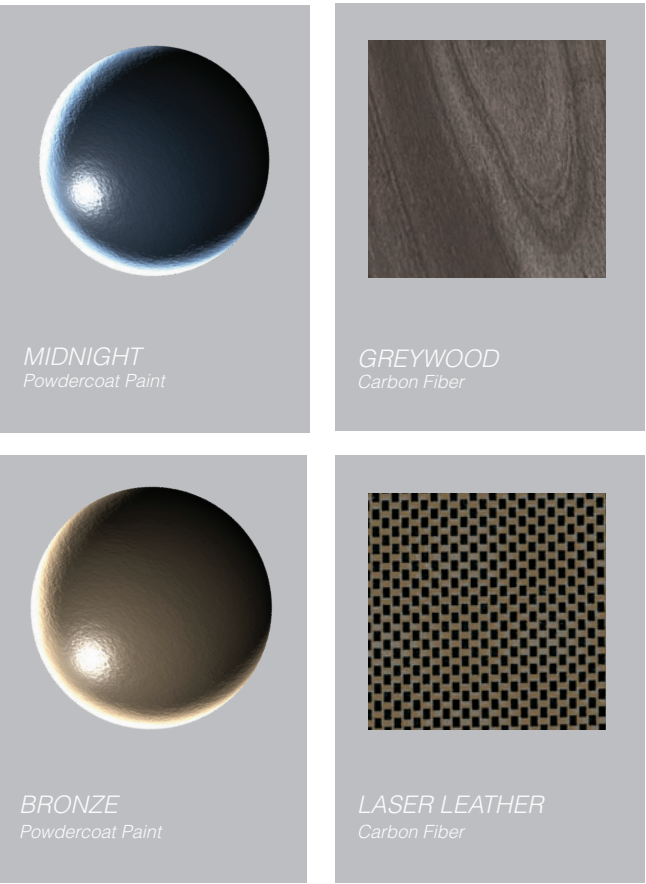
Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty



STANDARD IC

- Material: aluminum hinged lid 1/8" ABS plastic chassis with durable DR acrylic protective coating
Benefit: hinge allows quick access to controls and ancillary devices / lid and chassis are resistant to UV, impact and compatible surface cleaners / aluminum lid is strong and lightweight
- Material: electro-less nickel plated brass air and water Ts, routing valves, pneumatic control and distribution blocks
Benefit: 20 year limited warranty on all parts "under the hood" (on all control heads); refer to page 107 for warranty details
trouble-free handpiece operation
- Material: handpiece holder bar with glass-filled nylon injection molded handpiece holders
Benefit: allows for up to four handpiece positions, a syringe and ancillary touchpad at your desired angle and spacing
- Material: hinged aluminum control head lid
Benefit: available in nine two-coat powder coat paint colors
- Material: plastic chassis
Benefit: available in three carbon fiber look and two solid color finishes

EURO

- Material: whip arms rotate on nylon wheels with aluminum and brass components
Benefit: feather-touch whip movement / ultralight handpiece feel / long, trouble free life
- Material: control head swivels at its center point
Benefit: 240-degrees of control head rotation allows for beside-or over- the-patient placement
- Material: forward mounted aluminum whip arm system
Benefit: allows for greater reach of handpiece and syringe expanding control head placement options.
- Material: plastic control head lid
Benefit: available in three carbon fiber looks and two solid color finishes
- Feature: autoclavable silicone handpiece pad available in black and white
Benefit: blends with Designer Friendly® choices

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Standard IC Delivery Unit

Standard IC


Please use graph below to determine number of ancillaries and options possible for the Standard IC.



Standard Features

- Automatic control for three handpieces
- 3-way air/water syringe with tubing, 5'
- 4-hole handpiece tubing, 5'
- pneumatic master control block with diaphragm master on/off
- individual handpiece pressure adjustment
- individual water coolant adjustment
- non-retraction water coolant
- air/water coolant adjustment
- handpiece flush
- handpiece exhaust canister/oil collector
- handle brake release function

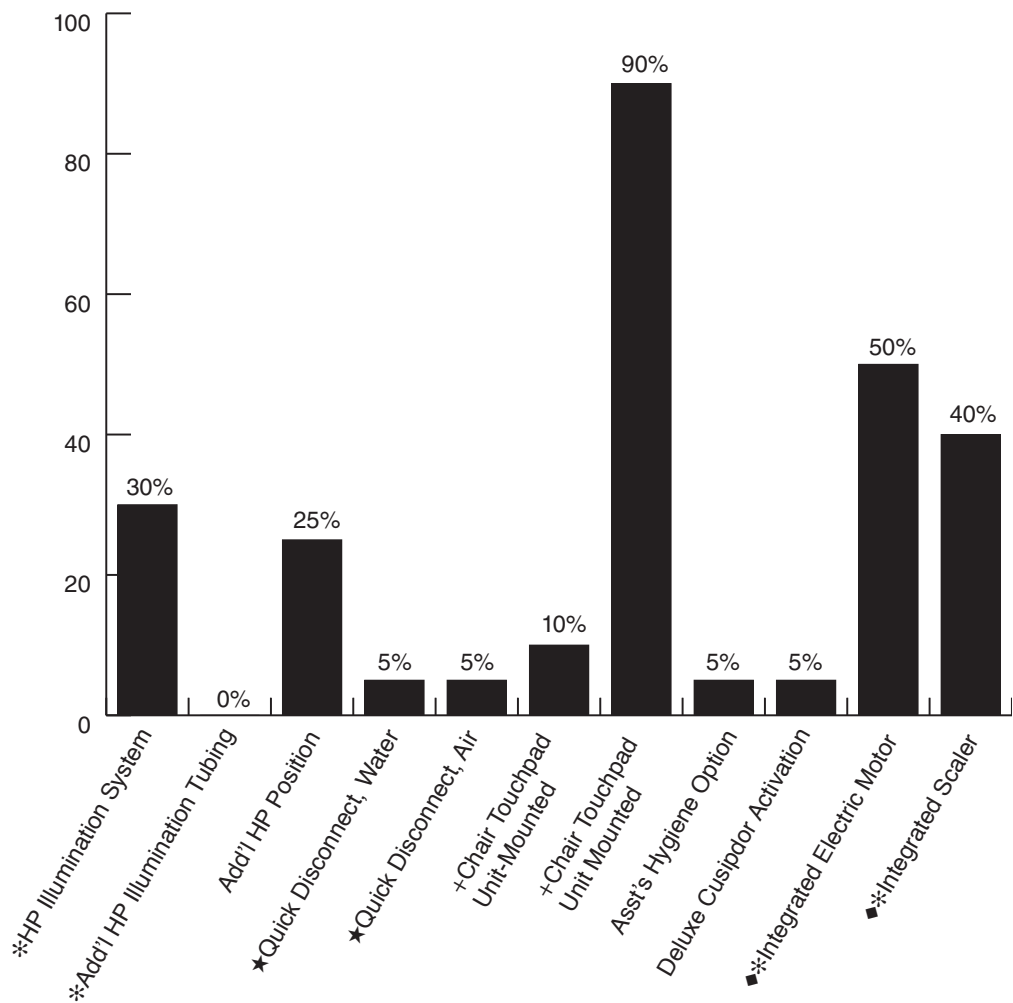% Standard IC Control Head Space Absorbed per Option



Please note:
If your ancillary and option choices exceed 100%, please call customer service for additional information.
* Uses existing handpiece position
★ Standard on SE and SI deliveries
+ 3900 Chair touchpad cable standard on PI and SI deliveries
◆ Two or more may require purchase of 0010-310 ancillary shelf (see page 72)

- Pivot & Fixed Chair Mount Packages
- Pivot & Fixed Chair Mount Package Ordering Guide
- Fusion Packages
- Fusion Package Ordering Guide
- Designer Friendly Colors & Surfaces
- Chairs & Stools
- Lights and Monitor Mounts
- Rear Mount Side Delivery Systems
- Ancillary Ordering & Options
- Control Head Accessories
- Simplified Cuspidors & Assistant's Instrumentation
- Simplified Rear & Side Delivery Systems
- Factory Options & Parts
- General Information
- Warranty

Euro Delivery Unit

Euro


Please use graph below to determine number of ancillaries and options possible for the Euro.



Standard Features

- Automatic control for four handpieces
- 3-way air/water syringe with tubing, 3'
- 4-hole handpiece tubing, 3'
- pneumatic master control block with diaphragm
- master on/off
- individual handpiece pressure adjustment
- individual water coolant adjustment
- non-retraction water coolant
- air/water coolant adjustment
- handpiece flush
- handpiece exhaust canister/oil collector
- dual handle brake release function
- swivel instrument tray

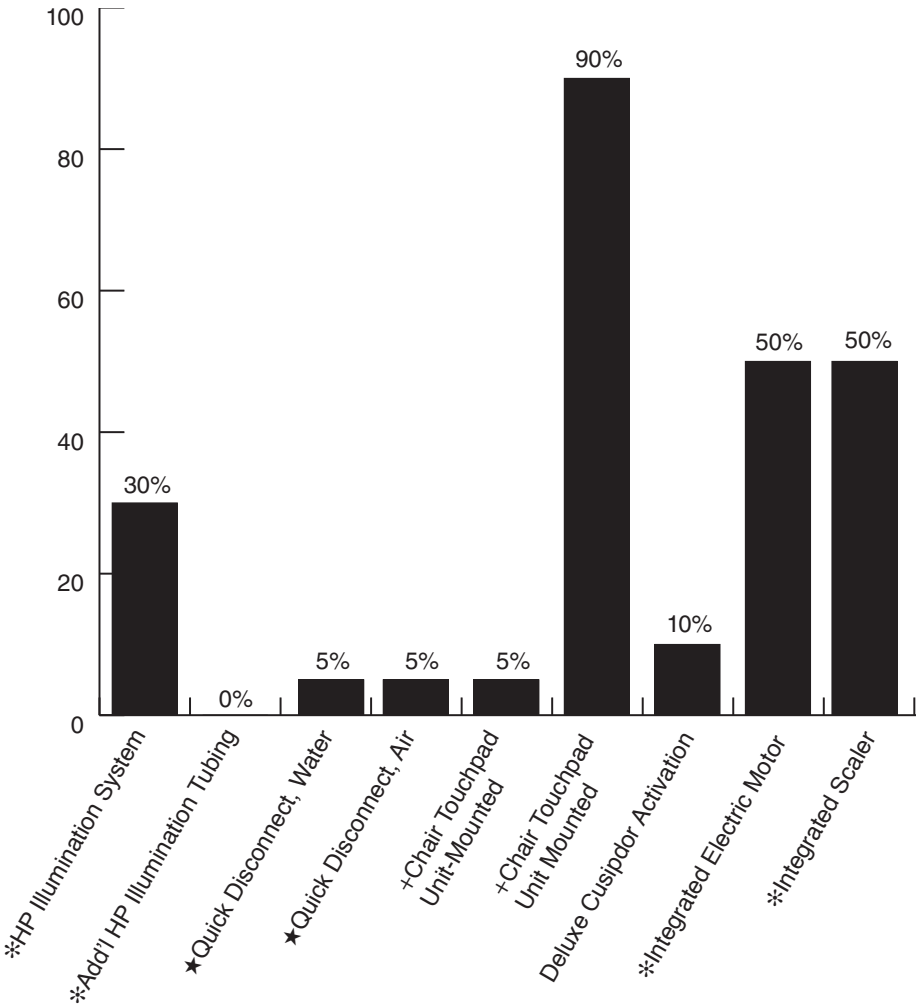% Euro Control Head Space Absorbed per Option



Please note:

If your ancillary and option choices exceed 100%, please call customer service for additional information.

- * Uses existing handpiece position
- ★ Standard on SE and SI deliveries
- + 3900 Chair touchpad cable standard on PI and SI deliveries

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Fusion Delivery Unit

Fusion


Please use graph below to determine number of ancillaries and options possible for the delivery system.



Standard Features

- Automatic control for three handpieces
- 3-way air/water syringe w/tubing, 5'
- 4-hole handpiece tubing, 5'
- pneumatic master control block w/diaphragm master on/off
- individual handpiece pressure adjustment
- individual water coolant adjustment
- non-retraction water coolant
- air coolant adjustment
- handpiece flush
- handpiece exhaust canister/oil collector
- handle brake release function

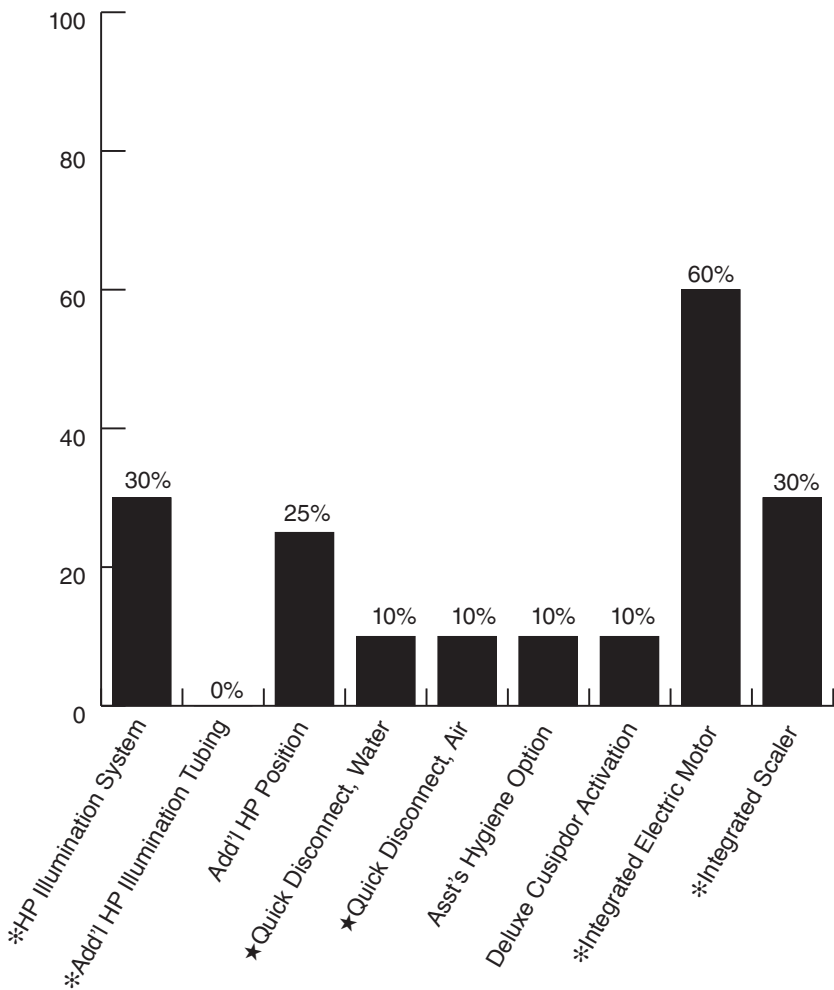% Control Head Space Absorbed per Option



Please note:
If your ancillary and option choices exceed 100%, please call customer service for additional information.

* Uses existing handpiece position

★ Standard on SI deliveries

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Handpiece Bundles

Save time and money with handpiece bundles

Our Star® handpiece bundles make it easy to be set up for success. Each bundle package is designed to complete a new operatory setup with easy ordering and tremendous savings.

Bundle Packages		Bundle Details	Part No	Quantity
Star 500 fiber optic Bundle		Star 500 Elite LubeFree, fiber optic	266040	3
Part Number	268925	Star 500 Swivel, fiber optic	266064	1
430 Torque fiber optic Bundle		430SWL Torque LubeFree, fiber optic	265890	3
Part Number	268927	LED Plug 'N Play Swivel	265646	1
430 Series fiber optic Bundle		430SWL Satin LubeFree, fiber optic	268840	3
Part Number	268928	LED Plug 'N Play Swivel	265646	1
430K Series fiber optic Bundle		430K LubeFree, fiber optic	262052	3
Part Number	268929			

Star fiber optic bundles require 6-pin tubing configuration and fiber optic light source
All other fiber optic bundles require fiber optic light source

NOTES:
Cannot be combined with other Star promotions or discounts.
*Savings calculated off the total retail price for the individual items in each package.

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Bundle Packages		Bundle Details	Part No	Quantity
Star Titan® Low-Speed Bundle	Part Number 268934	Star Titan 20K Motor	266271	1
		Titan 3 Motor To Angle Adaptor	264073	3
		Titan LubeFree Auto Latch Angle	264146	3
		Titan 3 Straight Nose Cone	257506	1
		Titan 4-Line Swivel	260144	1
Titan 3 Low-Speed Bundle	Part Number 268930	Titan 3 20K Motor	262262	1
		Titan 3 Motor to Angle Adaptor	264073	3
		Titan LubeFree Auto Latch Angle	264146	3
		Titan 3 Straight Nose Cone	257506	1
		Titan 4-Line Swivel	260144	1
Star 500 non-fiber optic Bundle	Part Number 268931	Star 500 Elite LubeFree, non-fiber optic	266041	3
		Star 500 Swivel, non-fiber optic	266065	1
430 Torque non-fiber optic Bundle	Part Number 268932	430SW Torque LubeFree	265891	3
		LED Plug ‘N Play Swivel	265646	1
430 Series non-fiber optic Bundle	Part Number 268933	430SW Satin LubeFree	268841	3
		LED Plug ‘N Play Swivel	265646	1

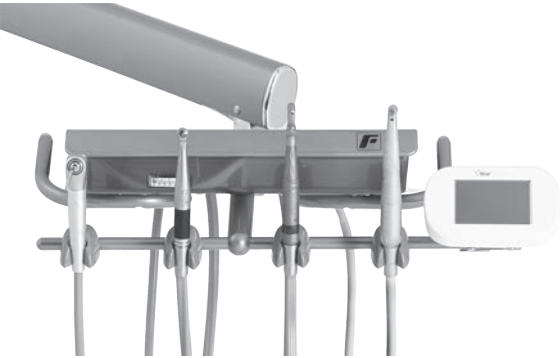
Star fiber optic bundles require 6-pin tubing configuration and fiber optic light source
All other fiber optic bundles require fiber optic light source

NOTES:
Cannot be combined with other Star promotions or discounts.
*Savings calculated off the total retail price for the individual items in each package.

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant’s Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Ancillary Products

Star ETorque™ Electric System




Factory Option

STAR-ETORQUE

Quiet and powerful, with a brushless motor reaching up to 200,000 rpm, this electric system can perform all procedures. It has a user-friendly display module that can be customized as needed and has a fully autoclavable motor. Handpiece attachments not included in price.

To order StarETorque and accessories, see the Star handpiece catalog.

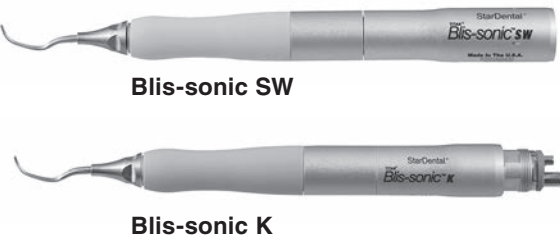Star® Titan® Blis-sonic™ Scaler



Blis-sonic SW

Blis-sonic K

Blis-sonic SW Scaler
(Includes 4 line swivel, Star PN 261547)

STAR-BLIS-SW

Blis-sonic K Scaler (Fixed Backed)

STAR-BLIS-K

Universal scaler tip

- Use for all Surfaces, medium and heavy calculus

Sickle scaler tip

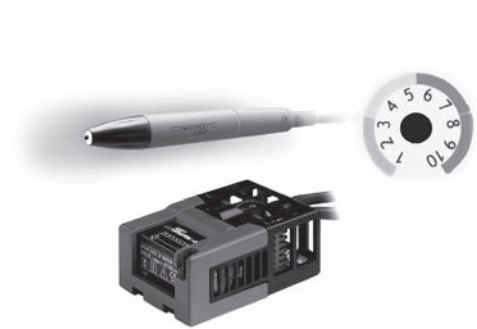
- Recommended for sensitive patients
- Less energetic than the Universal Tip
- Less suitable for deep root scaling

Perio scaler tip

- Recommended for lingual supragingival scaling of mandibular incisors and deep local pockets
- Provides same energy as Universal Tip

To order Titan Blis-sonic scalers and accessories, see the Star handpiece catalog.

Satelec SP



Factory option

SATELEC-S

Factory option

SATELEC-S-LED

- Uses existing handpiece position
- Satelec SP Newtron LED lit Piezo Ultrasonic Scaler w/power supply
- Color coding system ensures proper tip and power setting
- Three tips: #1, #10Z, H3 (Intro Scaling Kit)
- Tips available for periodontics, endodontics, scaling, restorative, and cavity preparation
- Cruise Control® automatically adjusts frequency amplitude
- On/off controlled by rheostat
- First ever autoclavable LED lit handpiece – 100,000 LUX
- Newtron technology means power, control and reliability
- Satelec warranty 1 year / tips 30 days from date of purchase
- Please call 800-289-6367 for warranty and technical support of the Scaler

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Ancillary Products

Cavitron Scaler



Factory option

CAVITRON

- Uses existing handpiece position Dentsply Cavitron Magnetostrictive Ultrasonic scaler
- Blue Zone® for patient comfort
- SPS® (Sustained Performance System)
constant balance between scaling efficiency and patient comfort
- Handpiece w/330 degree swivel and lavage control (reduces cord drag)
- Autoclavable Steri-mate® (infection control)
- PowerLINE 1000 30K insert
- Frequency 30 kHz
- Cavitron warranty 2 years from date of purchase
- Please call 800-289-6367 for warranty and technical support of the Scaler

HP Illumination System, 5 Hole with NEW integrated LED bulb

HP Illumination System, 6 Pin



Factory option

1700-021-CS

Factory option

1700-019-CS

- Uses existing handpiece position
- Dual channel outputs individually adjustable from 2.5VDC to 4.0VDC
- Accommodates one additional handpiece tubing
- Illumination delay (adjustable between 0-15 seconds) can be set individually for 2 HP
- Automatic activation of HP illumination w/the dental unit foot control
- Simple installation into any new or existing dental unit
- Indicator LEDs confirms correct operation of the control module
- Powers all brands of optic handpieces and air-motors
- ISO standard handpiece terminal with stainless steel nut
- Handpiece tubing, (4' to 5' installed)
- Power transformer, output 24VAC, 1.05 Amps
- Installation hardware and instructions

HP Illumination Tubing, 5 Hole with NEW integrated LED bulb



Factory option

0017-615-CS

- Used to replace power optic handpiece tubing
- ISO standard power optic handpiece terminal with stainless steel nut
- Handpiece tubing, (4' to 5' factory installed)

HP Illumination Tubing, 6 Pin



Factory option

0017-603-CS

- Used to replace power optic handpiece tubing
- ISO standard power optic handpiece terminal with stainless steel adapter and nut
- Handpiece tubing, (4' to 5' factory installed)

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

Ancillary Products

Chair Control, Touchpad, Unit-Mounted



Factory option 3914-998-CS

- 4 programmable preset positions

Additional Chair Control, Touchpad, Unit-Mounted (3900/3800)

Factory option 3914-999-CS

Assistant's Hygiene Instrumentation



Factory option 1408-CS

- Canister with 1/2" drain
- Mounting bracket
- Drain tubing
- Autoclavable standard HVE w/tubing, 5'
- Autoclavable saliva ejector w/tubing, 5'
- Standard vacuum holder
- Universal holder
- Installed on Standard IC and Value control heads (not available on Euro)

Ancillary Shelf



Factory option 0010-310

- Accommodates additional ancillary products
- Designed for easy access to installed devices
- Powder coated aluminum
- Installed on Standard IC control head only

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Forest

SIMPLIFIED



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

ORDERING EQUIPMENT MADE EASY

SIMPLIFIED

Too many options can sometimes be overwhelming. That’s why this select product offering has just two color choices. That will make upgrading a treatment room or matching existing equipment easy.



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Ordering Forest Equipment Made Easy

Simplified

Too many options can sometimes be overwhelming. That’s why this select product offering has just two color choices. You can make upgrading a treatment room or matching existing equipment easy.

Fixed Mount with Sidebox, External Umbilical

Choose Your Control Head



Fusion Control Head

Pnuematic 3 HP Syringe

Part Number 4295SE



Standard IC Control Head

Pnuematic 3 HP Syringe

Part Number 4485SE

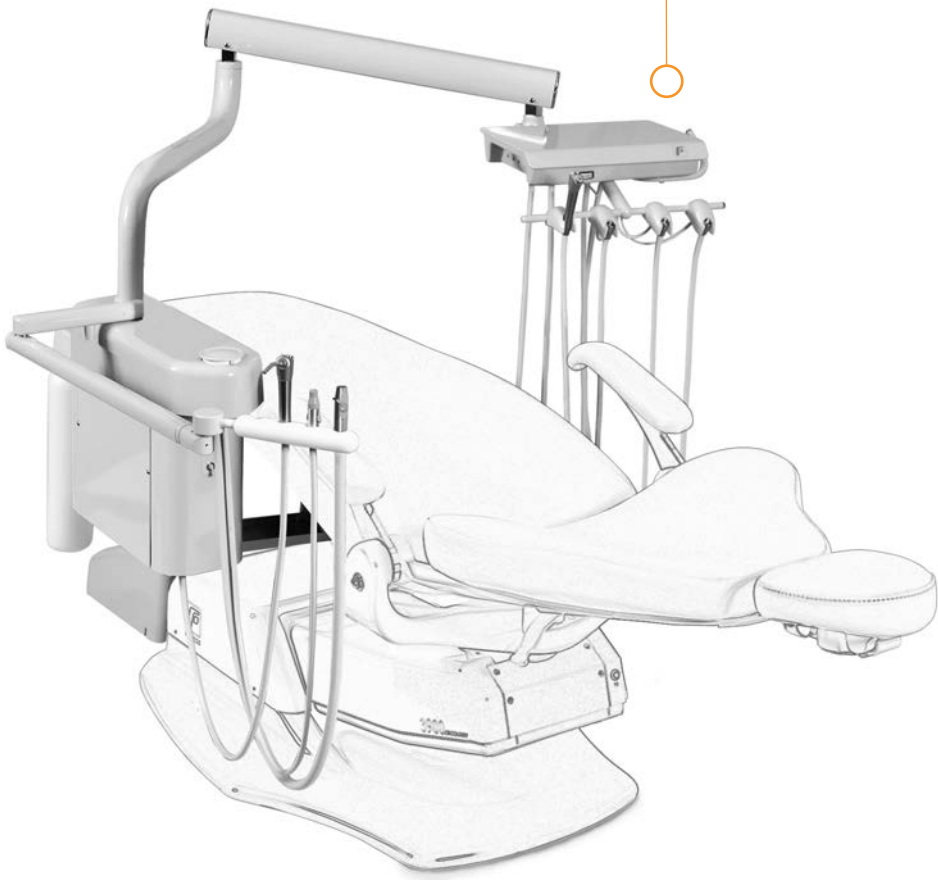


Euro Control Head

Pnuematic 4 HP Syringe

Part Number 4987SE

Delivery Systems can be mounted on any 2" fixed post mount brand of chair



Optional Equipment Choose One:

○ 1577 Cuspidor, Deluxe
(Sidebox comes with delivery)



○ 1679 Cuspidor, Basic
(Sidebox comes with delivery)



○ 9083 Light, Halogen
○ 1110-900 Light Post



○ 9078 Light, LED
○ 1110-900 Light Post



○ 9073 Light, LED
○ 1110-900 Light Post



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Fixed Mount without Sidebox, External Umbilical

Choose Your Control Head



Fusion Control Head

Pnuematic 3 HP Syringe

Part Number 4295FE



Standard IC Control Head

Pnuematic 3 HP Syringe

Part Number 4485FE



Euro Control Head

Pnuematic 4 HP Syringe

Part Number 4987FE

Delivery Systems can be mounted on any 2" fixed post mount brand of chair



Optional Equipment Choose One:

○ 1075FE Assistant's Telescoping Arm with Instrumentation



○ 1636FE Cuspidor, Basic



○ 1637FE Cuspidor, Basic with Assistant's Instrumentation



○ 9083 Light, Halogen



○ 1110-820 Light Post

○ 9078 Light, LED



○ 1110-820 Light Post

○ 9073 Light, LED



○ 1110-820 Light Post

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Standard Features	Factory Options	Part No	Pivot & Fixed Chair Mount Packages
Delivery System Choice of Control Head Automatic control for handpieces asepsis 3-way air/water syringe w/tubing, 5' (Euro tubing 3') handpiece tubing, 4-hole Midwest, asepsis style, 5' (Euro tubing 3') master control block master on/off individual handpiece pressure adjustment individual water coolant adjustment non-retraction water coolant air coolant adjustment handpiece flush handpiece exhaust canister/oil collector single handle brake release function self-contained clean water system w/ quick disconnect, single bottle, 1 1/2 liter wet/dry foot control Flexible control arm w/pneumatic brake vertical travel 28" / horizontal travel 62" External umbilical, 5' * Deluxe utility center air pressure regulator w/filter master air shut-off valve ! Chair adapter sold separately	DELIVERY SYSTEM Factory Installed STAR ETorque Electric System Handpiece illumination system One handpiece (5-hole), accommodates 1 add'l tubing Add'l handpiece illumination tubing (5-hole) One handpiece (6-pin), accommodates 1 add'l tubing Add'l handpiece illumination tubing (6-pin) ! Ancillaries use existing handpiece position (see page 71) Write in selection here _____ Write in selection here _____ ! Handpiece Bundle option (see pages 68-69) Write in selection here _____ Additional handpiece position, wet (4485) Additional handpiece position, air drive only (4485) Assistant's Hygiene Instrumentation (control head mount) Swivel tray, arm-mounted (over control head) Chrome accents Quick disconnect with flow control, 1/4" (specify air or water) Additional umbilical (maximum 9 feet) (per foot) Routing valve kit with water shutoff & regulator/factory installation Routing valve kit for use with cuspidor/factory installation 3800/3900 chair fixed 2" post mount chair adapter Δ non 3800/3900 2" post mount chair adapter, fixed (Only sold with unit) (specify chair make and model; non-cancellable/non-returnable)	MODEL CHOICE _____ * STAR-ETORQUE 1700-021-CS 0017-615-CS 1700-019-CS 0017-603-CS 1100-189-CS 1100-190-CS 1408-CS 1110-468 9000-998-CS 0016-455-CS 1115-999-CS 1115-315-CS 1115-316-CS 1115-195 1115-099-CA 1075FE 1400-113-CS 1400-113-CS 1110-820 1110-822 1300RM 1300TPM 1331-997	Pivot & Fixed Chair Mount Package Pivot & Fixed Chair Mount Package Ordering Guide Fusion Packages Fusion Package Ordering Guide Designer Friendly Colors & Surfaces Chairs & Stools Lights and Monitor Mounts Rear Mount Side Delivery Systems Ancillary Ordering & Options Control Head Accessories Simplified Cuspidors & Assistant's Instrumentation Simplified Rear & Side Delivery Systems Factory Options & Parts General Information Warranty
Assistant's Instrumentation 1075FE Central vacuum accessory package central vacuum canister w/disposable solids collector asepsis 3-way air/water syringe w/tubing, 5' autoclavable standard HVE lever valve w/tubing, 5' autoclavable saliva ejector lever valve w/tubing, 5' External umbilical, 5' Post mount, 2"	ASSISTANT'S INSTRUMENTATION Additional autoclavable standard HVE with tubing, 5' and autoclavable oral dry cup Hands-Free HVE, see page 91, Write in selection here _____ * CUSPIDOR, BASIC _____ with assistant's instrumentation Additional autoclavable standard HVE with tubing, 5' & autoclavable oral dry cup (1637FE) LIGHT, HALOGEN (mounting options, see page 45) OR LED OR LED •! LIGHT POST with 2" post clamp (specify light make/model) Light post, 35" (Forest) Light post, 35" (Pelton & Crane) MONITOR MOUNT ASSEMBLY (monitor screen not included) (see page 49) Ring mount (used with light post) Top post mount (used with monitor post) • Monitor post (for top post mount)	1075FE 1400-113-CS 1636FE 1637FE 1400-113-CS 9083 9078 9073 1110-820 1110-822 1300RM 1300TPM 1331-997	
Cuspidor, Basic 1636FE Gravity drain cuspidor w/swing away arm timed bowl rinse manual cup fill removable ceramic bowl Post mount adapter, 2" External 5' umbilical Water pressure regulator w/filter Master water shut-off valve	CUSPIDOR, BASIC Gravity drain cuspidor w/swing away arm timed bowl rinse manual cup fill removable ceramic bowl Post mount adapter, 2" External 5' umbilical Water pressure regulator w/filter Master water shut-off valve		
Cuspidor, Basic 1637FE Same as 1636FE above, w/additional features: Individual assistant's instrumentation holders Central vacuum accessory package central vacuum canister w/disposable solids collector asepsis 3-way air/water syringe w/tubing, 5' autoclavable standard HVE lever valve w/tubing, 5' autoclavable saliva ejector lever valve w/tubing, 5'	CUSPIDOR, BASIC Gravity drain cuspidor w/swing away arm timed bowl rinse manual cup fill removable ceramic bowl Post mount adapter, 2" External 5' umbilical Water pressure regulator w/filter Master water shut-off valve		
Equipment dimensions, see page 96 Shipping weight estimate: 4295FE & 4485FE Delivery Units 95 lb / 44 kg 4987FE Delivery Unit 100 lb / 46 kg	UNIT PLASTIC FINISHES: POWDER COAT PAINT FINISHES:	Choose only one post when sold with unit * When optional cuspidor added, standard utility center housing and cover is used to accommodate two umbilicals. ! Often overlooked Δ Lead time varies !	
	Simplified Surface Treatments		
	4295FE AVAILABLE IN FOREST GREY ONLY		
	UNIT PLASTIC FINISHES: POWDER COAT PAINT FINISHES:	Forest Grey White Forest Grey White	

Fixed Chair Post Mount, External Umbilical



Model 1636

Part Number 1636



Model 1637

Part Number 1637



Model 1075

Part Number 1075

Standard Features

Cuspidor, Basic 1636

- Gravity drain cuspidor w/swing away arm
- timed bowl rinse
- manual cup fill
- removable ceramic bowl
- Post mount adapter, 2"
- External umbilical, 5'

! Does not include post or chair adapter

Cuspidor, Basic 1637

- Same as 1636 above, w/additional features:
- Individual assistant’s instrumentation holders
- Central vacuum accessory package
 - central vacuum canister w/disposable solids collector
- 3-way air/water syringe w/tubing, 5'
- autoclavable standard HVE lever valve w/tubing, 5'
- autoclavable saliva ejector lever valve w/tubing, 5'

! Does not include post or chair adapter

Post Mount Vacuum Telescoping Arm 1075

- Post mount adapter, 2"
- Central vacuum accessory package
 - post mounted central vacuum canister w/disposable solids collector
- 3-way air/water syringe w/tubing, 5'
- autoclavable standard HVE lever valve w/tubing, 5'
- autoclavable saliva ejector lever valve w/tubing, 5'
- External umbilical, 8'

! Does not include post or chair adapter

Factory Options

CUSPIDOR, BASIC

Additional autoclavable standard HVE with tubing, 5' & autoclavable oral dry cup (1637) 1400-113-CS

Hands-Free HVE, see page 91, Write in selection here

Quick disconnect with flow control, 1/4" (specify air or water) 0016-455-CS

! Self-contained clean water system with quick disconnect 1700-310-CS

Standard utility center housing & cover 1115-020

► Standard utility center 1115-001

Additional umbilical (per foot) 1115-999-CS

! MOUNTING POSTS 2" diameter:

Control post with cap, 18" 1110-487

! LIGHT POSTS: 44" 2" diameter

Forest model 1110-794

Pelton & Crane 1110-785

ASSISTANT’S INSTRUMENTATION 1075

Additional standard autoclavable HVE with tubing, 5' & autoclavable oral dry cup 1400-113-CS

! Self-contained clean water system with quick disconnect 1700-323-CS

Assistant’s air shut-off & regulator utility center w/o housing & cover 1115-008

Standard utility center housing & cover 1115-020

Additional umbilical (per foot) 1115-998-CS

► Standard utility center includes: housing, cover, air & water pressure regulator with filter, master air and water shut-off valves.

! Often overlooked

⚠ **WARNING:** Isolated water sources (as opposed to city water hook-up) are highly recommended to allow full control of the source and quality of water circulating through the syringe and cuspidor water lines. If you do not have a water quality protocol, the optional self-contained water system allows you to isolate this product from the municipal water supply, however, it requires monitoring, maintenance and disinfection of the water lines, reservoir water bottle and source water.

Equipment dimensions, see page 96

Shipping weight estimate:

1636 17 lb / 8 kg

1637 20 lb / 9 kg

1075 15 lb / 7 kg

UNIT PLASTIC FINISHES:

☐ Forest Grey

☐ White

POWDER COAT PAINT FINISHES:

☐ Forest Grey

☐ White

Simplified Surface Treatments

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant’s Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Rear Beneath Shelf or Work Surface Mount



Model 1284 - Mounts under floating work surface

Shown with optional Clean Water System 1700-312 not included in price Assistant’s Utilities 1115-008 sold separately

Part Number 1284

Model 1286 - Mounts beneath shelf

Optional Clean Water System 1700-312 and Assistant’s Utilities 1115-008 sold separately

Part Number 1286

Standard Features	Factory Options	Part No
Horizontal Assistant’s Delivery System 1284/1286 Central vacuum accessory package, swing mount arm-mounted central vacuum canister w/disposable solids collector 3-way air/water syringe w/tubing, 5’ autoclavable standard HVE lever valve w/ tubing, 5’ autoclavable saliva ejector lever valve w/ tubing, 5’ Horizontal surface mount Link assembly, arm to plate (1286) External umbilical, 7’	DELIVERY SYSTEM Additional autoclavable standard HVE with tubing, 5’ & autoclavable oral dry cup Hands-free HVE, see page 91, Write in selection here Slide mount, 7 5/8” extension (1284) ! Standard utility center housing & cover ! Self-contained clean water system with quick disconnect, single bottle ! Assistant’s air shut-off & regulator utility center w/o housing & cover ! Often overlooked	MODEL CHOICE 1400-113-CS 1700-142 1115-020 1700-312 1115-008

⚠ WARNING: Isolated water sources (as opposed to city water hook-up) are highly recommended to allow full control of the source and quality of water circulating through the syringe and cuspidor water lines. If you do not have a water quality protocol, the optional self-contained water system allows you to isolate this product from the municipal water supply, however, it requires monitoring, maintenance and disinfection of the water lines, reservoir water bottle and source water.

! Simplified Surface Treatments

Mounting plate dimensions, see page 99
Shipping weight estimate:
1284 12 lb / 5 kg
1286 14 lb / 6 kg

- UNIT PLASTIC FINISHES:

☐ Forest Grey
- POWDER COAT PAINT FINISHES:

☐ Forest Grey

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant’s Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Simplified Rear and Side Delivery Systems

Side Cabinet & Wall Mount



Model 3235

12" Rigid Arm Cabinet Mount

Part Number 3235



Model 3230

22" Rigid Arm Wall Mount

Part Number 3230

Standard Features	Factory Options	Part No
Doctor Standard IC Delivery System 3235/3230 Standard IC control head automatic control for three handpieces 3-way air/water syringe w/tubing, 5' 4-hole handpiece tubing, 5' pneumatic master control block w/diaphragm master on/off individual handpiece pressure adjustment individual water coolant adjustment non-retraction water coolant air coolant adjustment handpiece flush handpiece exhaust canister/oil collector single handle brake release function self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter wet/dry foot control - Shadow powder coat standard Flexible control arm w/pneumatic brake (3235) cabinet vertical mount 12" rigid arm & 18" flex arm Flexible control arm w/pneumatic brake (3230) wall vertical mount 22" rigid arm & 18" flex arm External umbilical, 5' Standard utility center housing and cover (3230) Air pressure regulator w/filter Master air shut-off valve	DELIVERY SYSTEM Factory Installed STAR ETorque Electric System Handpiece Illumination System One handpiece (5-hole), accommodates 1 add'l tubing Add'l handpiece illumination tubing (5-hole) One handpiece (6-pin), accommodates 1 add'l tubing Add'l handpiece illumination tubing (6-pin) Assistant's hygiene instrumentation (control head mount) Quick disconnect with flow control, 1/4" (specify air or water) Swivel tray, arm-mounted (over control head) Standard utility center housing & cover (3235) Additional umbilical (per foot) Routing valve kit with water shutoff & regulator/factory installation Wall mount panel, 16" x 18"	STAR-ETORQUE 1700-021-CS 0017-615-CS 1700-019-CS 0017-603-CS 1408-CS 0016-455-CS 1110-468 1115-020 1115-999-CS 1115-315-CS 1110-586-CS
	! Often overlooked	

Simplified Surface Treatments

Mounting plate dimensions, see page 109
Shipping weight estimate:
3235, 3230 50 lb / 23 kg

- UNIT PLASTIC FINISHES:

☐ Forest Grey
- POWDER COAT PAINT FINISHES:

☐ Forest Grey

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Rear or Side Beneath Counter Mount



Model 2400

Beneath Counter Mount with swivel included
Countertop/work surface not included

Part Number 2400



Optional Beneath Counter Slide 1700-140



Optional Folding Arm Assembly 1110-115-CS
(Not available with STAR-ETORQUE option)

Standard Features

Doctor Standard IC Delivery System 2400

- Beneath counter swivel mount (work surface not included)
- Standard IC control head
 - automatic control for three handpieces
 - 3-way air/water syringe w/tubing, 5'
 - 4-hole handpiece tubing, 5'
 - pneumatic master control block w/diaphragm master on/off
 - individual handpiece pressure adjustment
 - individual water coolant adjustment
 - non-retraction water coolant
 - air coolant adjustment
 - handpiece flush
 - handpiece exhaust canister/oil collector
 - self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter
- wet/dry foot control - Shadow powder coat standard
- External braided umbilical, 7'
- Air pressure regulator w/ filter
- Master air shut-off valve

Factory Options

DELIVERY SYSTEM

- Factory Installed STAR ETorque Electric System
- Handpiece Illumination System
 - One handpiece (5-hole), accommodates 1 add'l tubing
 - Add'l handpiece illumination tubing (5-hole)
 - One handpiece (6-pin), accommodates 1 add'l tubing
 - Add'l handpiece illumination tubing (6-pin)
- Assistant's hygiene instrumentation (control head mount)
- Quick disconnect with flow control, 1/4" (specify air or water)
- Slide mount with swivel, 7 5/8" extension
- Folding arm assembly, 8" x 12" with wall mount
- Standard utility center housing & cover
- Routing valve kit with water shutoff & regulator/factory installation

MODEL CHOICE

Part No

- STAR-ETORQUE
- 1700-021-CS
- 0017-615-CS
- 1700-019-CS
- 0017-603-CS
- 1408-CS
- 0016-455-CS
- 1700-140
- 1110-115-CS
- 1115-020
- 1115-315-CS

! Often overlooked

Mounting plate dimensions, see page 109
Shipping weight estimate:
2400 28 lb / 13kg

- UNIT PLASTIC FINISHES: ☐ Forest Grey
- POWDER COAT PAINT FINISHES: ☐ Forest Grey

Simplified Surface Treatments

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Rear Beneath Counter Mount



Panel Mount is recessed in cabinet separate from duo articulating instrument delivery 2184BC.



Model 2184BC

Universal Heavy-Duty Mounting Plate 11" x 14 1/2" included (mounts to under counter)
Specify cabinet manufacturer

Part Number 2184BC



Model 2186BC

Universal Heavy-Duty Mounting Plate 11" x 14 1/2" included (mounts to under counter)
Specify cabinet manufacturer

Part Number 2186BC

Standard Features	Factory Options	Part No
Duo Delivery System 2184BC/2186BC Automatic panel control for three handpieces 3-way air/water syringe w/tubing, 5' 4-hole handpiece tubing, 5' pneumatic master control block w/diaphragm master on/off individual handpiece pressure adjustment individual water coolant adjustment non-retraction water coolant air coolant adjustment handpiece flush wet/dry foot control - Shadow powder coat standard Central vacuum accessory package, swing mount arm-mounted central vacuum canister w/disposable solids collector 3-way air/water syringe w/tubing, 5' autoclavable standard HVE lever valve w/tubing, 5' autoclavable saliva ejector lever valve w/tubing, 5' Articulating left/right instrument positioning Self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter Solid work surface, 19" x 24" ► Beneath counter rigid arm assembly w/heavy duty mounting plate, 11" x 14 1/2" External umbilical, 7' Air pressure regulator w/filter Master air shut-off valve	DELIVERY SYSTEM Factory Installed STAR ETorque Electric System Handpiece Illumination System One handpiece (5-hole), accommodates 1 add'l tubing Add'l handpiece illumination tubing (5-hole) One handpiece (6-pin), accommodates 1 add'l tubing Add'l handpiece illumination tubing (6-pin) Additional standard autoclavable HVE position with tubing, 5' & autoclavable oral dry cup Hands-Free HVE, see page 91, Write in selection here Quick disconnect with flow control 1/4" (for water) Quick disconnect with flow control 3/8" (for air) Storage organizer with bins & frame Swivel tray (requires purchase of 1110-179-CS) Swivel tray (mounts without storage organizer) Monitor mount assembly (monitor screen not included) (requires purchase of 1110-179-CS) Touchpad for 3900 chair (25 ft. Cable) Standard utility center housing & cover Routing valve kit with water shutoff & regulator/factory installation Standard laminate work surface in lieu of solid surface Special color laminate or solid work surface upgrade ! Often overlooked Δ Lead time varies	STAR-ETORQUE 1700-021-CS 0017-615-CS 1700-019-CS 0017-603-CS 1400-113-CS 0016-454-CS 0016-456-CS 1110-179-CS 1110-908-CS 1110-482-CS 1300-025-CS 1700-081-1-CS 1115-020 1115-315-CS 0023-003-SS-LAM 0023-003-X / 0023-003-SSV



► Beneath Counter Mount

Mounting plate dimensions, see pages 108
Panel cut-out, see page 108
Shipping weight estimate:
2184BC 85 lb / 36 kg
2186BC 85 lb / 36 kg

Simplified Surface Treatments

- UNIT PLASTIC FINISHES:
POWDER COAT PAINT FINISHES:
- ☐ Forest Grey
☐ Forest Grey
- SOLID WORK SURFACE:
- CORIAN®: ☐ Antarctica
- WILSONART®: ☐ Arctic Melange ☐ Baja Melange
- MEGANITE®: ☐ Cinder Granite ☐ Goldrush Gemstone ☐ Moonstone ☐ Platinum Granite
☐ Rocky Road Granite ☐ Silver Mist ☐ Starry Starry Night

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Simplified Rear and Side Delivery Systems

Panel Mounts



Model 2500

Part Number 2500

Optional Assistant's Hygiene Instrumentation 1408-CS

Model 2520

Part Number 2520

Standard Features	Factory Options	Part No
Doctor Control 2500 Automatic control for three handpieces 3-way air/water syringe w/tubing, 7' 4-hole handpiece tubing, 7' pneumatic master control block w/diaphragm master on/off individual handpiece pressure adjustment individual water coolant adjustment non-retraction water coolant air coolant adjustment self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter wet/dry foot control - Shadow powder coat standard Remote mount holder assembly w/control tubing, 6' External umbilical, 6' Air pressure regulator w/filter Master air shut-off valve	DELIVERY SYSTEM Handpiece Illumination System One handpiece (5-hole), accommodates 1 add'l tubing Add'l handpiece illumination tubing (5-hole) One handpiece (6-pin), accommodates 1 add'l tubing Add'l handpiece illumination tubing (6-pin) Standard utility center housing & cover Assistant's hygiene instrumentation with extended 10 1/2" holder bar (Drain Tubing 12') standard autoclavable HVE, saliva ejector with 2 holders and canister ! Quick disconnect for air or water - dealer installed (see page 94) ! Routing valve kit for city water - dealer installed (see page 94) ! Often overlooked	MODEL CHOICE 1700-021-CS 0017-615-CS 1700-019-CS 0017-603-CS 1115-020 1408-CS
Doctor Control 2520 Automatic control for two handpieces 3-way air/water syringe w/tubing, 7' 4-hole handpiece tubing, 7' pneumatic master control block w/diaphragm master on/off individual handpiece pressure adjustment non-retraction water coolant air coolant adjustment self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter wet/dry foot control - Shadow powder coat standard Remote mount holder assembly w/control tubing, 6' External umbilical, 6' Air pressure regulator w/filter Master air shut-off valve		

Wall/cabinet cut-out dimensions, see page 107
Shipping weight estimate:
2500 18 lb / 8 kg
2520 15 lb / 7 kg

! Simplified Surface Treatments

UNIT PLASTIC FINISHES:

☐ Forest Grey

POWDER COAT PAINT FINISHES:

☐ Forest Grey

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Rear or Side Mount Compact Controls



Model MT600

Work surface not included

Part Number MT600

Model MT610

Work surface not included

Part Number MT610

Standard Features		Factory Options	Part No
Doctor Control with 3 HP		DELIVERY SYSTEM	MODEL CHOICE
Top Shelf Mount	MT600TS	Handpiece Illumination System	
Beneath Shelf Mount w/Swivel	MT600BS	One handpiece (5-hole), accommodates 1 add'l tubing	1700-021-CS
Doctor Control with 2 HP		Add'l handpiece illumination tubing (5-hole)	0017-615-CS
Top Shelf Mount	MT610TS	One handpiece (6-pin), accommodates 1 add'l tubing	1700-019-CS
Beneath Shelf Mount w/Swivel	MT610BS	Add'l handpiece illumination tubing (6-pin)	0017-603-CS
All Units Include the Following Features:		Slide mount with swivel, 7 5/8" extension (order BS)	1700-140
Automatic control for three handpieces (MT600), or		Folding arm assembly, 8" x 12" with wall mount (order TS)	1110-115-CS
Automatic control for two handpieces (MT610)		Standard utility center housing & cover	1115-020
3-way air/water syringe w/tubing, 4'		Assistant's hygiene instrumentation with extended 10 1/2" holder bar (Drain Tubing 12') standard autoclavable HVE, saliva ejector with 2 holders & canister	1408-CS
4-hole handpiece tubing, 4'		! Quick disconnect for air or water - dealer installed (see page 94)	
pneumatic master control block w/diaphragm master on/off		! Routing valve kit for city water - dealer installed (see page 94)	
individual handpiece pressure adjustment			
non-retraction water coolant			
air coolant adjustment			
self-contained clean water system w/quick disconnect, single bottle, 1 1/2 liter		! Often overlooked	
wet/dry foot control - Shadow powder coat standard			
External umbilical, 12'			
Air pressure regulator w/filter			
Master air shut-off valve			

Mounting plate dimensions, see page 109
Shipping weight estimate:
MT600, MT610 20 lb / 10 kg

!

Simplified Surface Treatments

UNIT PLASTIC FINISHES:

☐ Forest Grey

POWDER COAT PAINT FINISHES:

☐ Forest Grey

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Simplified Rear and Side Delivery Systems

Mobile Carts



Model 5350

Part Number 5350



Shown with optional swivel tray,
arm mounted 1110-731



Model 5355

Part Number 5355

Standard Features

Doctor Standard IC Cart 5350

- Standard IC Control Head
 - automatic control for three handpieces
 - 3-way air/water syringe w/tubing, 4'
 - 4-hole handpiece tubing, 4'
 - pneumatic master control block w/diaphragm master on/off
 - individual handpiece pressure adjustment
 - individual water coolant adjustment
 - non-retraction water coolant
 - air coolant adjustment
 - handpiece flush
 - handpiece exhaust canister/oil collector
 - self-contained clean water system w/quick disconnect, single bottle 1 1/2 liter
 - wet/dry foot control - Shadow powder coat standard
- Adjustable height frame
- Twin wheel hard floor casters
- External 7' umbilical
- Standard utility center housing and cover
- Air pressure regulator w/filter
- Master air shut-off valve

Duo Standard IC Cart 5355

- Same as above, w/additional features:
- Central vacuum canister w/disposable solids collector
- autoclavable standard HVE lever valve w/tubing, 5'
- autoclavable saliva ejector lever valve w/tubing, 5'

Factory Options

DELIVERY SYSTEM	MODEL CHOICE	
Factory Installed STAR ETorque Electric System		STAR-ETORQUE
Handpiece Illumination System		
One handpiece (5-hole), accommodates 1 add'l tubing		1700-021-CS
Add'l handpiece illumination tubing (5-hole)		0017-615-CS
One handpiece (6-pin), accommodates 1 add'l tubing		1700-019-CS
Add'l handpiece illumination tubing (6-pin)		0017-603-CS
Additional standard autoclavable HVE with tubing, 5' & autoclavable oral dry cup (5355)		1400-113-CS
Hands-Free HVE, see page 91, Write in selection here		
Quick disconnect with flow control, 1/4" (specify air or water)		0016-455-CS
Swivel tray, arm-mounted		1110-731
Foot control support shelf, clamps to 2" base		0023-665
Additional umbilical (per foot)		1115-999-CS
Routing valve kit with water shutoff & regulator/factory installation		1115-315-CS

! Often overlooked

Equipment dimensions, see page 110
Shipping weight estimate:
5350 50 lb / 23 kg
5355 55 lb / 25 kg

! Simplified Surface Treatments

- UNIT PLASTIC FINISHES:

☐ Forest Grey
- POWDER COAT PAINT FINISHES:

☐ Forest Grey

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Standard Chair Adapters - for 2” Post

3900/3800 fixed accommodates external umbilical 1115-195

Monitor Posts - 2” Diameter with Oilite Bushing

(used with 1300TPM)

30” with dual slot, fixed sidebox mount 1331-998
30” with slot & post clamp, fixed mount 1331-997

Light Posts - 2” with Light Bushing

35” fixed (includes post clamp)
Forest 1110-820
Pelton & Crane 1110-822
44” fixed
Forest 1110-794
Pelton & Crane 1110-785

44” fixed with slot (for monitor cord)
Forest 1110-050

44” fixed with slot (sidebox access)
Forest 1110-900
Pelton & Crane 1110-890

44” fixed with dual slot (for monitor cord & sidebox access)
Forest 1110-170

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant’s Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty



250 Watt (Universal Power Supply)

Used when two or more ancillary devices are ordered with a Forest
Operator Light (excluding the 3900 soft-touch massager power source)

- Provides power for MULTIPLE electrical accessories such as: fiber-optics operator lights, scalers, curing lights, and electric micro motors
- Size: 7.5" - 5.5" - 2.5". Easily fits in all Forest floor boxes
- Quick disconnect connectors for all of Forest ancillary products
- Available in both 230 volt and 110 volt
- Auto-resettable Klaxon switch: inbuilt over current and over temperature protection, eliminating the need for manually changing fuses or circuit breakers

Please note: built in air-electric switch, power source turns on/off w/pilot air only



Handpiece Illumination System,
5 hole with NEW integrated LED bulb

Factory option 1700-021-CS

Handpiece Illumination System, 6 Pin

Factory option 1700-019-CS

- Dual channel outputs individually adjustable from 2.5VDC to 4.0VDC
- Accommodates one additional handpiece tubing
- Illumination delay (adjustable between 0-15 seconds) can be set individually for two handpieces
- Automatic activation of handpiece illumination w/the dental unit foot control
- Simple installation into any new or existing dental unit
- Indicator LED's confirms correct operation of the control module
- Powers all brands of optic handpieces and air-motors
- ISO standard handpiece terminal w/stainless steel nut
- Handpiece tubing, 7' long (4' to 5' installed)
- Power transformer, output 24VAC, 1.05 Amps
- Installation hardware and instructions



Handpiece Illumination Tubing,
5 hole with NEW integrated LED bulb

Factory option 0017-615-CS

Used to replace power optic handpiece tubing

- ISO standard power optic handpiece terminal w/stainless steel nut
- Handpiece tubing, 7' (4' to 7' factory installed)



Handpiece Illumination Tubing, 6 pin

Factory option 0017-603-CS

Used to replace power optic handpiece tubing

- ISO standard power optic handpiece terminal w/stainless steel adapter and nut
- Handpiece tubing, 7' (4' to 7' factory installed)



Standard Deluxe Air/Water Syringe with Quick-release Tip

Factory option	1109-010
with Sheathed Tubing Assembly, 5'	
Sterling	1107-755
Shadow	1107-755-SH
Black	1107-755-BLK

Straight deluxe Air/Water Syringe with Quick-release Tip

Factory option	1109-025
with Sheathed Tubing Assembly, 4'	
Sterling	1109-585-STG
Shadow	1109-585-SHW
Black	1109-585-BLK

Standard Deluxe One Button Syringe with Quick-release Tip

Factory option	1109-035
with Sheathed Tubing Assembly, 5'	
Sterling	1109-485-STG
Shadow	1109-485-SHW
Black	1109-485-BLK




Additional HVE with Tubing 5' and Autoclavable Oral Dry Cup

Factory option	1400-113-CS
• Standard vacuum valve	
Oral dry cup only	1400-010



Hands-Free HVE Chair-Mounted Flexible HVE Holder w/ HVE and 6' Tubing

Factory options	
White Holder Arm w/Sterling HVE Tubing	1400-279-STG
White Holder Arm w/Grey HVE Tubing	1400-279-GRY
Black Holder Arm w/Shadow HVE Tubing	1400-279-SHW
Black Holder Arm w/Black HVE Tubing	1400-279-BLK
Oral dry cup (pictured above)	1400-010

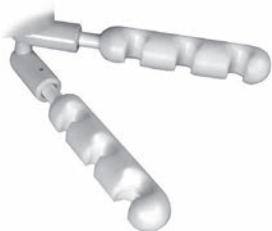

Dual Articulating Quad Holder Assembly

Factory option	0024-227-CS
4 standard, 4 universal, White	0024-227
4 standard, 4 universal, Forest Grey	0024-227-G
4 standard, 4 universal, Designer Friendly® (specify color)	0024-227-C
8 universal, White	0024-224
8 universal, Forest Grey	0024-224-G
8 universal, Designer Friendly (specify color)	0024-224-C



Support Shelf

Factory option	0023-665
Forest Grey	0023-665-G

- Provides adequate support for use of foot control while allowing cart mobility.

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty



Contoured Utility Center Housing and Cover

Factory option	3900-101
Forest Grey	3900-101-G
Shadow	3900-101-SH
Designer Friendly® (specify color)	3900-101-C

Complete Contoured Utility Center

Factory option	3900-102
Forest Grey	3900-102-G
Shadow	3900-102-SH
Designer Friendly (specify color)	3900-102-C
<ul style="list-style-type: none">Air and water pressure regulator w/filterMaster air and water shut-off valves	



Deluxe Utility Center Housing and Cover

Forest Grey	1115-007-G
Shadow	1115-007-SH
Black	1115-007-BLK
White	1115-007-WH
<ul style="list-style-type: none">Rear umbilical access	
Forest Grey	1115-009-G
Shadow	1115-009-SH
Black	1115-009-BLK
White	1115-009-WH
<ul style="list-style-type: none">Right umbilical access	

Complete Deluxe Utility Center

Forest Grey	1115-005-G
Shadow	1115-005-SH
Black	1115-005-BLK
White	1115-005-WH
<ul style="list-style-type: none">Rear umbilical accessAir and water pressure regulator w/filterMaster air and water shut-off valves	



Standard Utility Center Housing and Cover

Factory option	1115-020
Forest Grey	1115-020-G
Designer Friendly (specify color)	1115-020-C

Complete Standard Utility Center

Factory option	1115-001
Forest Grey	1115-001-G
Designer Friendly (specify color)	1115-001-C
<ul style="list-style-type: none">Air and water pressure regulator w/ filterMaster air and water shut-off valves	



Recessed utility center housing and cover

1115-305

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty



**Tray Holder Assembly, Use with 1095/1295
Rear Pivot Assistant's Instrumentation**

Factory option	1110-680
Forest Grey	1110-680-G
Designer Friendly® (specify color)	1110-680-C
Includes:	
• Pivoting powder coated tray holder, 9 3/4" x 13 1/2"	
• Clamp assembly adapts to 1 1/2" horizontal or vertical post	



**Slide with Swivel 7 5/8" Extension,
use w/1284, 2184, 2400, MT's**

Factory option (2400/MT's)	1700-140
Forest Grey (2400/MT's)	1700-140-G
Factory option (1284, 2184) (Does not include swivel bracket)	1700-142



Folding Arm Assembly, use with 2400, MT's

Factory option	1110-115-CS
Forest Grey	1110-115-G
Includes:	
• 8" folding arm with 1/2" pin and umbilical retainer	
• 4" x 4" wall mounting plate which accepts 1/2" pin	
• 12" folding arm with 1/2" pin and umbilical retainer	
• Control mount with 1/2" pin	



Universal Heavy Duty Mounting Plate

10" x 24" beneath counter	0010-084-BLK
11" x 14 1/2" beneath counter	0010-086-BLK
10" x 24" in lieu of 11" x 14 1/2"	0010-084-CS



26" Rigid with 24" Flexible Arm In Lieu of Top Post Mount

Factory option	1110-771-RM
----------------	-------------

22" Rigid with 18" Flexible Arm In Lieu of Top Post Mount

Factory option	1110-536-RM
----------------	-------------

26" Rigid With 18" Flexible Arm In Lieu Of Top Post Mount

Factory option	1110-537-RM
• Designed for use with 2" diameter post.	



Floating Work Surface Arm Assembly, 20"

Forest Grey	0023-815-G
Designer Friendly (specify color)	0023-815-C



Flow Control

Quick Disconnects

1/4" Female and Male QD with Flow Control	
Factory Installed	0016-454-CS
1/4" Female QD with flow control, black/grey/shadow/white	
Dealer Installed	0016-454-BLK /-G /-SH /-WH
1/4" Female and Male QD With Flow Control, Arm and Control Head Mount	
Factory Installed	0016-455-CS
1/4" Female QD with Flow Control, Arm and Control Head Mount, Black /Grey	
Dealer Installed	0016-455-BLK /-G
3/8" Female and Male QD With Flow Control, Factory Installed	
Factory Installed	0016-456-CS
3/8" Female QD with Flow Control, Black /Grey/Shadow/White	
Dealer Installed	0016-456-BLK /-G /-SH /-WH



Water Shut-Off
and Regulator

Routing Valve Kits

— Check with Local Authorities About Installation Requirements	
Factory installed with water shut-off and regulator	1115-315-CS
Factory installed without water shut-off and regulator	1115-316-CS
Dealer Installed Includes:	
Toggle Valve/Bracket/Tubing/Water Shut-Off/Regulator	
black / grey / shadow / white	1115-310-BLK / -G / -SH / -WH
Dealer Installed Includes: Toggle Valve, Bracket, And Tubing	
black / grey / shadow / white	1115-312-BLK / -G / -SH / -WH







High Density Bottle

1 1/2" liter with quick disconnect adapter	1100-431
1 1/2" liter with quick disconnect upgrade kit	1100-423
For use with clean water system	
• Opaque bottle	
• High Density Polyethylene (HDPE)	



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

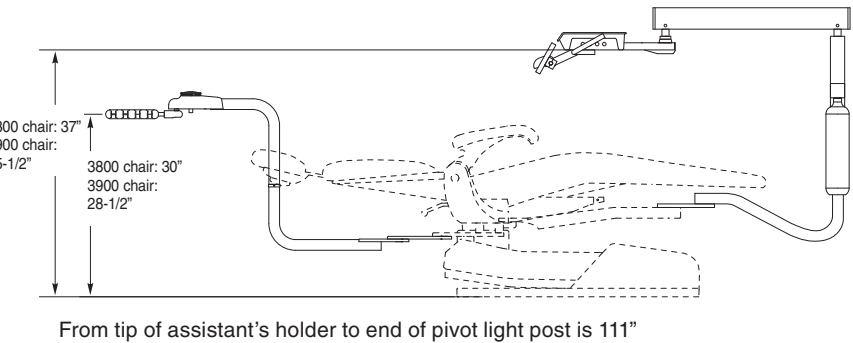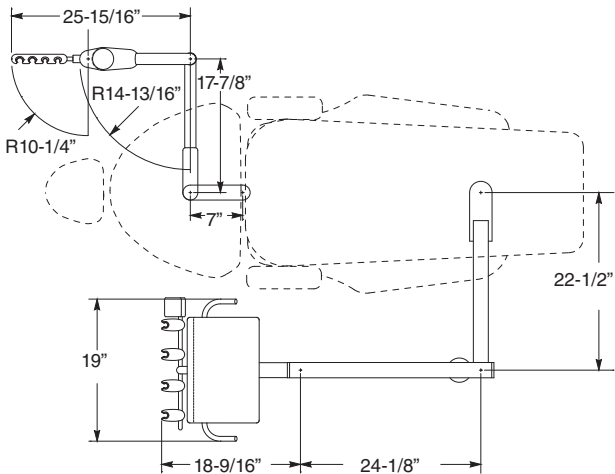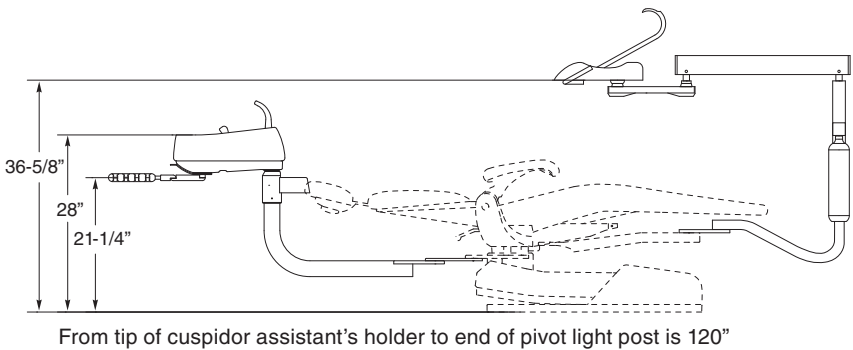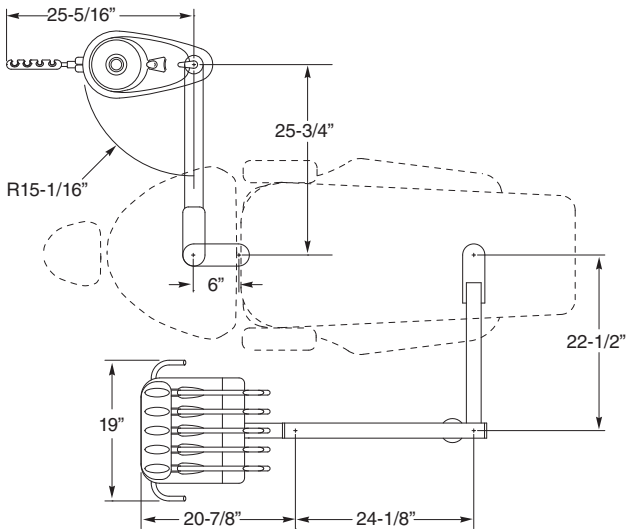Pivot Delivery Systems




Fixed with Sidebox Delivery Systems

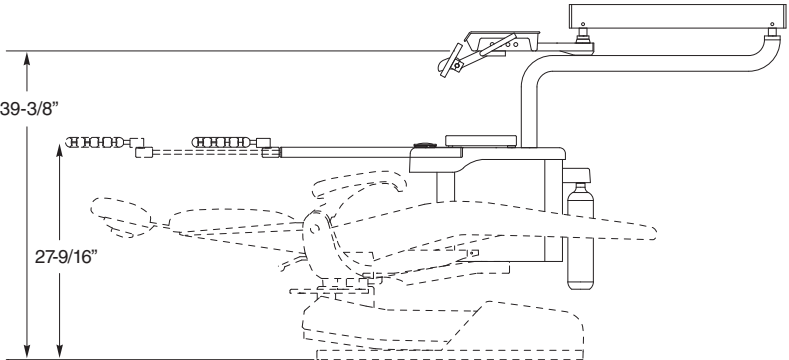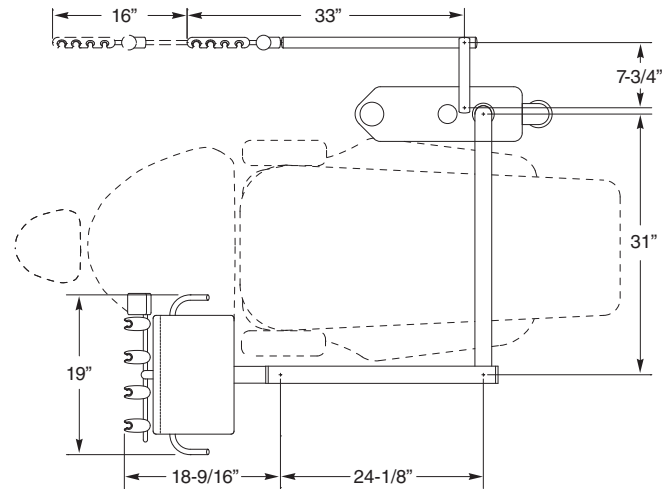


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

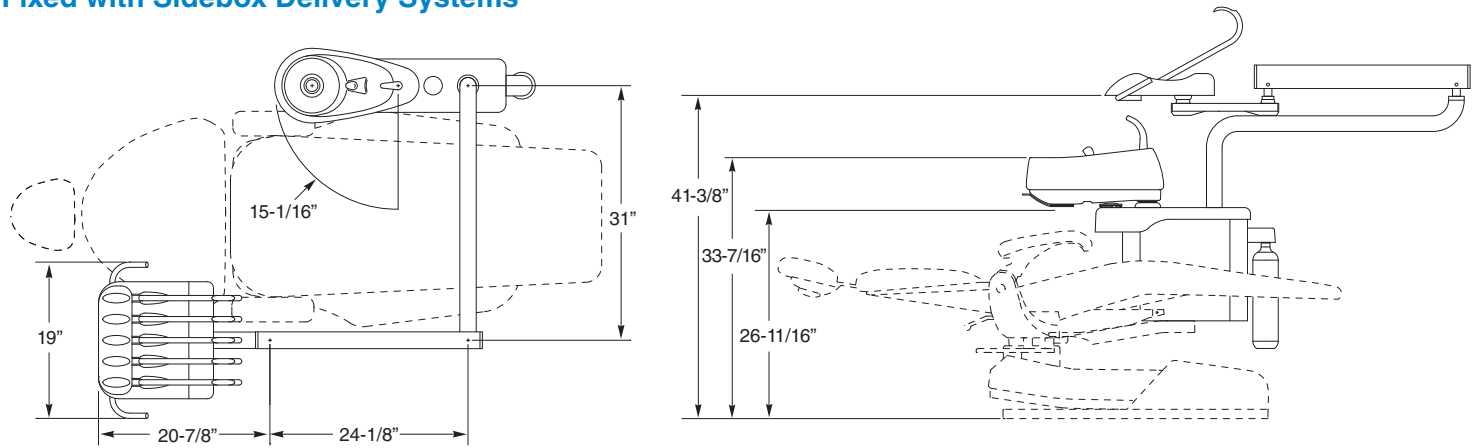
Fixed with Sidebox Delivery Systems




Fixed without Sidebox Delivery Systems

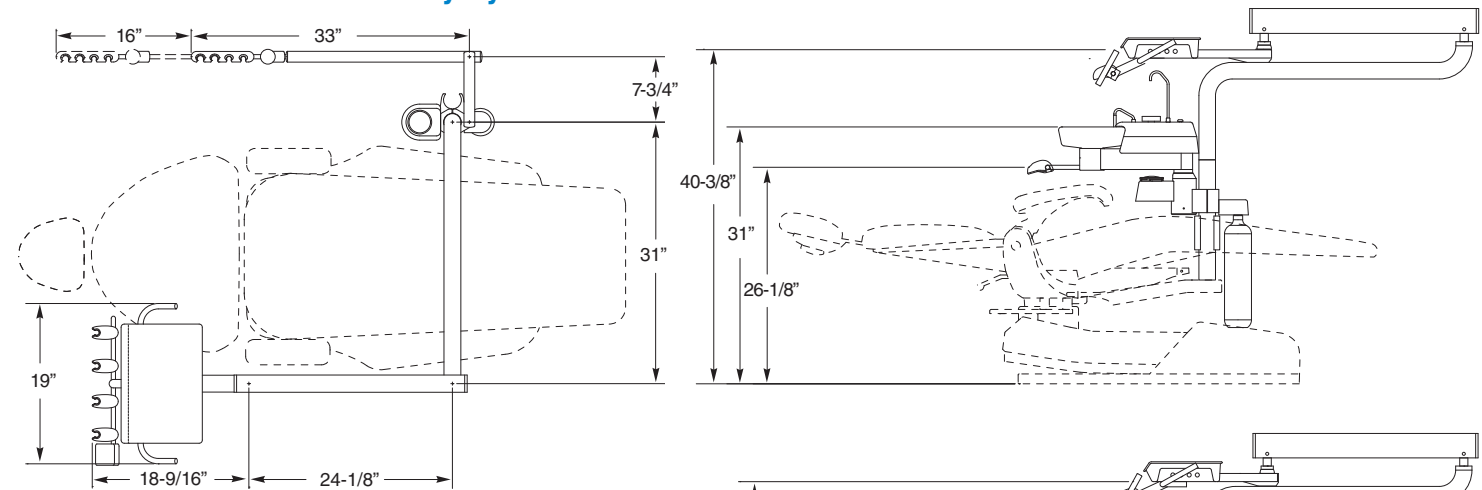

1636/1637 Basic Cuspidor

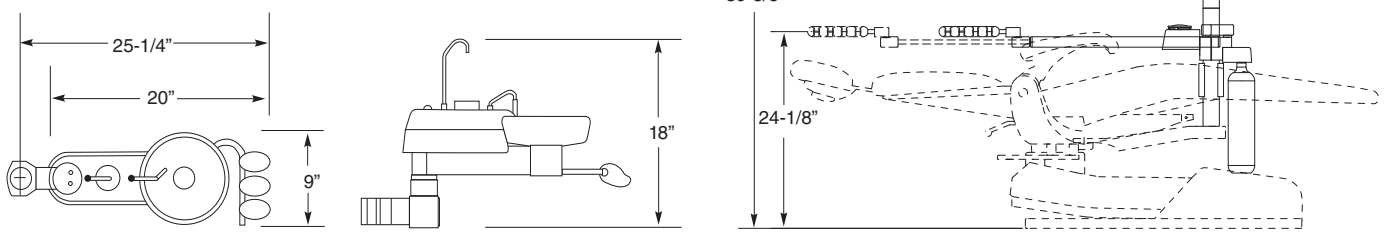

3800 Chair

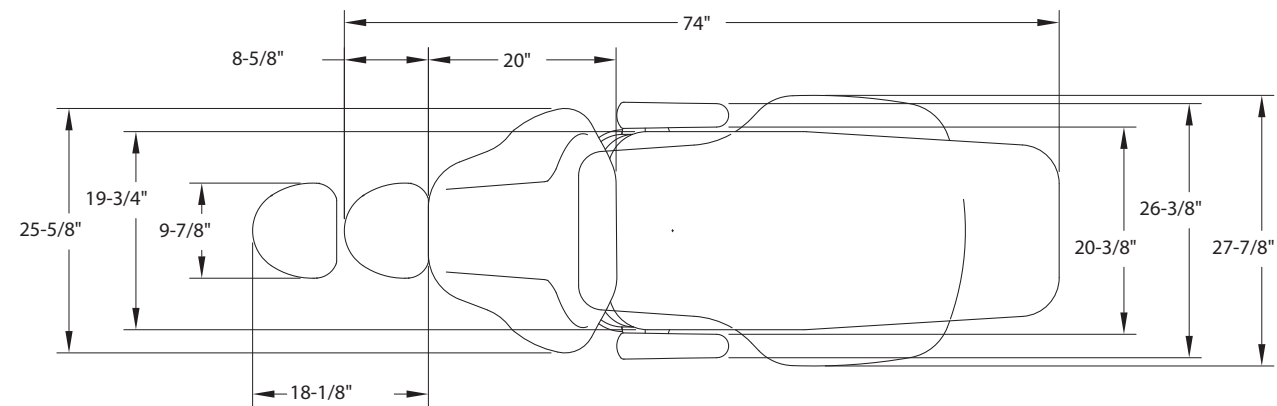

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

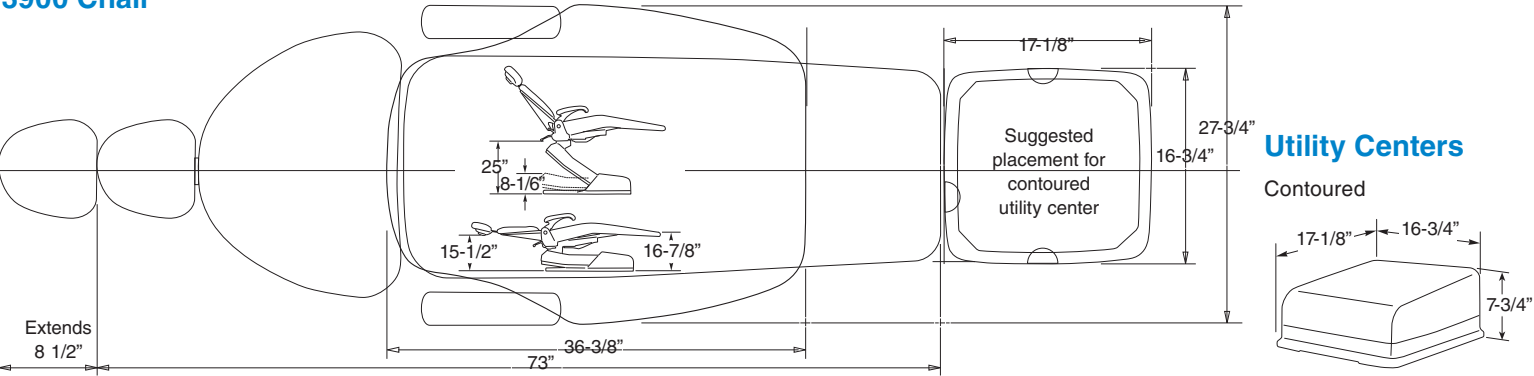
General
Information

Warranty


General Information

Dimensions


3900 Chair




6176 Doctor's Comfort Stool

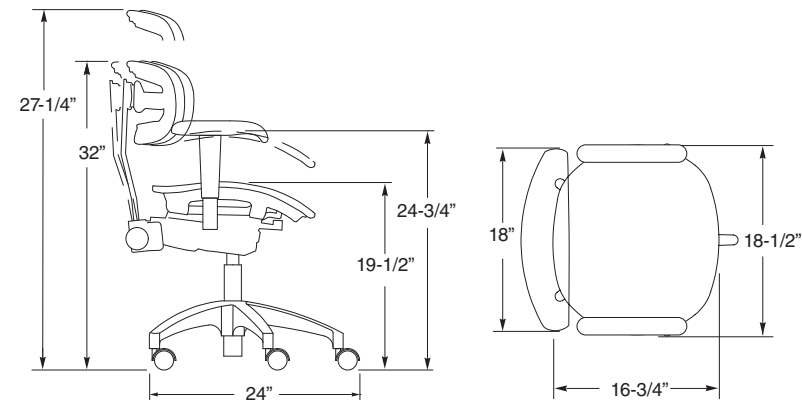

6159C Doctor's Deluxe Stool

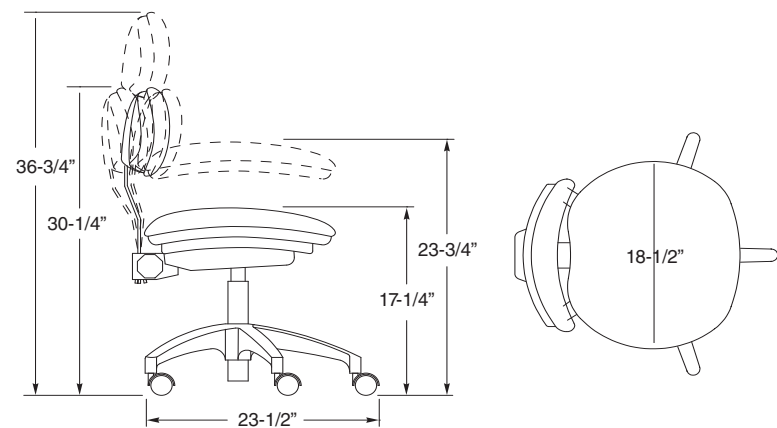


6156C Doctor's Economy Stool

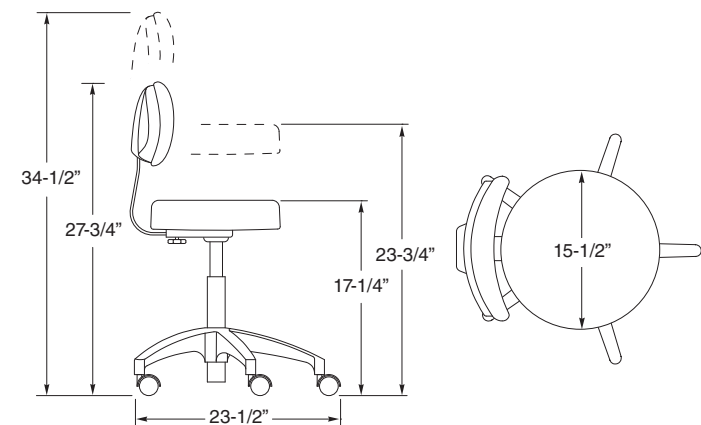


Deluxe

Standard

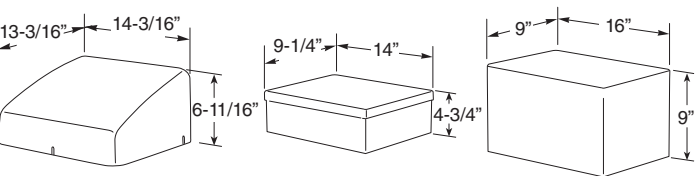
Recessed




6187 Assistant's Comfort Stool

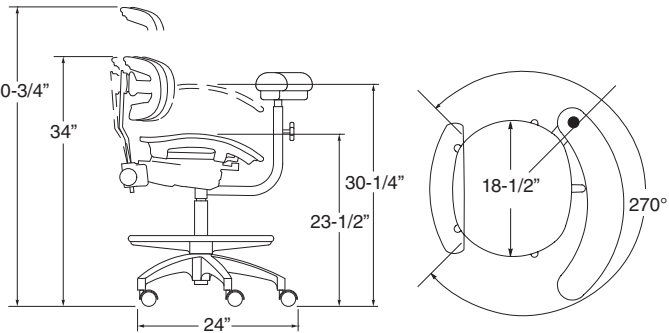

6160C/6161C Assistant's Deluxe Stool

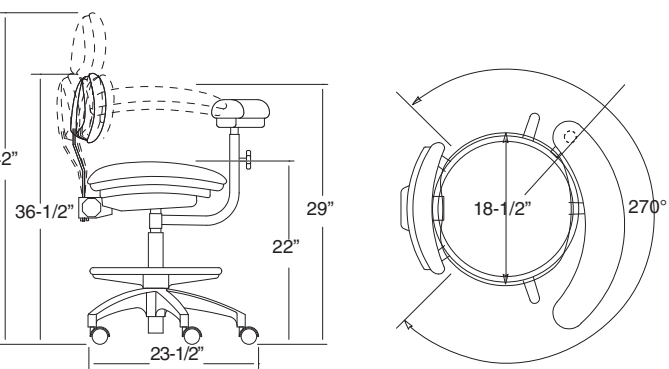

6157C Assistant's Economy Stool

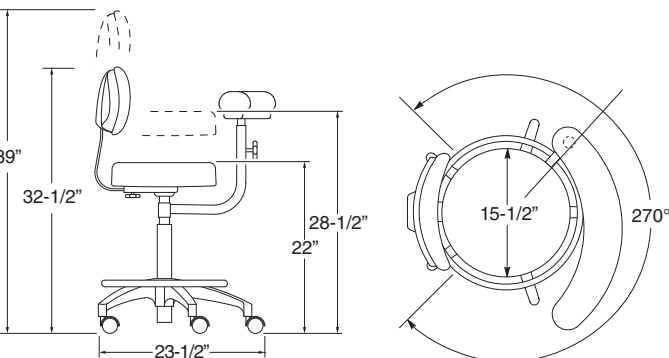

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

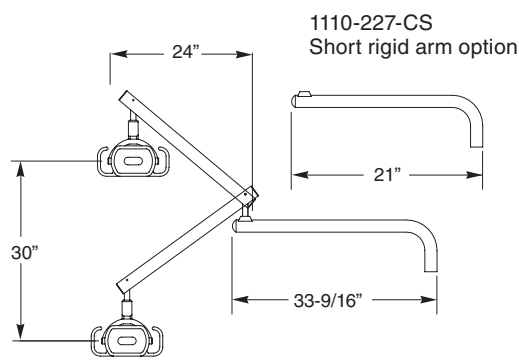
Factory
Options & Parts

General
Information


Warranty

Dimensions


Chair-Mounted Light



Ceiling-Mounted Light

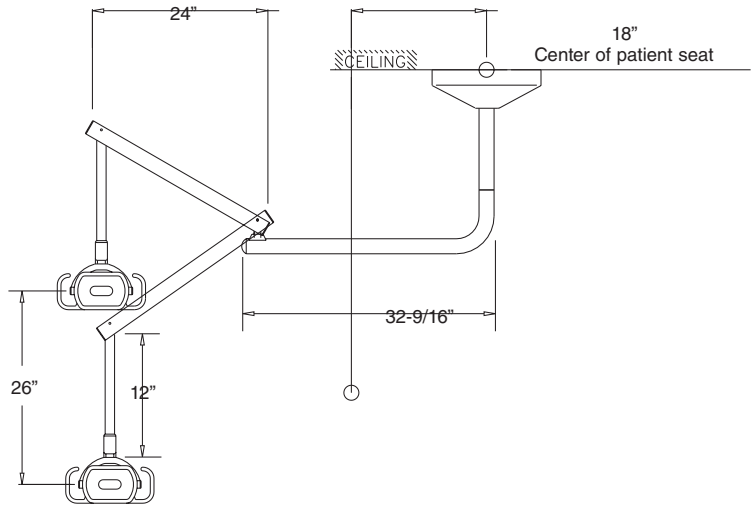


Track Mount Light/Monitor

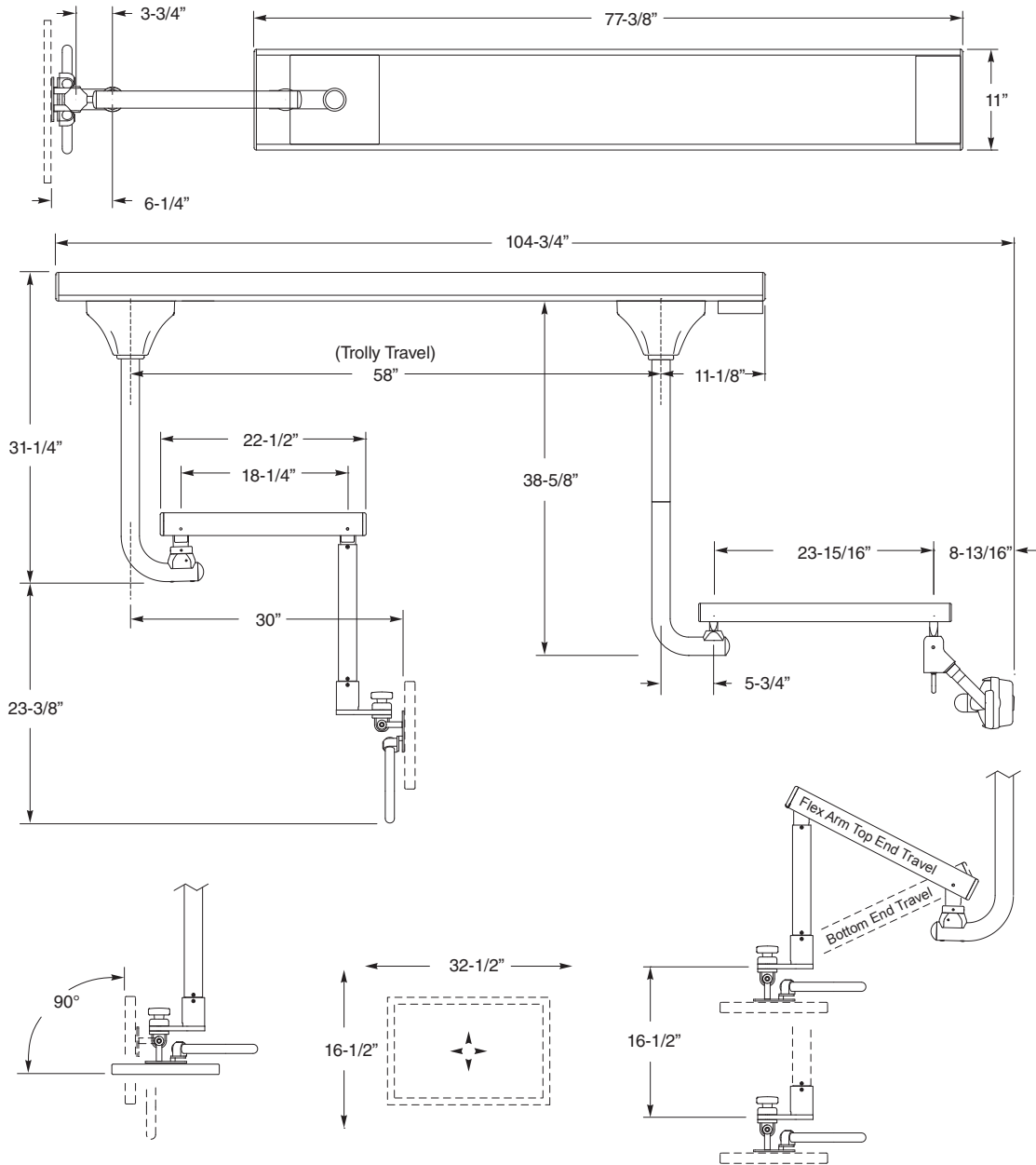


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

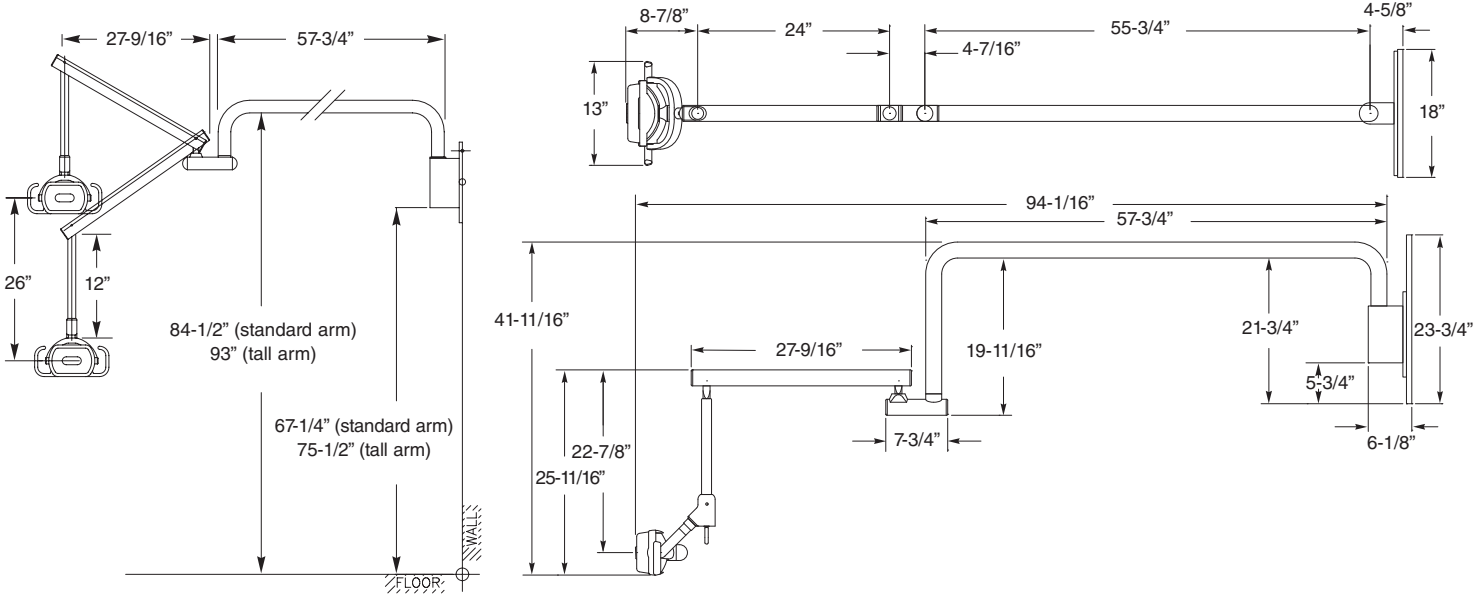
Standard & Tall Wall & Cabinet Mounted Light Tall Wall & Cabinet Mounted Light




Wall & Cabinet Mount

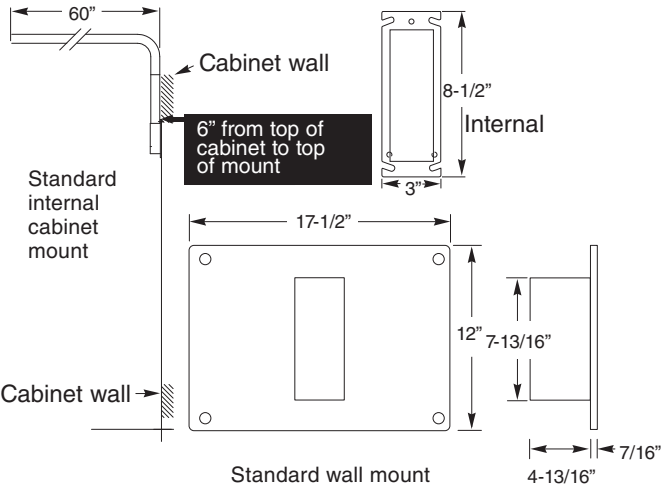


Optional Cabinet Mount

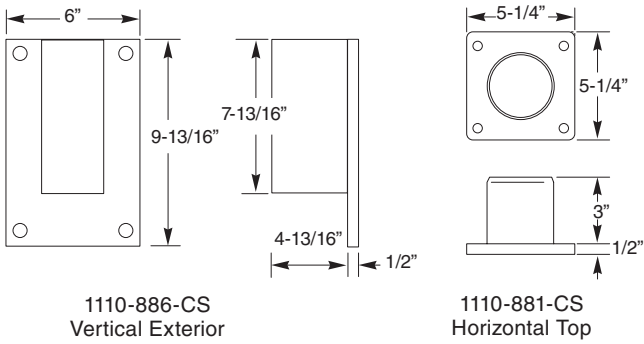

1300RM

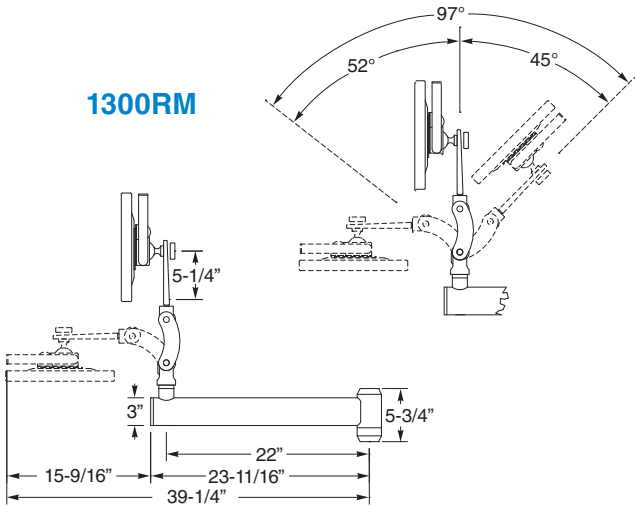



1300TPM

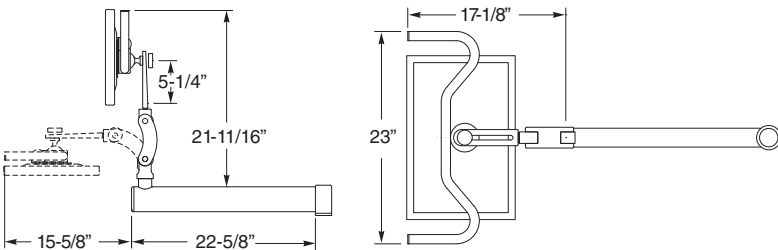



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

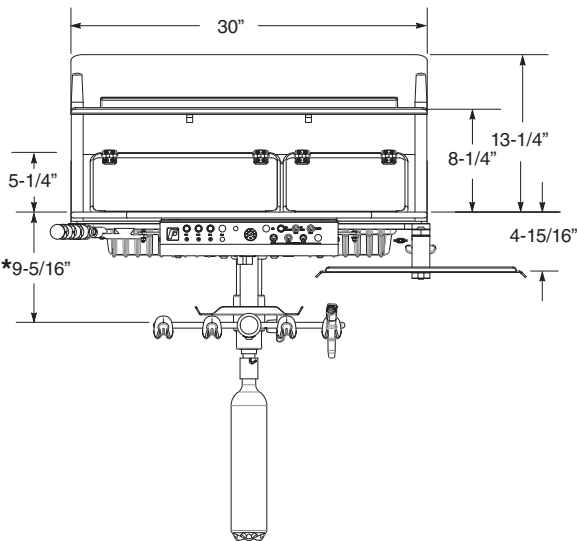
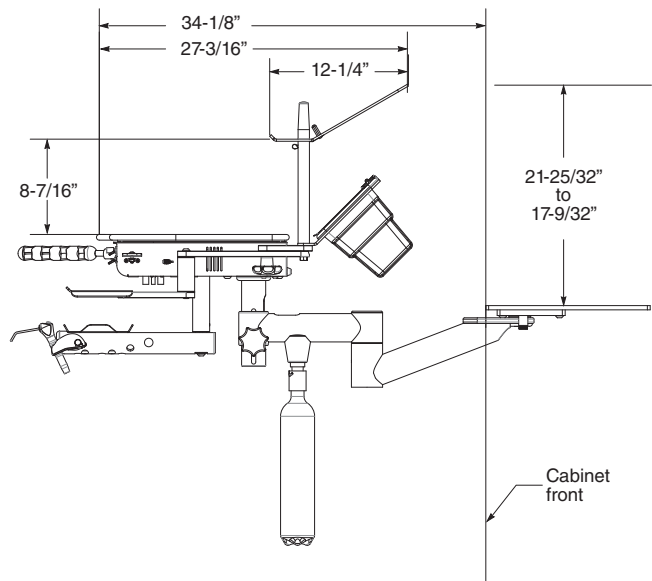
7020PRO-BC Beneath Counter




7021PRO-BC Beneath Counter

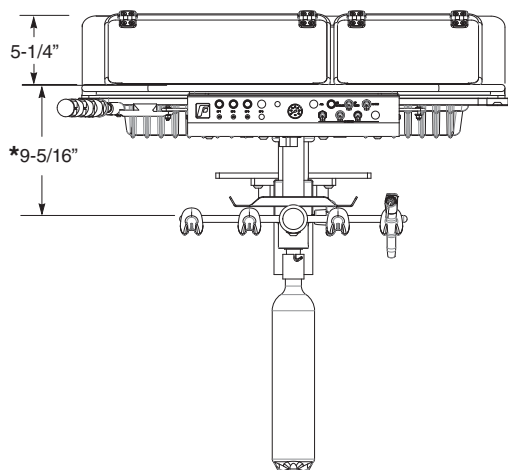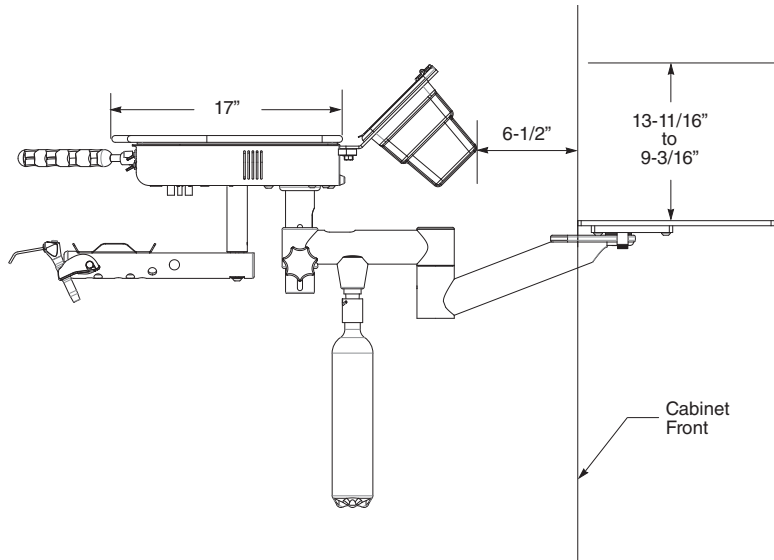

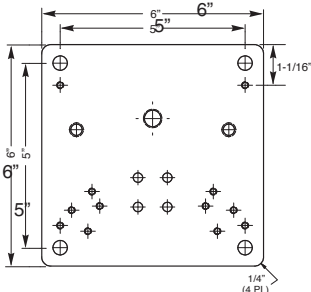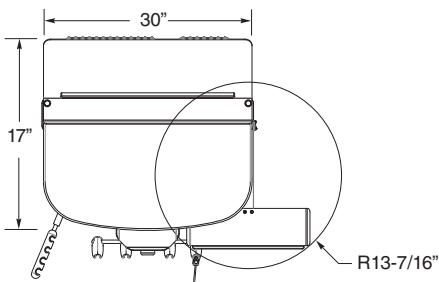- Large contoured work surface standard for 7020/7021
- Work surface 9/16" thick
- Shown with organization features
- Height adjustable work surface, 5"



0010-086

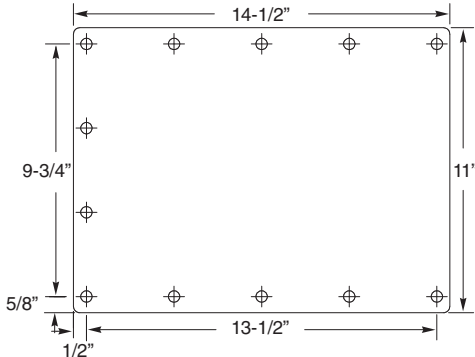

7020/7021PRO-BC Mounting Plates

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

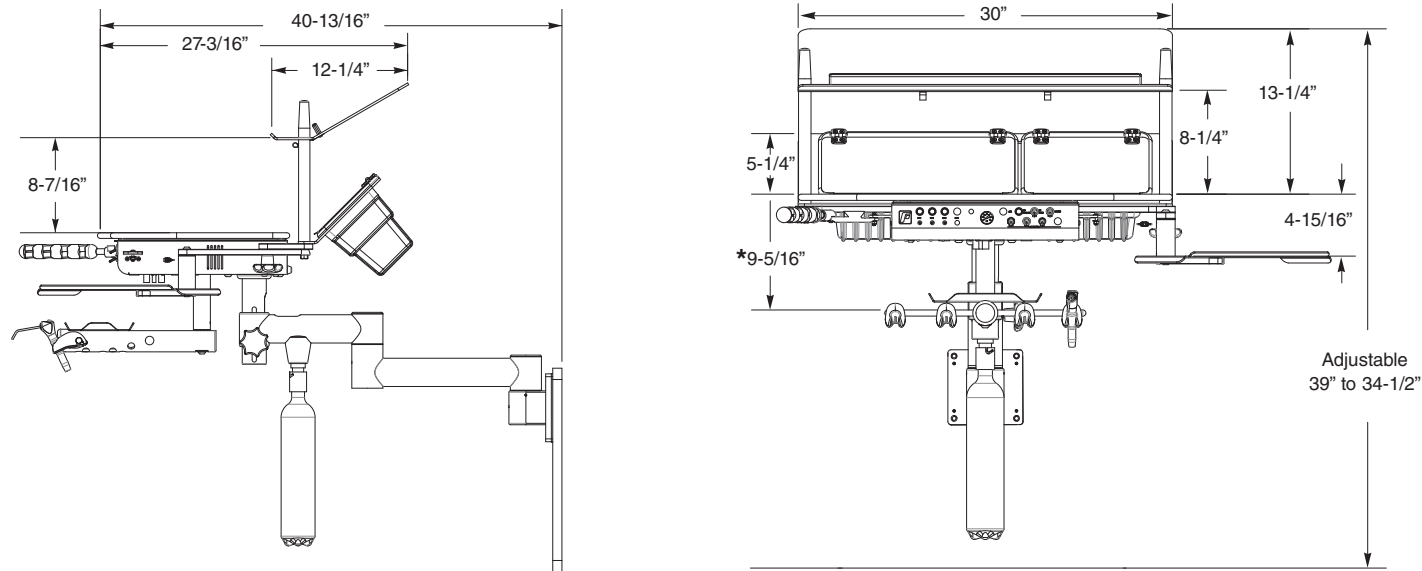
Factory
Options & Parts

General
Information


Warranty

Dimensions


7020PRO-WM Wall Mount




7021PRO-WM Wall Mount

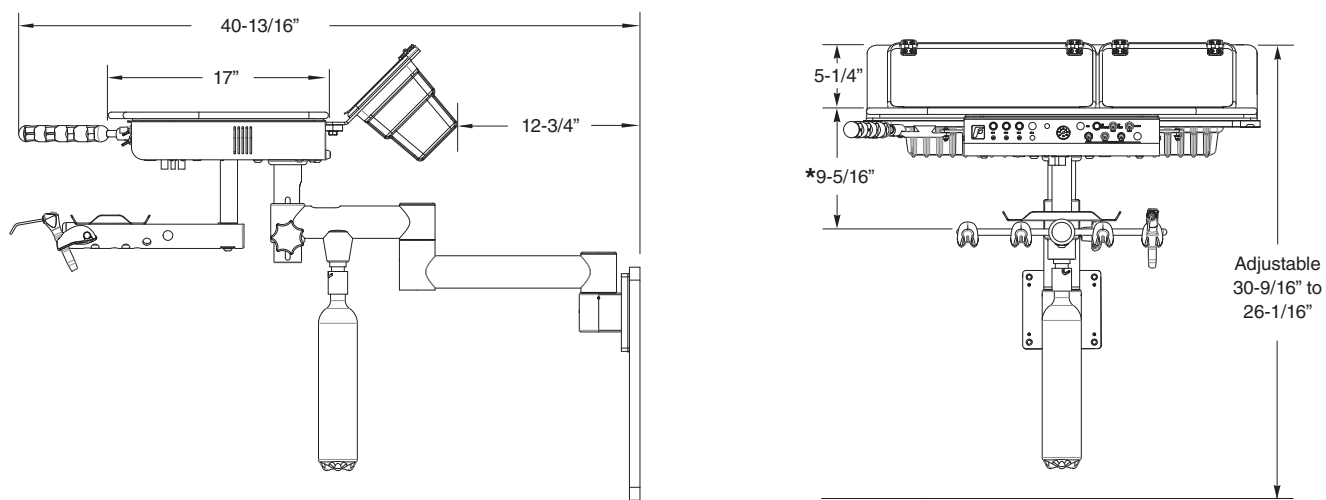

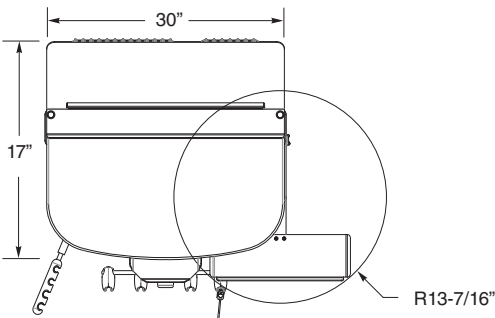
- Large contoured work surface standard for 7020/7021
- Work surface 9/16" thick
- Shown with organization features
- Height adjustable work surface, 5"



7020/7021 WM Mounting Plate

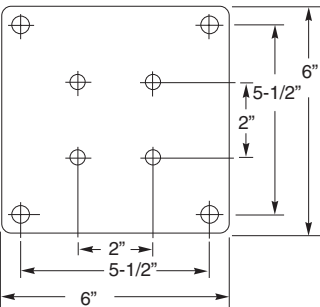



7020/7021 WM Leveling Mounting Board

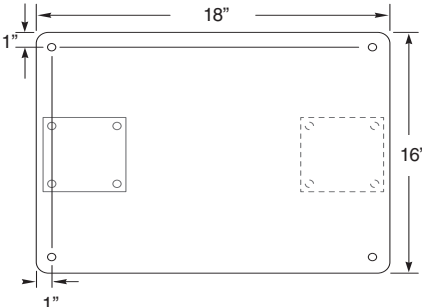



Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

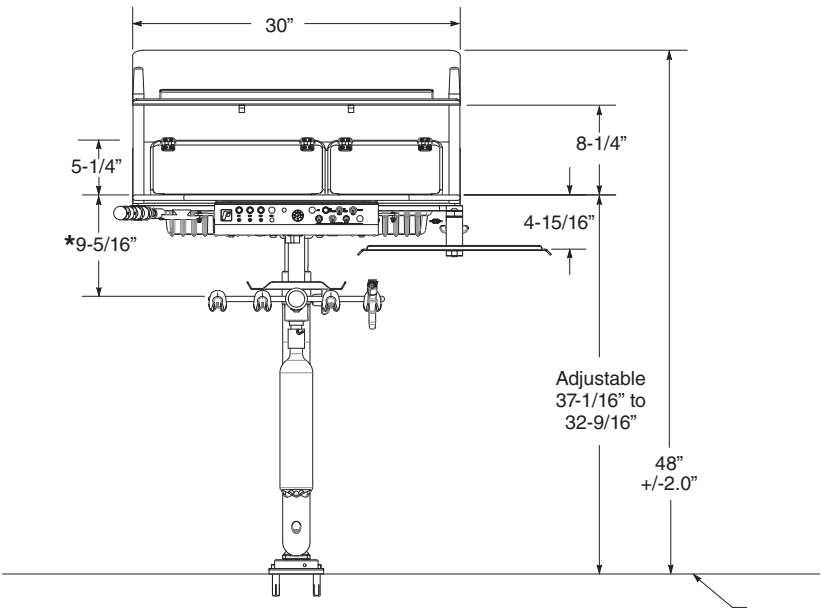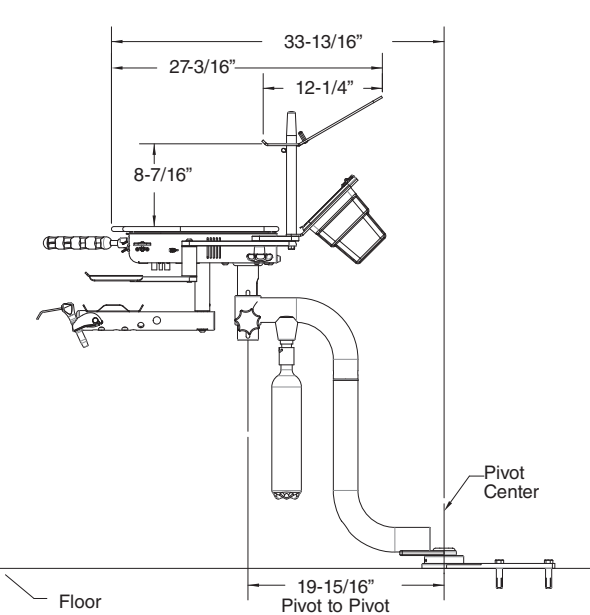7020PRO-UC Under Cabinet Floor Mount






7021PRO-UC Under Cabinet Floor Mount

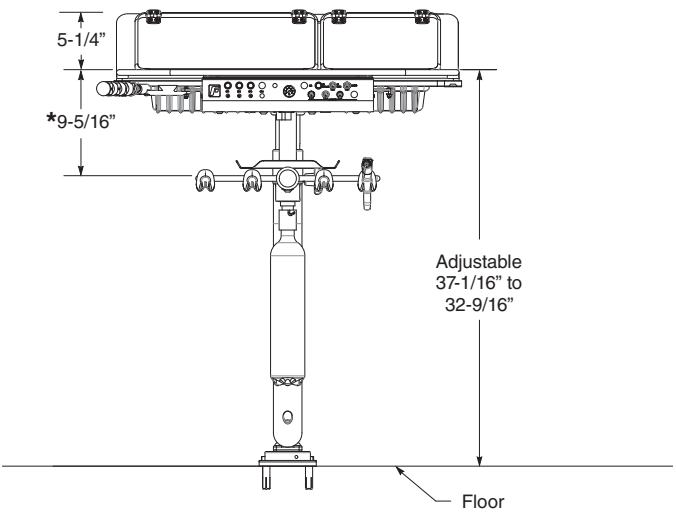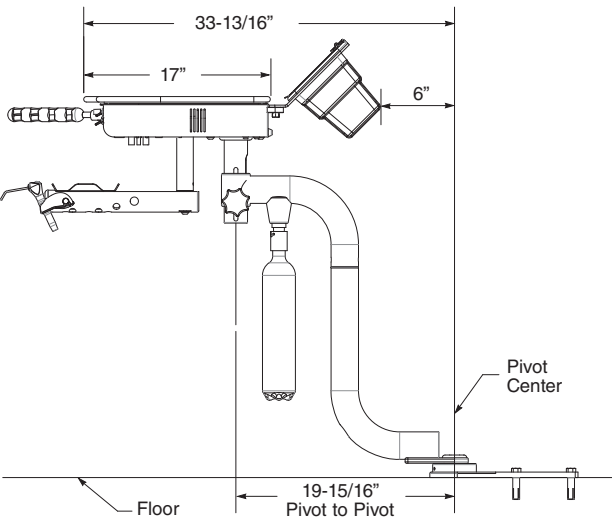

- Large contoured work surface standard for 7020/7021
- Work surface 9/16" thick
- Shown with organization features
- Height adjustable work surface, 5"

7020/7021PRO-UC Under Cabinet Floor Mount Leveling/Mounting Plate

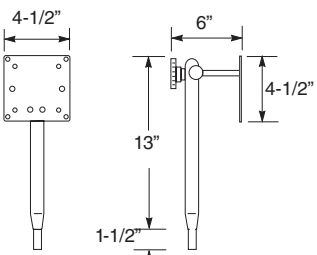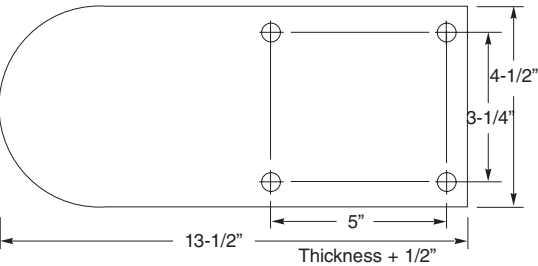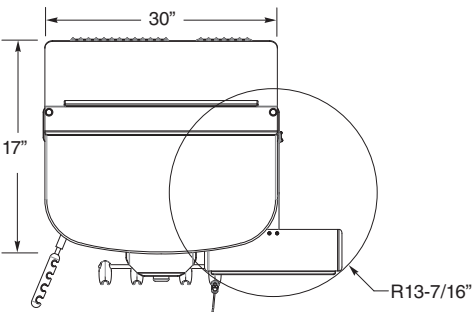1300-025-CS Monitor Assembly Option



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

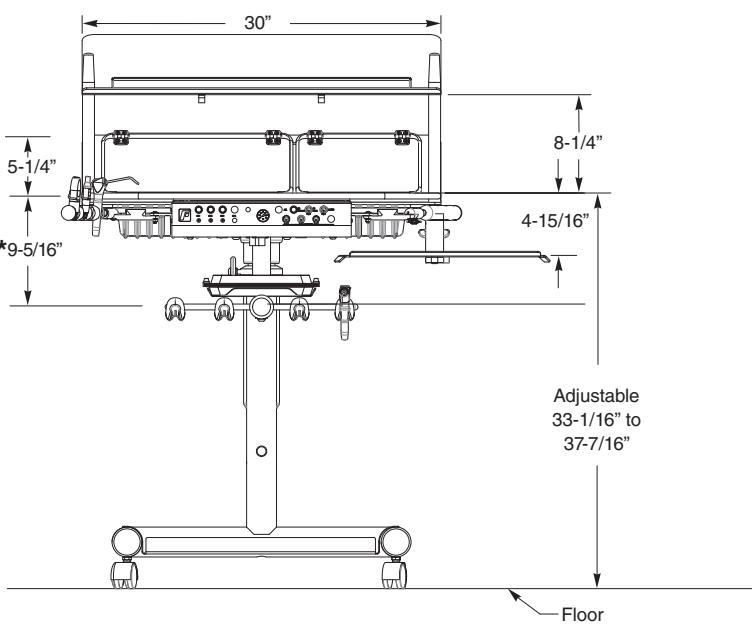
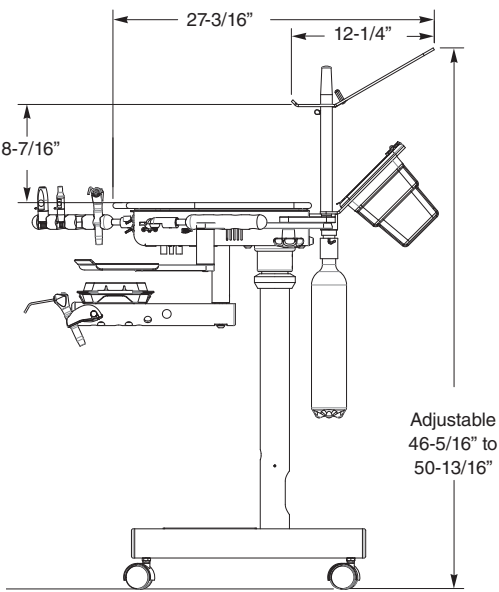
Factory
Options & Parts

General
Information



Warranty

Dimensions


7020PRO-MC Duo Cart




7021PRO-MC Duo Cart

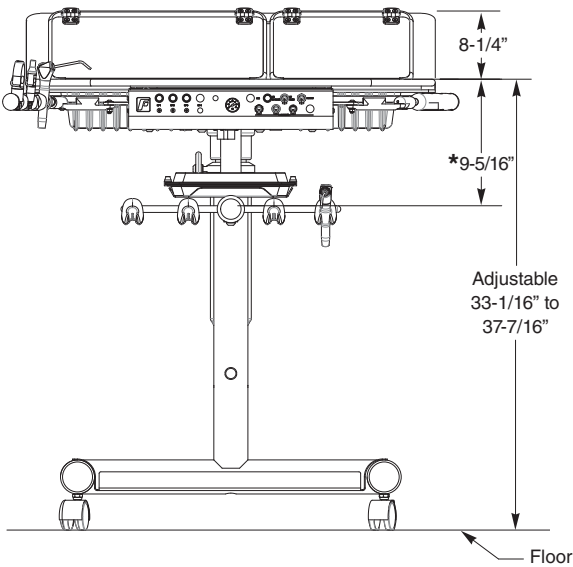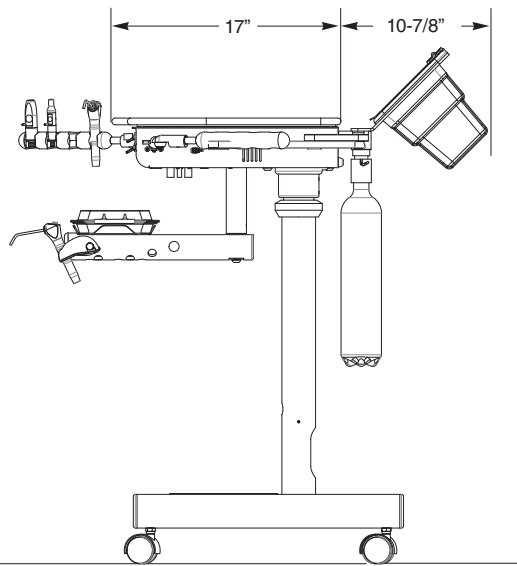

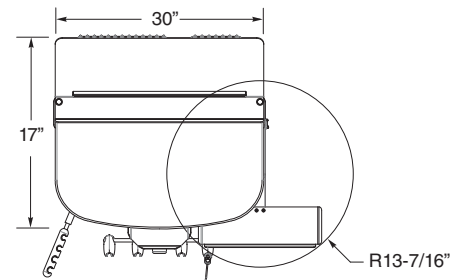- Large contoured work surface standard for 7020/7021
- Work surface 9/16" thick
- Shown with organization features
- Height adjustable work surface, 5"



1300-025-CS Monitor Assembly Option

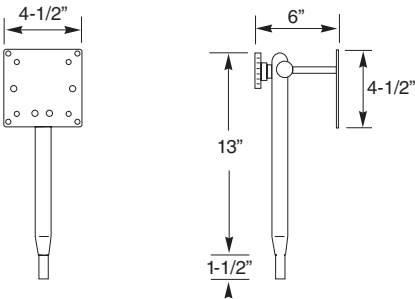

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

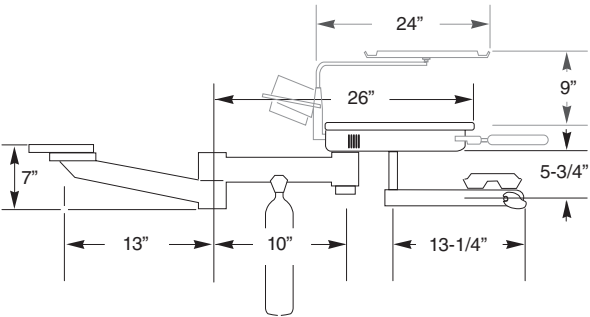
Factory
Options & Parts

General
Information

Warranty


Dimensions

7010BC Beneath Counter




- Work surface 9/16" thick (semi-round), 7010
- Shown with optional storage organizer and swivel tray
- Shown with additional optional asst's quad holder assembly
- Height adjustable work surface, 5"

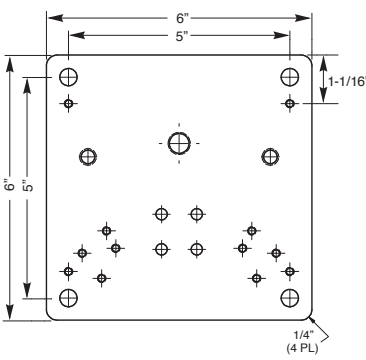
7010BC Mounting Plates




0010-086

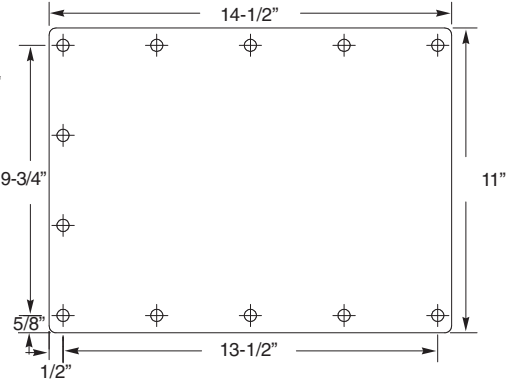


0010-084-CS Option

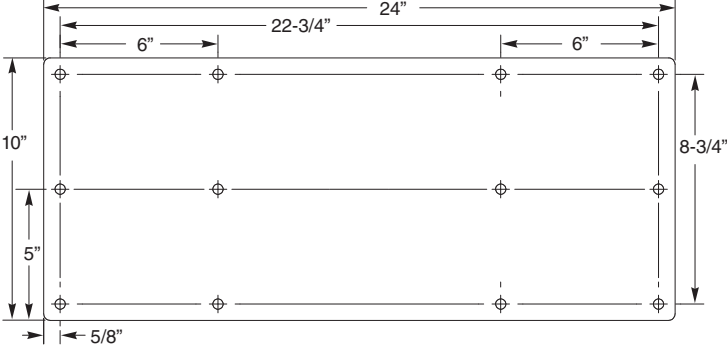


7010WM Mounting Plate

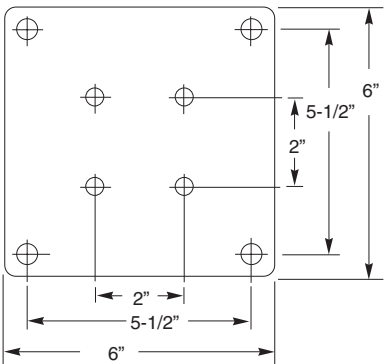


7010WM Leveling/Mounting Board

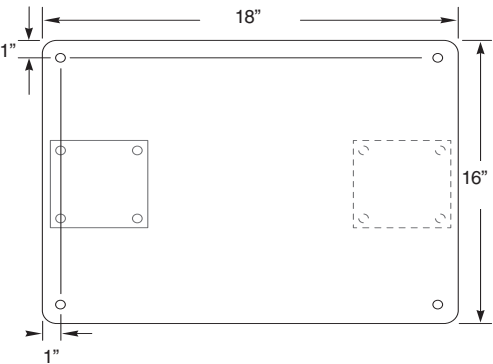

7010WM Wall Mount

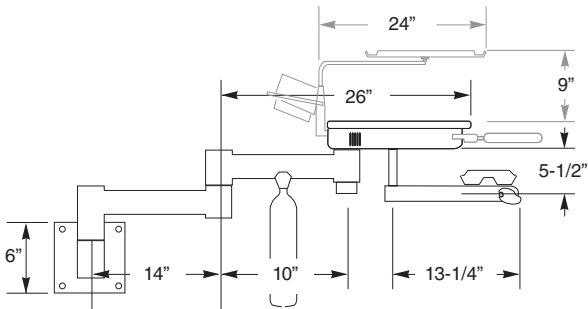


1300-025-CS Monitor Assembly Option

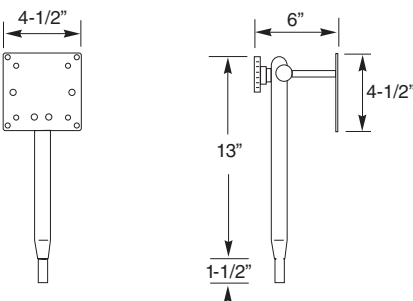


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

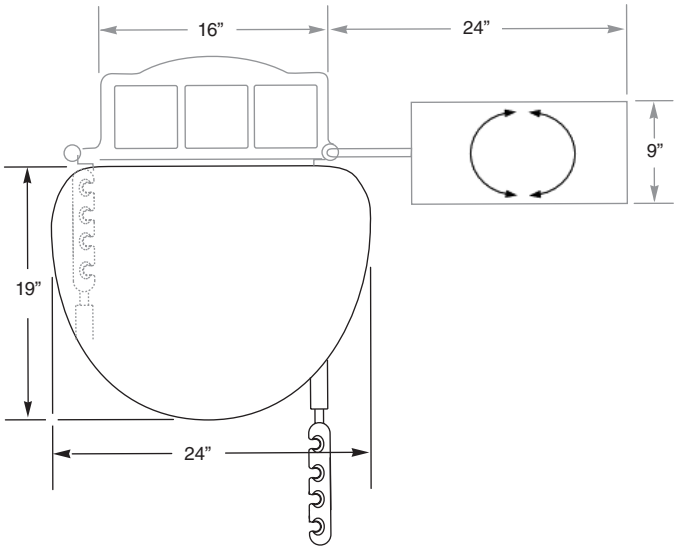
- Work surface 9/16" thick (semi-round), 7010
- Shown with optional storage organizer and swivel tray
- Shown with additional optional asst's quad holder assembly
- Height adjustable work surface, 5"




7010UC Under Cabinet Floor Mount

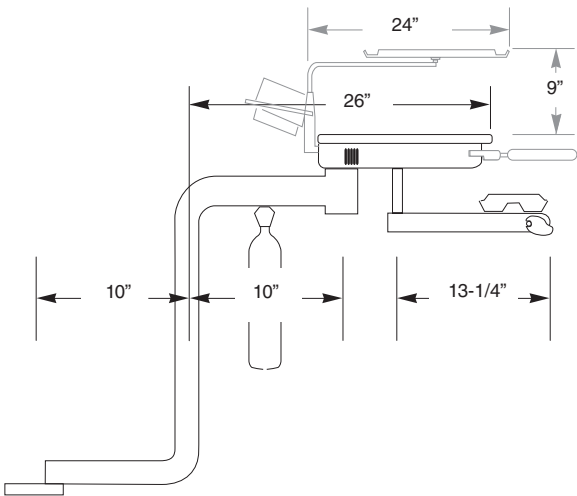


1300-025-CS Monitor Assembly Option

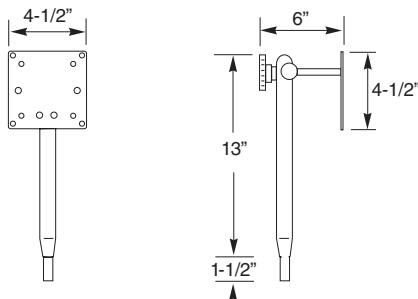

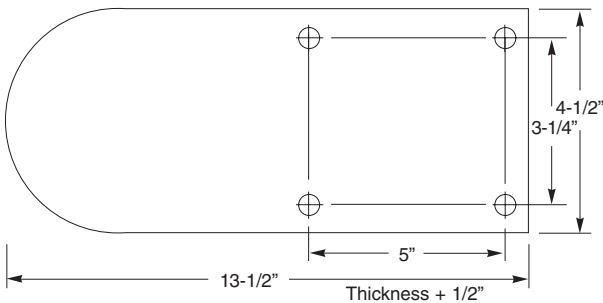7010UC Under Cabinet Floor Mount
Leveling/Mounting Plate



7010MC Duo Cart

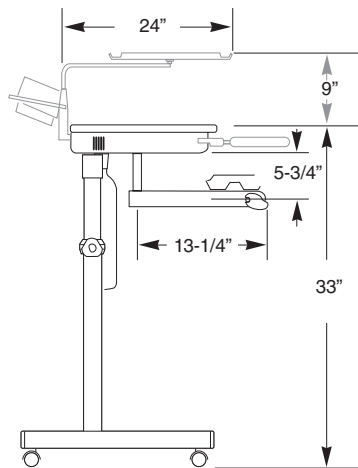

Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

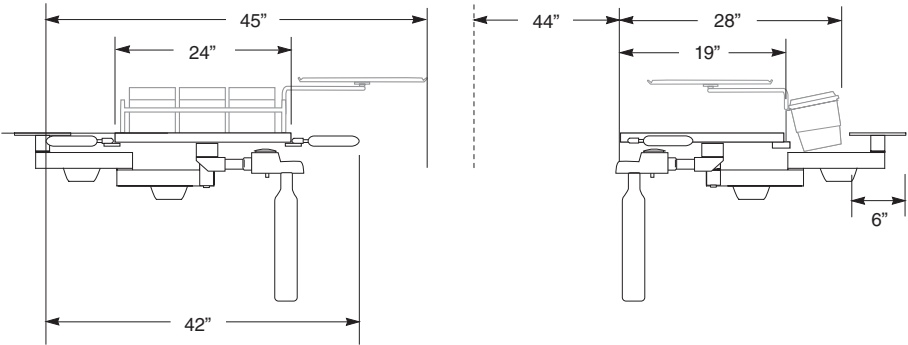
Factory
Options & Parts

General
Information


Warranty

Dimensions


3486BC and WM Rear Beneath Counter Mount




3486BC Heavy Duty Mounting Plate

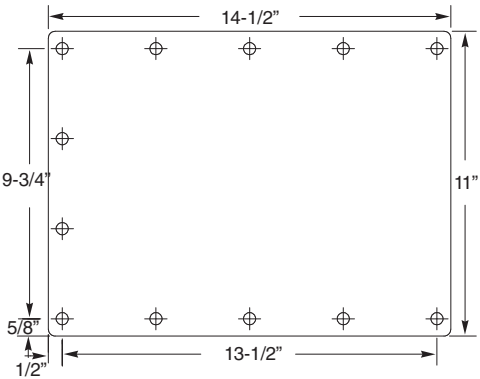


3486WM Wall Mount Hub

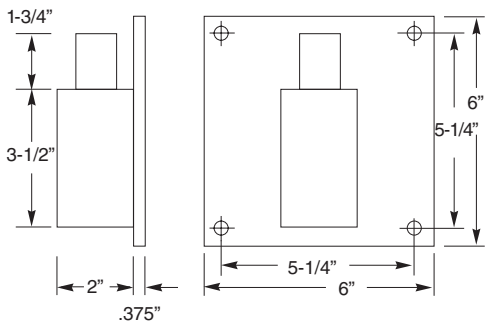


3486WM Leveling/Mounting Board

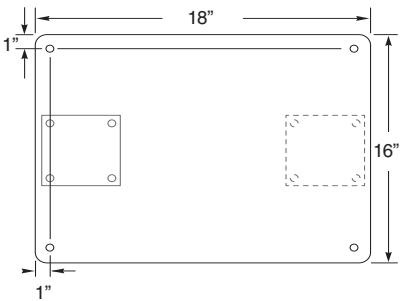


3486BC and WM Mounting Arm

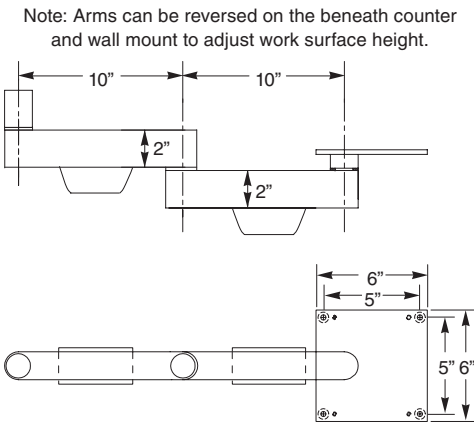

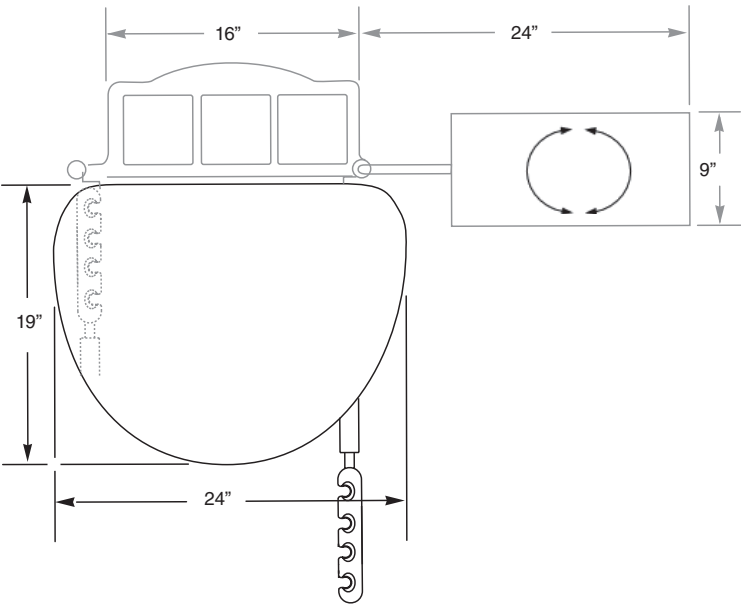- Semi-round work surface standard for 2184BC/2186BC/3486
- Work surface 9/16" thick (semi-round and rectangular)
- Shown with optional storage organizer and swivel tray
- Shown with additional optional asst's quad holder assembly




3486MC Mobile Cart

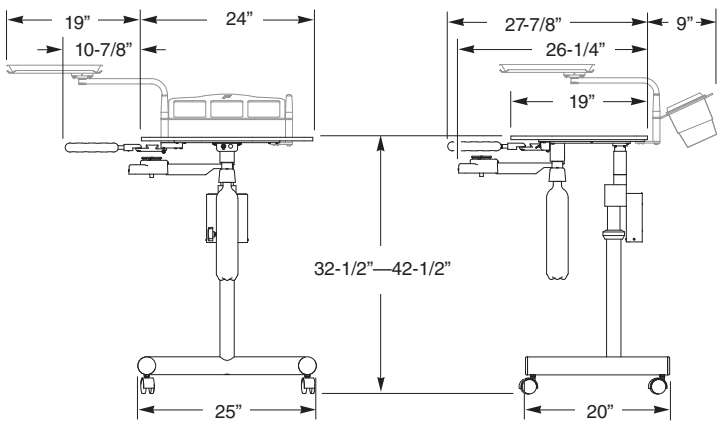


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

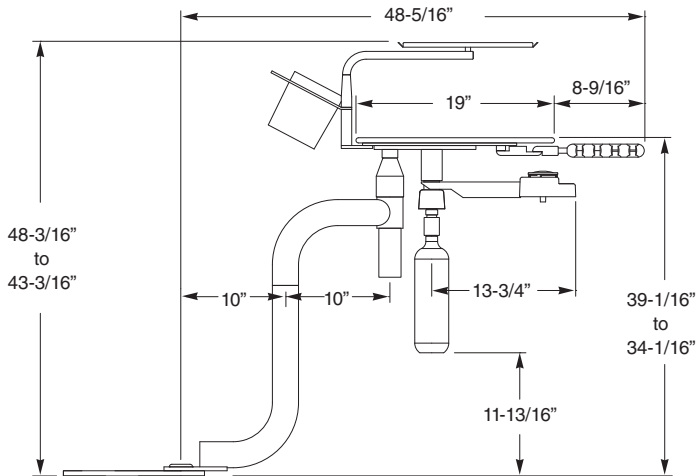
3486UC Under Cabinet Floor Mount Assistant's



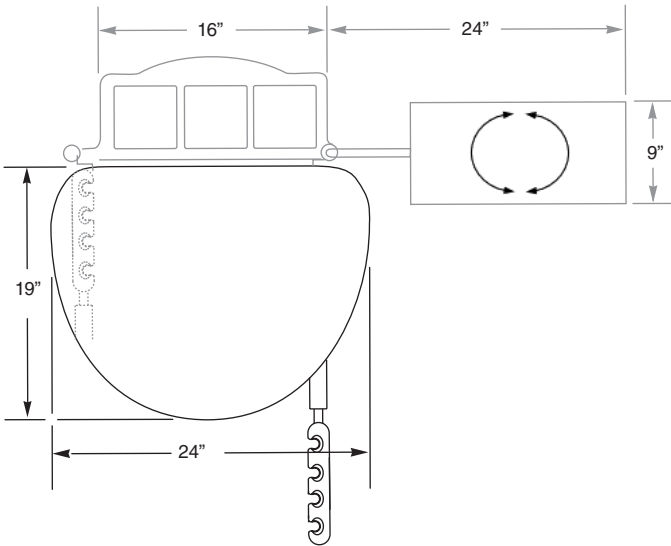
- Semi-round work surface standard for 2184BC/2186BC/3486
- Work surface 9/16" thick (semi-round and rectangular)
- Shown with optional storage organizer and swivel tray
- Shown with additional optional asst's quad holder assembly




2184 Rear Beneath Work Surface Duo Unit

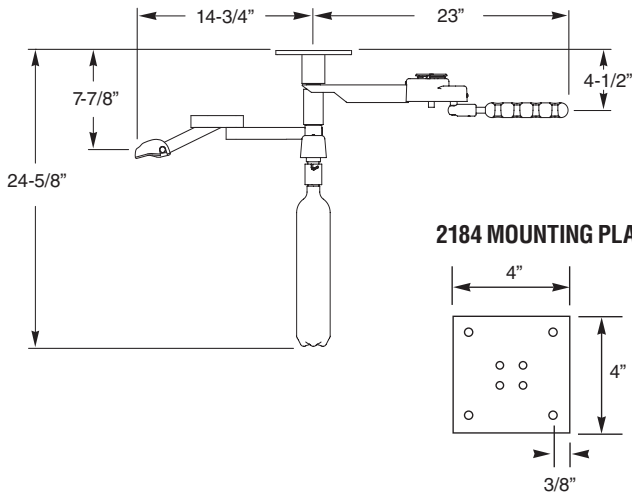

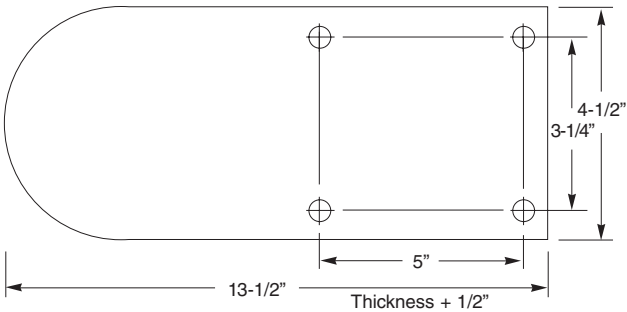
3486UC Under Cabinet Floor Mount Leveling/Mounting Plate




1300-025-CS Monitor Assembly Option

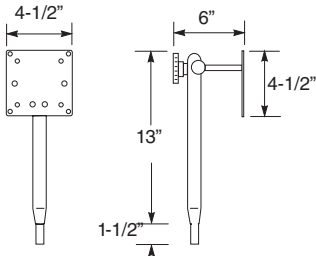

2184/2500 Wall Cut-Out

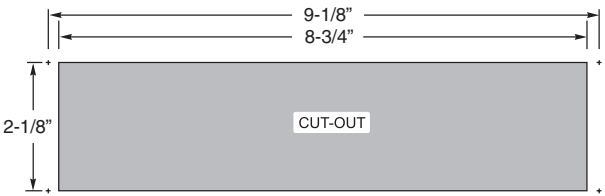


1700-142 SLIDE MOUNT OPTION, 1284/2184

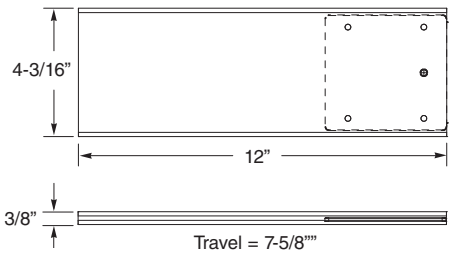


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

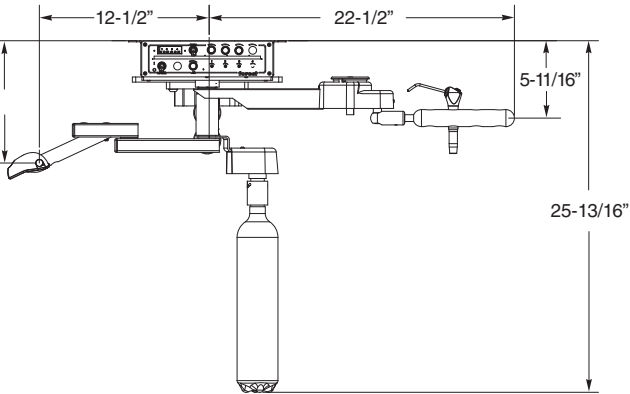
2186 Rear Beneath Work Surface Duo Unit




2186 Control Mounting

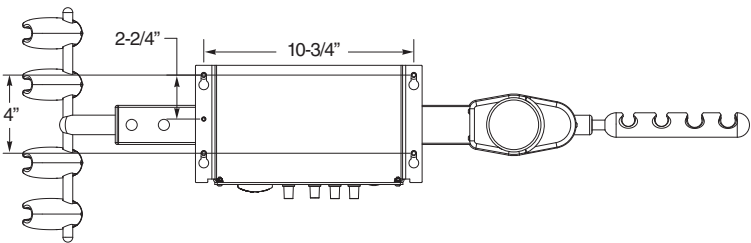


2184BC Beneath Counter Duo Unit

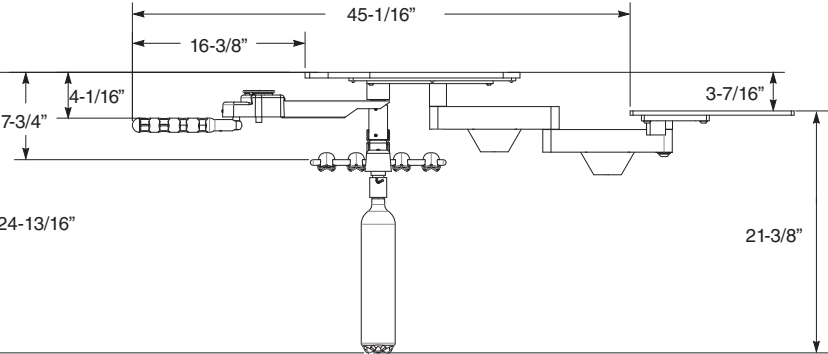


2184BC/2186BC Beneath Counter Duo Unit

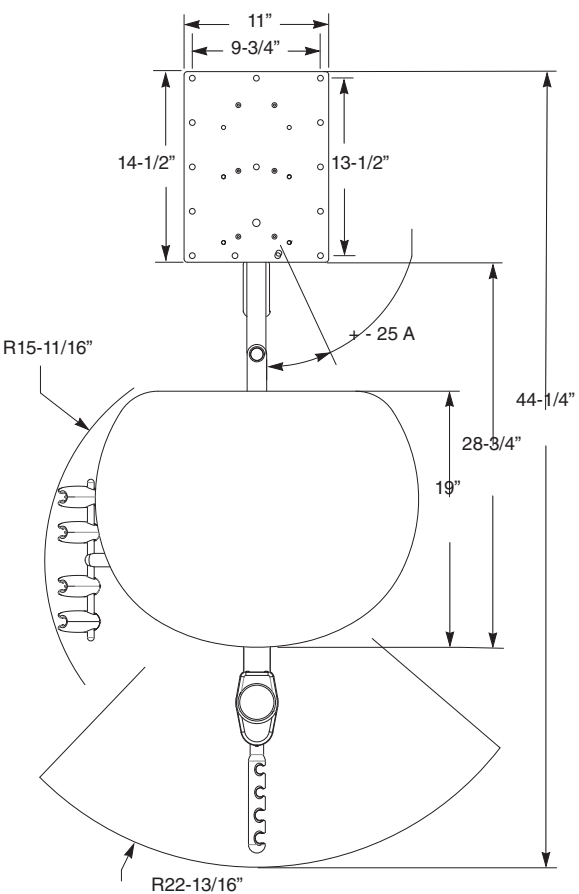


2186BC Rear Beneath Counter Duo Unit

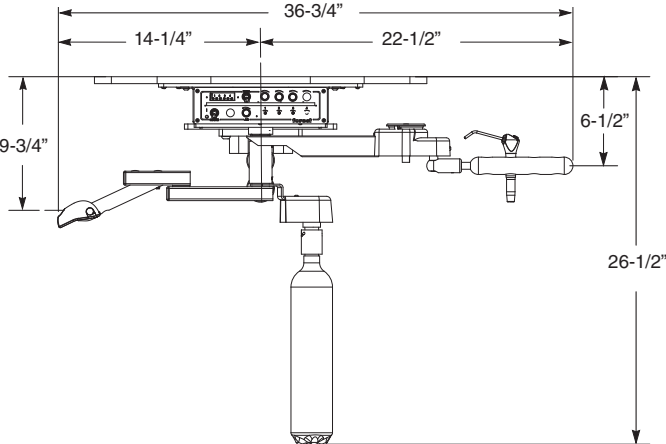


2500 PANEL MOUNT

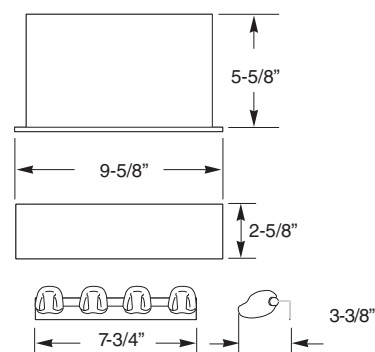

2520 PANEL MOUNT

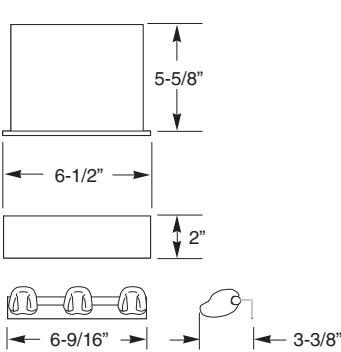


2520 Wall Cut-Out

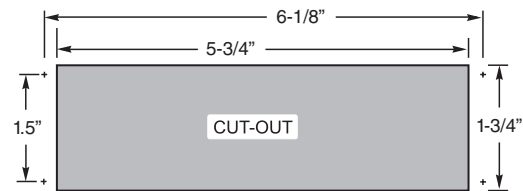


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

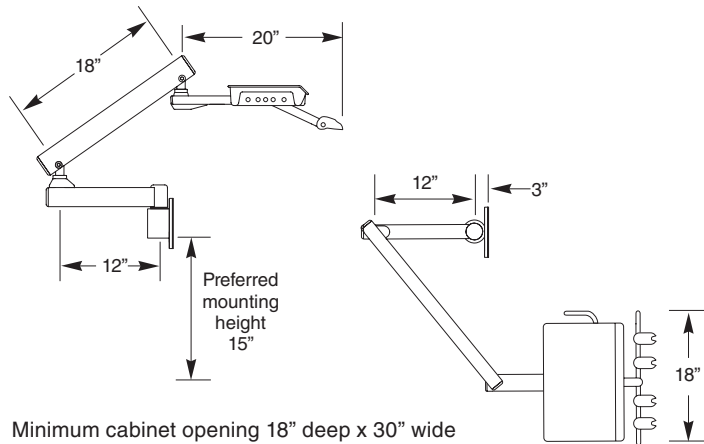
3235 Cabinet Side Delivery




3230 Wall Side Delivery

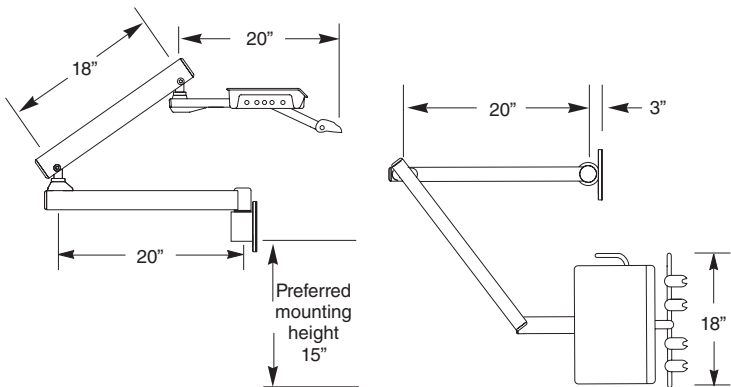


2400 Rear Mount Beneath Counter

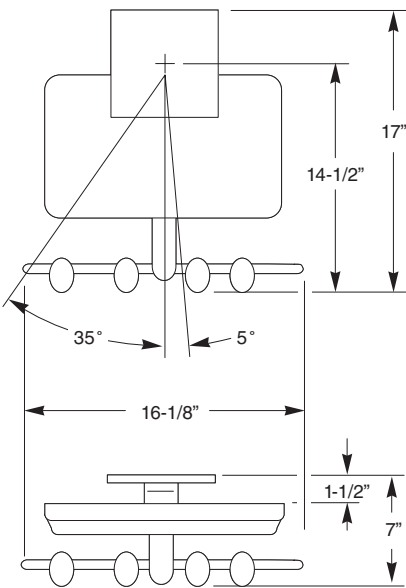


2400 Mounting Plate

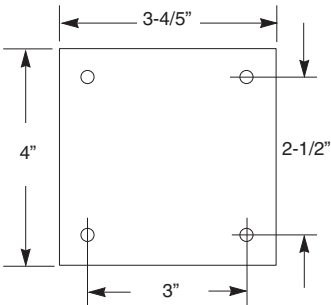


3230/3235 Leveling Mounting Plate

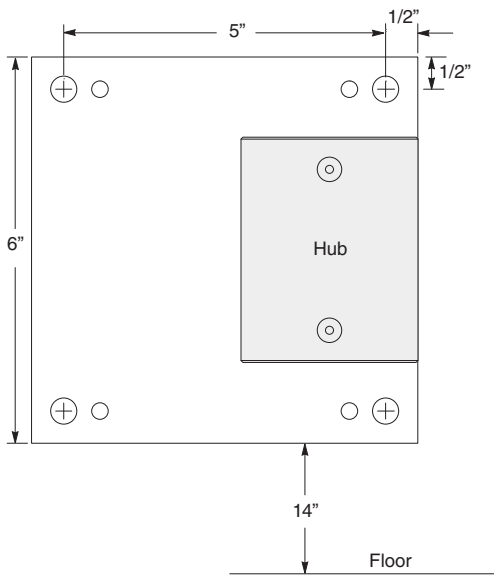


1110-115-CS
Folding Arm System Option

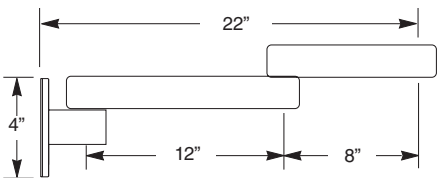


MT600/MT610
Compact Controls

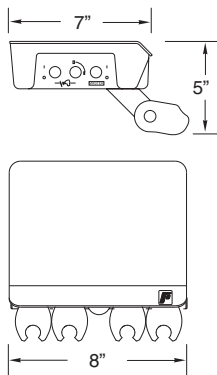

MT600/610BS
Leveling/Mounting Plate

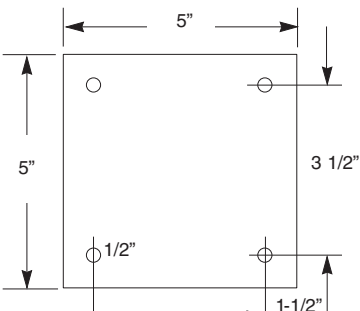


1700-140 Slide Mount
With Swivel Option

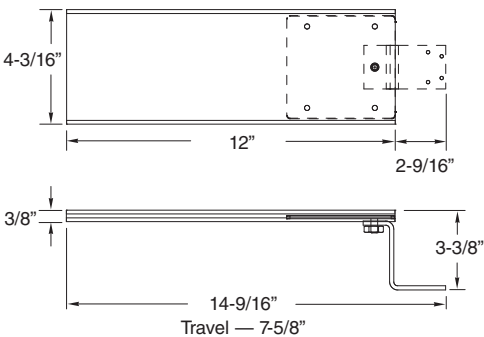


Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information
Warranty

Dimensions

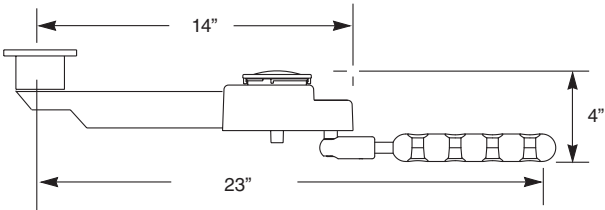
1284 Rear, Beneath Work Surface




1286 Rear, Beneath Shelf

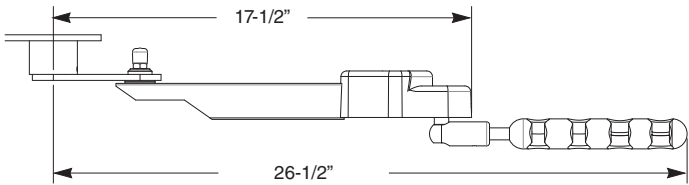


1700-142 Slide Mount Option, 1284/2184

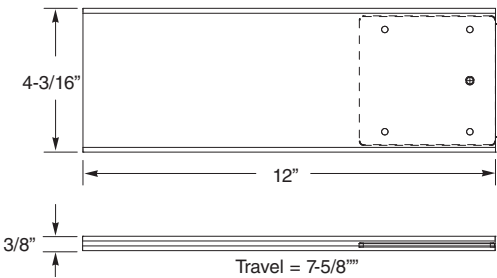


1284/1286 Mounting Plate

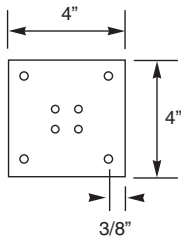

5350/5355 Standard IC Cart

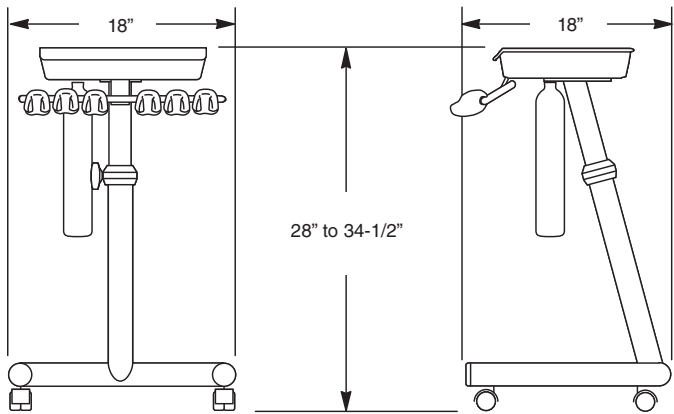

5350/5355 Work Surface

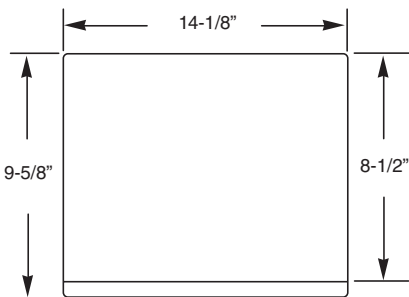

Proper maintenance will prolong the life of the dental equipment. It is important to adhere to routine practices for patient safety and avoid premature aging of the dental equipment. Following the step by step details outlined in the Forest Operator Guides included with the product will provide valuable information to achieve optimum success maintaining patient safety and the lifespan of the treatment room equipment.

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant's Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts

DENTALEZ and its employees are proud of the products we provide to the dental community. We stand behind these products with a warranty against defects in material and workmanship as provided below.

In the event that you experience difficulty with the application or operation of any of our products, please contact our technical service department at our expense at (800)423-3555.

If we cannot resolve the issue by telephone, we will arrange for a representative to contact you or suggest that the product be repaired by a dealer service technician.

If product return is required, we will provide you with a Return Authorization number and shipping instructions to return the product to the proper facility. If the product is under warranty we will ask you to provide proof of purchase such as a copy of your invoice. Please be sure to include the Return Authorization number on the package you are returning. It is the responsibility of the purchaser/end user to immediately cease using the product upon discovery that the product is not in proper working order and to dispose of any defective product that is not returned to DENTALEZ.

Freight costs for product returns are the responsibility of the customer. Products under warranty will be repaired or replaced, at our sole discretion, and returned at our expense. We are not responsible for shipping damages. We will, however, help you file a claim with the freight carrier.

DENTALEZ warrants all equipment and parts to be free of defects in material and workmanship, under normal usage, under the following terms:

Forest	Limited Warranty Period*
Patient chairs / Provider stools	5 years from date of installation
Chair / Stool upholstery	2 years from date of installation
Delivery systems	5 years from date of installation
7000 series delivery system chassis	20 years from date of installation
Handpiece illumination systems	6 months from date of installation
Power suppliers, Syringe buttons, HVE and SE valves	1 year from date of installation
Water bottles, gas springs, and handpiece tubing	2 years from date of installation
Air and Water tubing connectors (T's), Distribution block, Master control block, Flow valves	20 years from date of installation
Operatory lights	5 years from date of installation

Please note the following additional terms of our warranty and return policy:

- This warranty covers manufacturing defects only and does not cover defects resulting from abuse, accidents, misuse, improper handling, cleaning, care or maintenance. Damage, including, without limitation, discoloration of upholstery, resulting from contact, such as dye transfer, from garments worn by dental staff/team, the use of chemicals, cleaners, disinfectants or sterilization is not covered under this warranty nor is normal wear and tear. Failure to follow or observe our operating/use, maintenance and/or installation instructions voids this warranty as does the use of non-authorized parts and repairs made by an unauthorized repair facility.
- **Liability is limited to repair or replacement of the defective part(s) or products at our sole discretion. All other liabilities, in particular liability for damages, including, without limitation, special, exemplary, consequential or incidental damages are excluded.**
- THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. NO EMPLOYEE, REPRESENTATIVE OR DEALER IS AUTHORIZED TO CHANGE THIS WARRANTY IN ANY WAY OR TO GRANT ANY OTHER WARRANTY.

CANCELLATION POLICY:

Standard orders: Written cancellation request required within 3 days of order entry date, otherwise a 25% restock fee will be charged. Expedited ship date requests are non-cancellable.

CUSTOM ORDERS:

- Custom configurations are those that deviate from Forest’s published product offering; these orders are non-cancellable and non-returnable.
- May be subject to an additional “Deviation” charge and/or a non-refundable engineering charge.
- Custom powder coat, custom plastic surface treatments and/or custom treatments such as multiple color combinations or custom colors may result in an upcharge of 10-15% per model and are non-cancellable and non-returnable.

WARRANTY REPLACEMENTS:

Parts replaced on a product that is in warranty will be warranted for the duration of that product’s original warranty

PRODUCT RETURN:

Opened products or product returns more than a year old cannot be returned for credit. There will be a 25% (\$25.00 minimum) restocking charge on all items authorized for return.

Ultraleather is a registered trademark of Ultrafabrics, Inc., Formica is a registered trademark of The Diller Corporation, Corian is a registered trademark of DuPont, Wilsonart is a registered trademark of Wilsonart LLC., Meganite is a registered trademark of Lai, James, Naugahyde is a registered trademark of United States Rubber Products, Inc., and Active tilt is a registered trademark of Crown Seating. Designer Friendly is a registered trademark of DTE Oregon, Inc.

Pivot & Fixed Chair Mount Packages
Pivot & Fixed Chair Mount Package Ordering Guide
Fusion Packages
Fusion Package Ordering Guide
Designer Friendly Colors & Surfaces
Chairs & Stools
Lights and Monitor Mounts
Rear Mount Side Delivery Systems
Ancillary Ordering & Options
Control Head Accessories
Simplified Cuspidors & Assistant’s Instrumentation
Simplified Rear & Side Delivery Systems
Factory Options & Parts
General Information

Heritage

EQUIPMENT

EMBRACING OUR TRADITION

DENTALEZ has been providing high-quality equipment to the dental industry for over 60 years. The DENTALEZ Heritage product line continues the legacy of DENTALEZ equipment in a dedicated online catalog while supplies last. The future of DENTALEZ equipment will be re-imagined through the lens of the Forest equipment portfolio.



Pivot & Fixed
Chair Mount
Packages

Pivot & Fixed
Chair Mount
Package
Ordering Guide

Fusion
Packages

Fusion
Package
Ordering Guide

Designer
Friendly
Colors &
Surfaces

Chairs
& Stools

Lights and
Monitor Mounts

Rear Mount
Side Delivery
Systems

Ancillary
Ordering
& Options

Control Head
Accessories

Simplified
Cuspidors &
Assistant's
Instrumentation

Simplified
Rear & Side
Delivery
Systems

Factory
Options & Parts

General
Information

Warranty

[DENTALEZ.COM/HERITAGE](https://www.dentalez.com/heritage)

Contact your local representative for more information.

	Pivot & Fixed Chair Mount Packages
	Pivot & Fixed Chair Mount Package Ordering Guide
	Fusion Packages
	Fusion Package Ordering Guide
	Designer Friendly Colors & Surfaces
	Chairs & Stools
	Lights and Monitor Mounts
	Rear Mount Side Delivery Systems
	Ancillary Ordering & Options
	Control Head Accessories
	Simplified Cuspidors & Assistant's Instrumentation
	Simplified Rear & Side Delivery Systems
<p>Visit DSXsales.com to find your local representative. For DENTALEZ, Forest, Ramvac, and Imagine + orders and/or questions please call our Customer Experience Representatives at 866.DTE.INFO.</p> <p>Regarding Radiance cabinet customer service, please contact ILS at ilsdental.com/contact.</p> <p>2024 DENTALEZ, Inc. All rights reserved. DENTALEZ, Star and Dentoform are registered trademarks and Imagine + and Nevin are trademarks of DENTALEZ INC. Ramvac is a registered trademark of Ramvac Dental Products Inc.</p>	
General Information	
Warranty	